



Experimental facilities

Heinz Maier-Leibnitz Zentrum (MLZ)

MLZ is a cooperation between:

Bayerisches Staatsministerium für
Bildung und Kultus, Wissenschaft und Kunst



SPONSORED BY THE



Federal Ministry
of Education
and Research

The Heinz Maier-Leibnitz Zentrum (MLZ):

The Heinz Maier-Leibnitz Zentrum is a leading centre for cutting-edge research with neutrons and positrons. Operating as a user facility, the MLZ offers a unique suite of high-performance neutron scattering instruments. This cooperation involves the Technische Universität München, the Forschungszentrum Jülich and the Helmholtz-Zentrum Geesthacht. The MLZ is funded by the German Federal Ministry of Education and Research, together with the Bavarian State Ministry of Education, Science and the Arts and the partners of the cooperation.

The Forschungs-Neutronenquelle Heinz-Maier-Leibnitz (FRM II):

The Forschungs-Neutronenquelle Heinz-Maier-Leibnitz provides neutron beams for the scientific experiments at the MLZ. The FRM II is operated by the Technische Universität München and is funded by the Bavarian State Ministry of Education, Science and the Arts.

Welcome to the MLZ!

The Heinz Maier-Leibnitz Zentrum (MLZ) offers beam time at high-class instruments hosted at the Research Neutron Source FRM II at Garching. The source provides the complete spectrum of neutrons: cold, thermal, hot neutrons as well as fast neutrons from the converter facility and in addition, the world's highest flux of mono-energetic positrons. An ultra cold source is under construction. Structure research, large scale structures, spectroscopy, imaging and analysis as well as particle physics are topics our instrument suite is equipped for.

Continuous improvements on our instrument suite make it necessary to revise this brochure regularly - in your hands you are holding the third edition! Find all details about our instruments, laboratories, and sample environment here, get in touch with our scientists, discuss your ideas, and become part of our user community: Just submit your proposal to the next deadline after which it will be reviewed by our international referees based solely on the scientific merit of the proposed experiment. The mayor part of the available beam time is distributed to external users.

On February 21st, 2013, the MLZ was inaugurated in order to present a cooperation that had been established in 2011. The partners are the TUM and neutron centres of the Helmholtz Association in Jülich and Geesthacht. They receive additional funds for the cooperation from the Federal Ministry of Education and Research. The source FRM II itself is funded by the Bavarian State Ministry of Sciences, Research and the Arts via the Technische Universität München (TUM).

At the MLZ, not only the cooperation partners assemble, but also institutes of the Max-Planck-Gesellschaft and many university institutes, merging their knowledge and substantial resources with the intention to serve the German and international neutron community the best. A list of them is given at this brochure's end.

The MLZ looks forward to your compelling ideas and experiments!

Yours,



Prof. Dr. Winfried Petry

Scientific Director MLZ
Scientific Director FRM II
and
Chair of Institute for Functional Materials
Technische Universität München



Prof. Dr. Thomas Brückel

Scientific Director MLZ
Director of the Jülich Centre for
Neutron Science (JCNS)
and
Peter Grünberg Institute (PGI)
Forschungszentrum Jülich GmbH

Neutron Source

The neutron source FRM II.....	8
Secondary neutron sources.....	10
Neutron guides.....	12

Structure Research

BIODIFF diffractometer for large unit cells.....	16
HEIDI single crystal diffractometer on hot source.....	18
MIRA dual wavelength band instrument.....	20
POLI polarised hot neutron diffractometer.....	22
RESI thermal neutron single crystal diffractometer.....	24
SPODI high resolution powder diffractometer.....	26
STRESS-SPEC materials science diffractometer.....	28

Large Scale Structures

KWS-1 small angle scattering diffractometer.....	32
KWS-2 small angle scattering diffractometer.....	34
KWS-3 very small angle scattering diffractometer with focusing mirror.....	36
MARIA magnetic reflectometer with high incident angle.....	38
NREX neutron reflectometer with X-ray option.....	40
REFSANS horizontal time-of-flight reflectometer with GISANS option.....	42
SANS-1 small angle neutron scattering.....	44

Spectroscopy

DNS diffuse scattering neutron time-of-flight spectrometer.....	48
J-NSE neutron spin echo spectrometer.....	50
PANDA cold three axes spectrometer.....	52
PUMA thermal three axes spectrometer.....	54
RESEDA resonance spin echo spectrometer.....	56
SPHERES backscattering spectrometer.....	58
TOFTOF cold neutron time-of-flight spectrometer.....	60
TRISP three axes spin echo spectrometer.....	62

Imaging & Analysis

ANTARES cold neutron radiography and tomography station.....	66
MEDAPP and NECTAR fission neutron facility for scientific, medical, and industrial applications.....	68
PGAA prompt gamma and in-beam neutron activation analysis facility.....	70

Particle Physics

EDM neutron electric dipole moment measurement	74
MEPHISTO facility for particle physics with cold neutrons	76

Positrons

NEPOMUC neutron induced positron source Munich	80
CDBS coincident Doppler-broadening spectrometer.....	82
PAES positron annihilation induced Auger-electron spectrometer	83
PLEPS pulsed low energy positron system.....	84
SPM scanning positron microscope	85

Irradiation facilities

Irradiation facilities.....	88
-----------------------------	----

Sample environment (SE), Laboratories (LAB), and User facilities (UF)

SE	
magnetic field	94
high pressure.....	97
specialised equipment.....	98
low temperatures	100
high temperatures.....	106
LAB	
Biology Laboratory.....	108
Chemistry Laboratory	108
Materials Science Laboratory.....	109
Sample Preparation Laboratory.....	109
UF	
TEM: JEOL 200 kV JEM-FS2200.....	110
Thin Film Lab: MBE	110
X-ray Tomography: v tome x.....	111

User Office

From idea to publication - General information for users	114
Access to the MLZ.....	116
Partner institutions.....	118
Imprint.....	122



The neutron source FRM II



Secondary neutron sources



Neutron guides

Neutron Source

The neutron source FRM II



Powerful neutron source

The neutron source Heinz Maier-Leibnitz (FRM II) is a multipurpose research reactor with a particular focus on beam tube experiments. It has been designed for providing neutrons to scientific users from all over the world as well as for medical and industrial applications. The FRM II is operated as a Corporate Research Center by the Technische Universität München (TUM) in Garching near Munich, Germany. Its first criticality was achieved in March 2004.

The FRM II is the most powerful neutron source in Germany and reaches the highest neutron flux density ($8 \cdot 10^{14} \text{ n/cm}^2\text{s}^{-1}$, max. undisturbed flux density) relative to its thermal power (20 MW) throughout the world. More than 30 experimental facilities are operated by scientific teams from German universities, research institutes of the Helmholtz Association, and the Max Planck Society at the neutron source. Today, 26 beam tube facilities are operational. Further seven irradiation facilities mainly for medical and industrial application are in service, an irradiation facility for the production of the medical isotope ^{99}Mo is under construction. The FRM II is equipped with cold, thermal, hot, and fast neutron sources and thus covers a broad range of applications, including a device for the continuous production of an intense positron beam.

Fuel element

The FRM II was designed for an exclusive purpose: the production of neutrons for basic research and applied physics. Its high performance is based on

the concept of a compact core: a single, cylindrical fuel element with a diameter of just 24 cm is sufficient for 60 days of reactor operation. The fuel zone measures 70 cm in height and contains about 8 kg of uranium in the form of U_3Si_2 . Like other high-performance neutron sources around the world, the FRM II uses highly enriched uranium.

The fuel element is located in the centre of a moderator tank filled with heavy water (D_2O). The tips of the beam tubes are placed in the region of the maximum thermal neutron flux density. Various vertical irradiation channels are arranged in the moderator tank. The beam tubes guide the neutrons to the experiments in the Experimental Hall and Neutron Guide Hall West.

The Experimental Hall hosts scientific instruments exhibiting high neutron flux densities and gives access to the positron beam lines, whereas the Neutron Guide Hall West is connected via six neutron guides to the cold neutron source. A second guide hall will be connected to the reactor building soon in order to extend the number of available instruments.

Safety first

The highest priority is always given to safety at the FRM II. The inherent safety stems from its principle design, with a compact fuel element built into the centre of the moderator tank filled with heavy water. Fuel element, moderator tank, and beam tubes are located in the reactor pool filled with 700 m³ of highly purified water. While passing the fuel element, the temperature of the cooling water only increases from 37°C to a maximum of about 53°C. Neither steam nor high pressures are produced. Three subsequent cooling circuits guarantee the safe dissipation of the 20 MW via the air path.

Redundant safety installations (i.e. multiple, independently constructed units) are a key feature of the safety concept of the FRM II. The single central control rod inside the fuel element is sufficient to regulate and shut down the reactor. Additionally, a redundant set of five shut-down rods is available. Each of these systems is constructed such that the

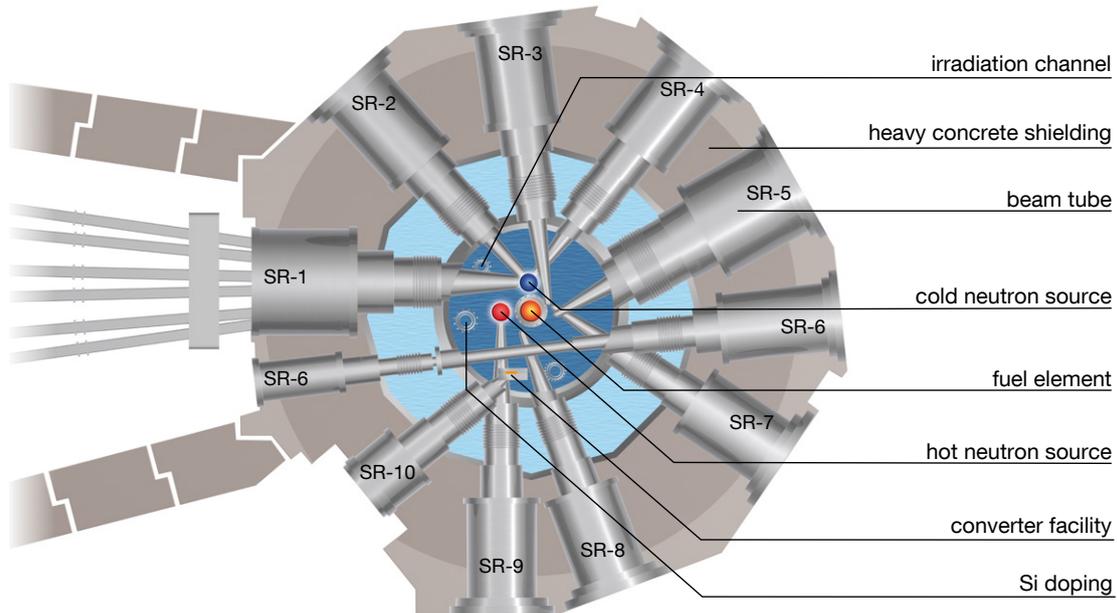


Figure 1: Horizontal section of the reactor pool showing the beam tubes, fuel element, as well as cold and hot neutron source. Beam tubes SR-1, -2, and -4 are fed by the cold source, SR-9 by the hot source. The through-going beam tube SR-6 will be used by the ultra cold neutron source. The converter plate for fast neutrons supplies beam tube SR-10 with the tumour treatment facility. The remaining beam tubes are placed into the highest neutron flux taking up the thermal neutrons.

reactor can be shut down in a fast and durable manner, completely independently.

The 1.8 m thick outer concrete wall of the reactor building protects the reactor against all impacts from outside. It was designed to resist the crash of a fast military as well as a civilian aircraft. This feature has

been approved by independent experts. Furthermore, the building can withstand earthquakes up to level VI ½ (MSK), which is beyond the strength of possible earthquakes in the region, or a high flood-water from the nearby river Isar with a height that might occur once every 10,000 years.

Technical Data

Reactor main parameters

- Thermal power: 20 MW
- Max. undisturbed thermal neutron flux density: $8 \cdot 10^{14} \text{ n/cm}^2 \text{ s}^{-1}$
- 10 horizontal; 2 tilted beam tubes
- D₂O moderator
- H₂O cooling water

Staff and money

- 435 Mio€ construction costs (2003)
- ~ 400 employees on-site (MLZ + FRM II)

Instruments

- 26 instruments in routine operation (2015)
- 6 instruments under construction

Fuel element

- Dimensions:
height: 133 cm; outer diameter: 24 cm;
active zone: 70 cm
- 8 kg HEU in U₃Si₂ alloy in 113 fuel plates
- Enrichment: 92.5% in ²³⁵U
- 60 days in a row/ fuel element – typical 240 days of operation per year

Technical Director Dr. Anton Kastenmüller

Email: anton.kastenmueller@tum.de
Office: +49.(0)89.289.12154



www.frm2.tum.de/en/the-neutron-source

Secondary neutron sources



Figure 1: The cold source vessel surrounded by the beam tubes SR-1, -2, -4. View into the moderator tank of the FRM II before operation.

The different instruments at the neutron source FRM II are supplied by various secondary sources slowing down or speeding up the neutrons after thermalisation in the moderator tank. This allows a large variety of applications. Furthermore, beam tube SR-11 hosts the positron source NEPOMUC, which is described in the chapter *Positrons*.

The cold source

At FRM II, a major part of the experiments is carried out using low energy (“cold”) neutrons. For this purpose the beam tubes SR-1, -2, and -4 are fed from a special secondary source, the so-called cold neutron source, which is located within the moderator tank close to the maximum of the thermal neutron flux density.

The major component of the cold neutron source is the moderator chamber containing about 12 l of liquid deuterium (D_2) at a temperature of approximately 25 K. Neutrons interacting with the liquid D_2 in the moderator chamber are slowed down at this low ambient. Consequently their energy spectrum is shifted considerably into the low energy range (see fig. 2).

Within the moderator tank, the deuterium is available in a closed circuit driven by natural convection. The liquefaction of the deuterium is achieved by means of a He refrigerator acting as the heat sink of a heat exchanger which is installed above the

moderator chamber. The liquefied D_2 rinses into the moderator chamber where it is evaporated again due to neutron moderation and gamma heating.

The moderator chamber and the heat exchanger form the inpile section of the cold neutron source. Important components outside the moderator tank are the refrigerator, a buffer tank and a metal hydride storage unit used to keep the deuterium during periods when the cold neutron source is out of operation and warmed to room temperature e.g. during maintenance periods of the reactor. The subsequent preparation of the cold source for low temperature operation takes approximately one week.

The energy distribution of the neutrons moderated by the cold neutron source is shown in fig. 2. It has its maximum at 40 meV corresponding to a wavelength of 1.4 Å. From the figure it is clearly visible that the corresponding intensity of neutrons exhibiting very long wavelengths is considerably increased as compared to the thermal spectrum generated in the moderator tank at an ambient temperature of about 70°C.

The hot source

Neutrons of short wavelengths in the range of 0.1 eV to 1 eV are used to investigate the structure of condensed matter. As only a small fraction of this spectral range is present in the thermal neutron distribution, the neutrons are moderated upwards from about 300 K to 2200 K. This spectrum shift is performed by the hot neutron source.

The hot moderator consists of a graphite block thermally insulated and positioned in the moderator tank next to the maximum thermal neutron flux density. The graphite cylinder is heated by gamma radiation from the reactor. A double-skinned zircaloy container with interposed insulating graphite felt insulates the graphite block, ensuring a secure containment of the hot graphite. At a reactor power of 20 MW, the temperature inside the container rises to about 2000°C. The hot source provides neutrons to beam tube SR-9, which supplies the single crystal diffractometers HEiDi and POLI.

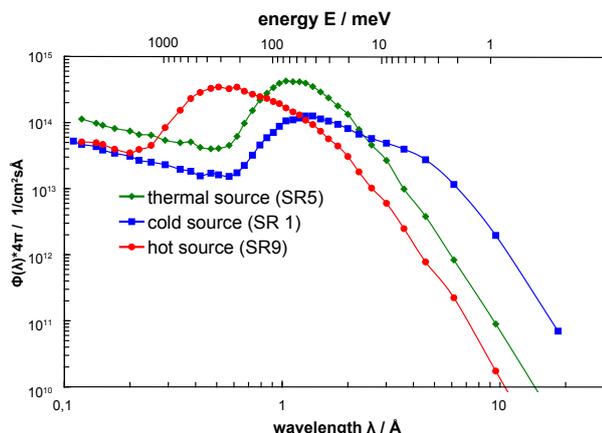


Figure 2: Neutron spectra at the entrance of the beam tubes at 20 MW reactor power.

Converter facility for fast neutrons

In order to obtain a high-intensity neutron beam with an unmoderated fission spectrum, an arrangement of uranium plates is inserted as a secondary source (so-called converter) in front of the tip of beam tube SR-10. It supplies the tumour treatment facility MEDAPP and the radiography and tomography station NECTAR with fast neutrons. Slow neutrons induce nuclear fission in the uranium plates causing the emission of neutrons with an average energy of 1.9 MeV. Without moderation, the fast neutrons are led through a horizontal beam tube to the experiments. The two converter plates deliver a thermal power of about 80 kW.

Technical Data

Cold source

- T = 25 K
- Volume of moderator chamber: 25 l
- Volume of liquid D₂: 12 l
- Mass of D₂ in cold source: 2.4 kg
- Cold moderator pressure: 1.5 bar
- 18 instruments fed by cold source

Hot source

- Mass of graphite in hot source: 14 kg
- Temperature at 20 MW reactor power: ~ 2000°C



Figure 3: The moderator of the hot source in the reactor pool of the FRM II during installation.

Converter facility

- Mass of uranium in converter plates: 540 g
- Thermal power: 80 kW
- Degree of enrichment: 92.5% in ²³⁵U


Dietmar Päthe

Email: dietmar.paethe@frm2.tum.de
Office: +49.(0)89.289.12122



www.frm2.tum.de/en/secondary-sources

Neutron guides

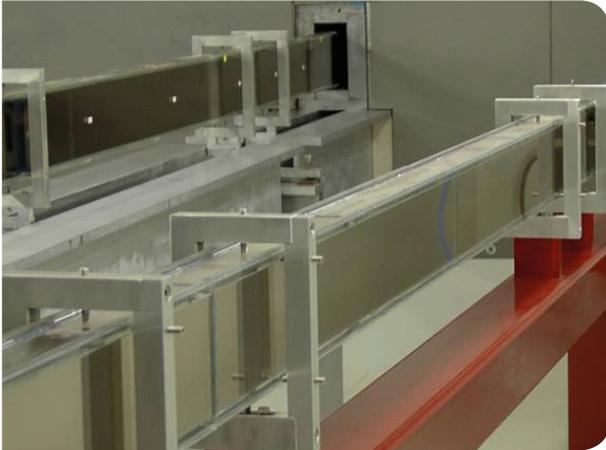

The MLZ makes extensive use of modern neutron guides to transport and distribute the neutrons over large distances in the Experimental Hall (SR-2, -5, -8) as well as in the Neutron Guide Halls. Adapted to the needs of the instruments with respect to wavelength distribution and angular dispersion the guide elements are coated with 58Ni or supermirror coatings with m values up to 3.0; on focussing sections up to $m = 3.6$.

Cold neutron guides

Beam tube SR-1 facing the cold neutron source delivers the neutron beams for the entire Neutron Guide Hall West. The in-pile unit of SR-1 consists of a mirror box with $m = 2.2$ supermirror coating on Al-plates and a dividing section of 2.1 m length, where the beam is divided into the six principal neutron guides. Fig. 1 shows schematically how the neutron guides are further split in order to serve a maximum number of instruments, especially with terminal positions.

Besides SR-1, the cold neutron three axes spectrometer PANDA at beam tube SR-2 has also supermirror inserts in the the primary beam's in-pile section. In the near future, SR-4b will be equipped with a neutron guide for the nuclear and particle physics beam line MEPHISTO in the Neutron Guide Hall East. Until now, almost 500 m of cold neutron guides have already been installed.

Thermal neutron guides

Beam tubes SR-8a and SR-8b are equipped with supermirror guides with coatings up to $m = 3$. At SR-5b a polarising supermirror bender provides the instrument TRISP with polarised neutrons. In total, roughly about 50 m of thermal guides are installed. In the near future, elliptical focussing thermal guides will provide neutron beams for the instruments in the Neutron Guide Hall East.

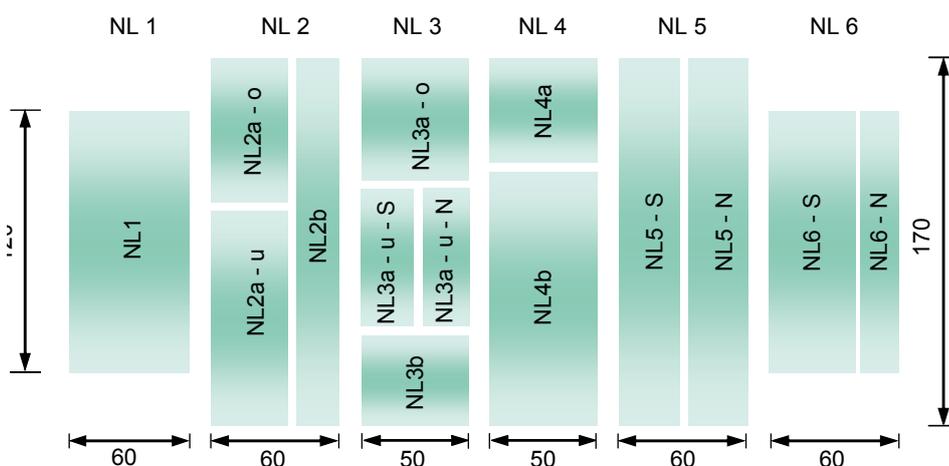


Figure 1: The sectioning of the six principal neutron guides of SR-1. The subsections are indicated by the letters a, b, u, o, S, N.

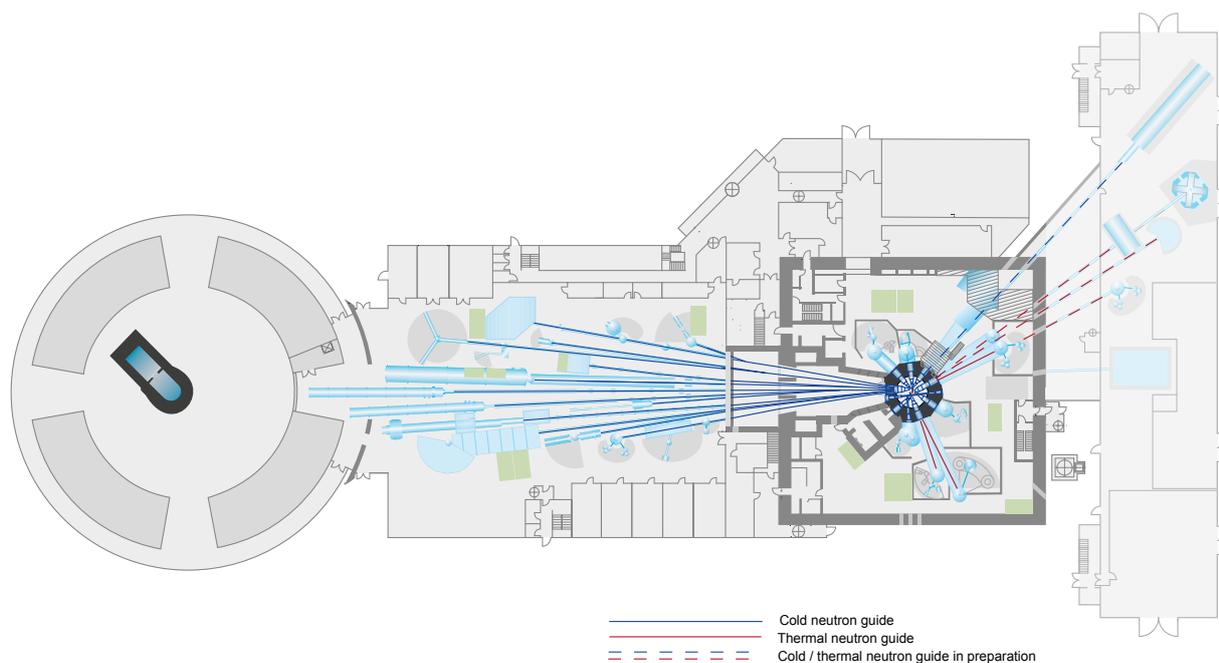


Figure 2: Floorplan of the FRM II with the neutron guide system.

Technical Data

Guide	NL1	NL2a-o	NL2a-u	NL2b	NL3a-o	NL3a-uS	NL3a-uN	NL3b
Length (m)	40	48	60	57	46	30	29	51
Section (mm ²)	60 × 120	44 × 60	44 × 100	12 × 170	50 × 50	10 × 56	38 × 56	50 × 45
Coating up to	m = 2.5	m = 3.0	m = 2.0	m = 2.0	m = 3.0	m = 3.0	m = 3.0	m = 2.0
Curvature (m)	1000	160	2000	400	460	30	460	800
Instrument(s)	BIODIFF NREX	J-NSE	TOFTOF	REFSANS	KWS-2	KWS-3	unused	KWS-1

Guide	NL4a	NL4b	NL5-S	NL5-N	NL6-S	NL6-N
Length (m)	34	52	70	35	54	35
Section (mm ²)	50 × 50	50 × 110	29 × 170	29 × 170	60 × 120	10 × 120
Coating up to	m = 2.0	m = 3.0	m = 2.0	m = 2.0	m = 2.2	m = 2.0
Curvature (m)	480	390	1640	400	1000	84
Instrument(s)	SANS-1	PGAA	RESEDA TREFF	MARIA	MIRA-2 DNS SPHERES	MIRA


Dr. Peter Link

Email: peter.link@frm2.tum.de
Office: +49.(0)89.289.14622





BIODIFF
diffractometer for large unit-cells

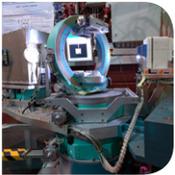

HEIDI
single-crystal diffractometer with hot neutrons




MIRA
dual wavelength band instrument




POLI
polarised hot neutron diffractometer



RESI
thermal neutron single crystal diffractometer

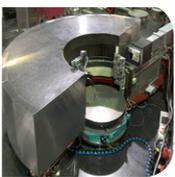

SPODI
high resolution powder diffractometer

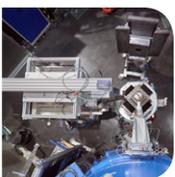

STRESS-SPEC
materials science diffractometer

Structure Research

BIODIFF

diffractometer for large unit cells



The monochromatic single crystal diffractometer BIODIFF is a joint project of the FRM II (TUM) and JCNS (Forschungszentrum Jülich).

BIODIFF is designed to handle crystals with large unit cells and is dedicated to the structure determination of biological macromolecules. In biological macromolecules, like proteins and nucleic acids, hydrogen atoms play an important role. Hydrogen atoms take part in the substrate binding process and are essential for proton transfer reactions during the catalysis in many enzymes. Therefore, the knowledge about the protonation states of amino acid residues in the active centre of proteins is often crucial for the understanding of their reaction mechanisms. However, hydrogen atoms, especially rather flexible ones, are barely detectable in X-ray structure determinations of proteins. On the other hand, hydrogen atoms are clearly visible in neutron crystallography experiments even at moderate resolutions ($d_{\min} < 2.0 \text{ \AA}$).

BIODIFF is the first instrument along the cold neutron guide NL1 and is positioned in a distance of about 32.5 m from the cold source. Using a pyrolytic graphite monochromator PG(002) the diffractometer covers a wavelength range of 2.4 Å to about 5.6 Å. Higher order wavelength contaminations are removed by a neutron velocity selector. The main

detector of the diffractometer consists of a neutron imaging plate system in a cylindrical geometry to cover a large solid angle. A fast LiF/ZnS scintillator CCD camera is foreseen for additional detection abilities. The main advantage of this instrument is the possibility to adapt the wavelength to the size of the sample crystal's unit cell while operating with a clean monochromatic beam that keeps the background level low.

Typical Applications

The main field of application is the neutron structure analysis of proteins, especially the determination of hydrogen atom positions. Typical questions in this field of interest are:


- Enzymatic mechanism (protonation states of amino acids)
- Ligand binding mediated by hydrogen bonds
- Investigation of the hydration shell of proteins
- H/D-exchange pattern as a monitor of structural stability/flexibility

Sample Environment

Besides standard sample environment BIODIFF provides:

- Oxford Cryosystems Cryostream 700 plus with a temperature range of 90 K to 500 K
- Closed cycle cryostat 3.5 - 325 K

- ① Neutron guide NL1
- ② Monochromator
- ③ Selector
- ④ Variable slits
- ⑤ Sample chamber
- ⑥ Cylindrical image plate detector
- ⑦ CCD detector
- ⑧ Beam stop

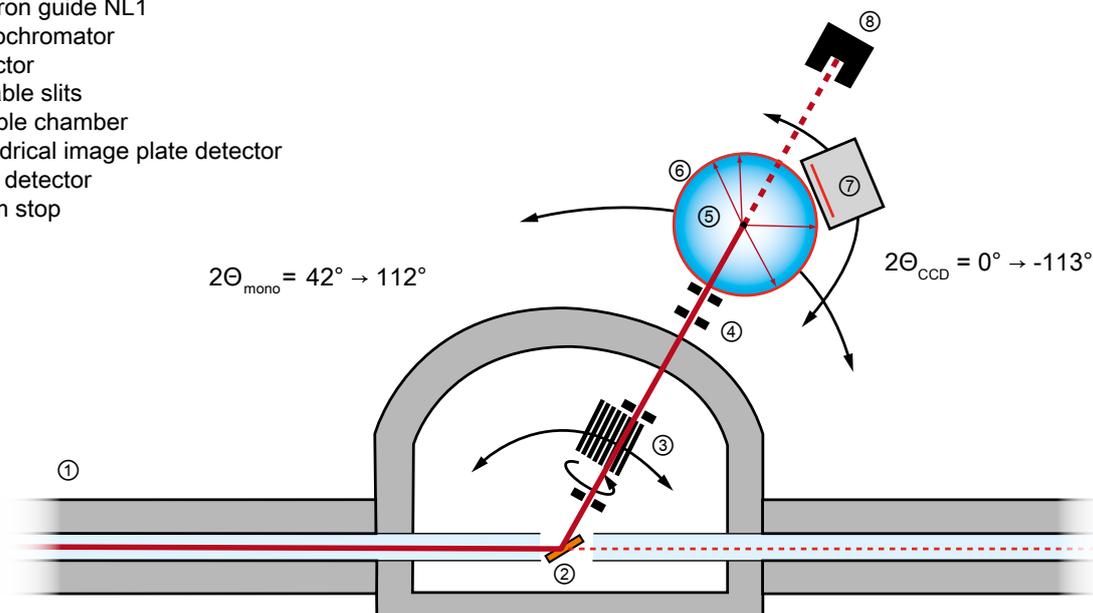

Technical Data

Primary beam

- Neutron guide NL1; supermirror $m = 2$
- Monochromator:
PG(002) mosaicity: $0.4 - 0.5^\circ$
- Higher order filter:
Astrium type velocity selector
transmission 87 % for 2.4 Å
- Wavelength range:
2.4 – 5.6 Å with selector
2.4 – 6.1 Å without selector
- Collimation by adjustable slits down to
 $\varnothing = 1 \text{ mm}$

Beam properties at the sample position

- Wavelength resolution at sample position:
 $\Delta\lambda/\lambda = 2.9\%$ at 2.4 Å
- Beam divergence (no slits)
0.8° FWHM horizontal
0.7° FWHM vertical

Main detector

- Neutron image plate (cylindrical)
- BaFBr:Eu²⁺ mixed with Gd₂O₃
 - Dimensions:

radius	200 mm
angular range	$\pm 152^\circ$ horizontal $\pm 48^\circ$ vertical
 - Pixel size (quadratic) 125, 250, 500 μm
 - Readout time (with erasing):
5 min (for 500 μm pixel size)

Auxiliary detector

- CCD camera with scintillator
- ZnS mixed with ⁶LiF
 - Dimensions:

Active scintillator area (flat)	200 x 200 mm ²
Distance to sample	100 mm
 - 2Θ-angle around sample position $0^\circ - 113^\circ$
 - CCD chip with 2048 x 2048 pixels
 - Pixel size: 13.5 x 13.5 μm²
 - Overall spatial resolution $\approx 300 \times 300 \mu\text{m}^2$
(limited by scintillator thickness)
 - Minimum readout time
 $\approx 1 \text{ sec}$ (full resolution); $< 1 \text{ sec}$ (binning mode)

Dr. Tobias Schrader

Dr. Andreas Ostermann


Email: t.schrader@fz-juelich.de andreas.ostermann@frm2.tum.de
 Office: +49.(0)89.289.10743 +49.(0)89.289.14702
 Instrument: +49.(0)89.289.14565



www.mlz-garching.de/biodiff
 DOI: 10.17815/jlsrf-1-19

HEIDI

single crystal diffractometer on hot source

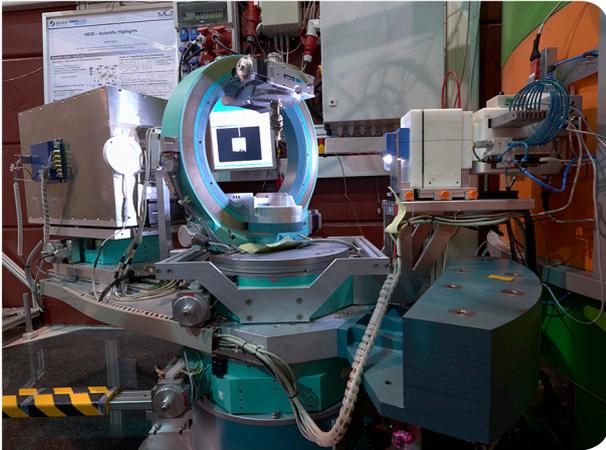

The single crystal diffractometer HEiDi is designed for detailed studies on structural and magnetic properties of single crystals using unpolarised neutrons and Bragg's Law:

$$2 d_{hkl} \sin(\Theta) = \lambda$$

Because of the large variety of short wavelengths and resolutions, HEiDi is suitable for studies on a lot of crystalline compounds – many of them of potential interest for energy or data storage technologies – like:

- HT superconductors (e.g. cuprates, FeAs-pnictides)
- Multiferroics (e.g. manganates) and other complex ferro- and antiferromagnetic compounds (e.g. olivines)
- Ionic conductors (e.g. nickelates)
- Ferroelectrics (e.g. KDP family)
- Mixed crystals (e.g. AsSe compounds)
- Highly absorbing compounds (e.g. with Gd, Sm, Eu, Dy)
- Frustrated magnetic materials (e.g. pyrochlores)

Applications (in general)


- Structure analysis
- Hydrogen bonds
- Static and dynamic disorder
- Harmonic and anharmonic mean square displacements
- Twinning
- Magnetic structure and order
- Structural and magnetic phase transitions
- Incommensurate structures

Applications (in detail)

- Studies of atomic positions and bond distances in compounds with heavy and light elements or elements of similar electron shells
- Temperature dependent studies for determination of phase transitions
- Studies of order – disorder phase transitions, e.g. H bonds by determination of anisotropic mean square displacements using large Q range up to $\sin(\Theta)/\lambda > 1$
- Structure determination of compounds with highly absorbing elements (Gd, Sm, Cd, Dy) with short wavelengths
- Studies on magnetic phase transitions and T dependencies (ferri, ferro and antiferro magnets, multiferroics)
- Studies on HT superconductors (e.g. cuprates, FeAs pnictides)
- Sample characterisation by profile analysis
- Determination of sample orientation, e.g. for preparation of experiments on three axes instruments
- Presentation of fundamentals of crystallography and structure analysis for education

Sample Environment

- Closed cycle cryostat (2 K – RT)
- Mirror furnace (RT – 1500 K)
- Micro furnace (RT – 500 K)
- Uniaxial pressure cell (from PUMA)

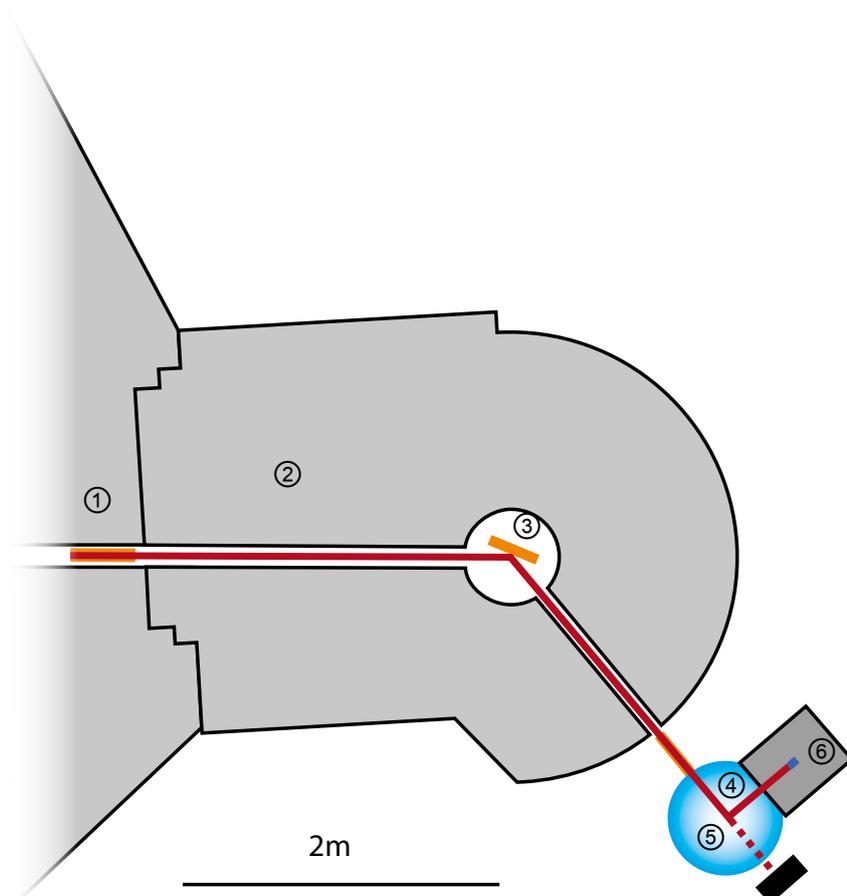

- ① Primary beam collimator
- ② Shielding
- ③ Monochromator
- ④ Sample
- ⑤ Diffractometer
- ⑥ Detector

Technical Data

Beam-tube

- SR-9b (hot source)
- Flux at sample $1.4 \cdot 10^7 \text{ cm}^{-2}\text{s}^{-1}$ ($\lambda \approx 1.17 \text{ \AA}$)
- Gain by hot source $\times 10$ ($\lambda \approx 0.6 \text{ \AA}$)

Wavelength

$2\Theta_M$	Ge(311)	Cu(220)	Ge(422)	Cu(420)
20°	0.503	0.443	0.408	0.280
40°	1.168	0.870	0.793	0.552
50°	1.443	1.079	0.993	0.680

Q-range

$2\Theta_M$	Ge(311)	Cu(220)	Ge(422)	Cu(420)
20°	1.46	1.95	2.12	3.09
40°	0.74	0.99	1.09	1.57
50°	0.60	0.80	0.87	1.27

Optical components

- Single detector optimised for small wavelengths (sensitivity > 90% at 0.3 \AA)
- Analyzer PG(002); optional for studies of purely elastic scattering and background suppression
- Neutron filters for suppression of $\lambda/2$ - or $\lambda/3$ -contamination of the monochromatised beam

Dr. Martin Meven

Dr. Andrew Sazonov

Email: martin.meven@frm2.tum.de andrew.sazonov@frm2.tum.de
 Office: +49.(0)89.289.14727 +49.(0)89.289.11764
 Instrument: +49.(0)89.289.14828



www.mlz-garching.de/heidi
 DOI: 10.17815/jlsrf-1-20

MIRA

dual wavelength band instrument




MIRA is a multipurpose instrument. With its two beam ports, namely MIRA-1 and MIRA-2, it provides neutrons over a wide range of wavelengths $3.5 \text{ \AA} < \lambda < 20 \text{ \AA}$. The instrument can easily be moved from one port to the other without changing the sample environment. The large variety of different options offered can also be combined in most cases. This allows for a fast realisation of experiments in a very flexible way using a number of available options:


- Cold diffraction
- Cold three axes spectroscopy for extreme environments
- Small angle neutron scattering (SANS)
- MIEZE
- Reflectometry
- 3D-Polarimetry

All options can be measured in polarised mode. Using the finger detector, the instrument has a very low background of less than 0.1 cps. For MIRA-2 a q-range up to 2.5 \AA^{-1} with an q-resolution of 0.01 \AA^{-1} can be reached. Vertical and horizontal B-fields up to 2.2 T and vertical B-fields up to 7.5 T are available. Temperatures from 50 mK to 1500 K can be applied using the standard sample environment.

Typical Applications

- Dynamics of magnetic excitations
- Determination of magnetic structures, especially large scale structures, i.e. helical spin density waves or magnetic lattices
- Quasi-elastic measurements in magnetic fields with high resolution
- Determination of structures and dynamics in extreme environments, like pressure
- Determination of layer thickness of films, for instance in polymer physics
- Reflectometry from magnetic multilayers
- Polarisation analysis

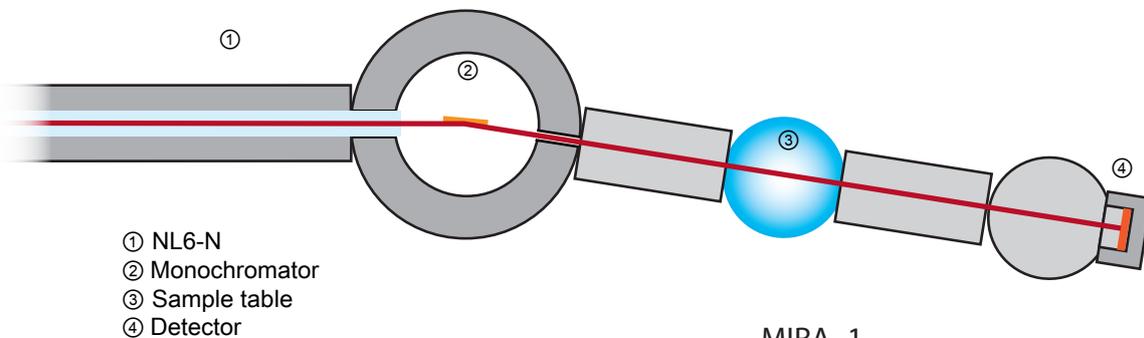

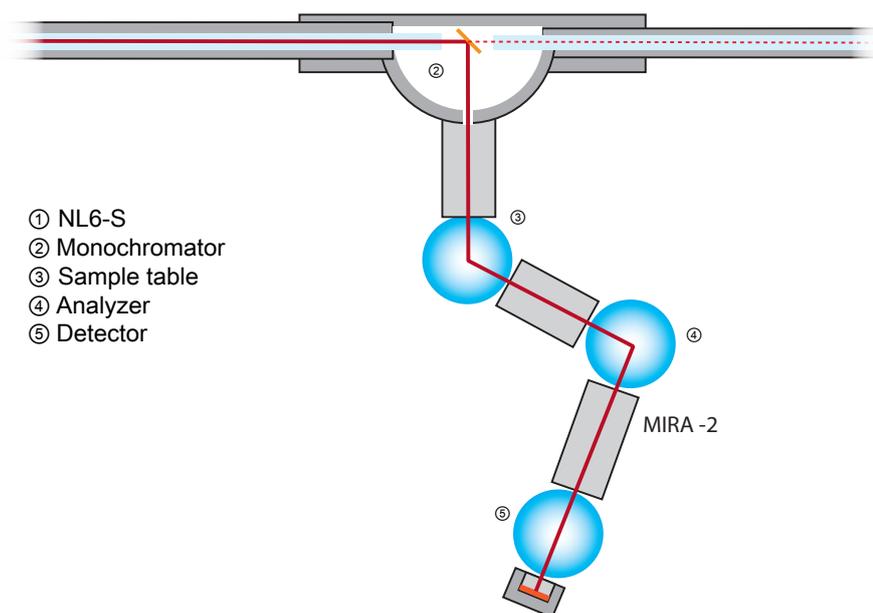

Technical Data

MIRA-1

Primary beam

- Neutron guide: NL6-N
- Dimensions: 10 x 120 mm² (width x height)
- Curvature: 84 m
- Coating: sides m = 2, top/bottom m = 2

Monochromator

- Intercalated HPGO $\Delta\lambda/\lambda = 2\%$
- Multilayer $\Delta\lambda/\lambda \approx 3\%$ (5% polarised)
- $6 \text{ \AA} < \lambda < 20 \text{ \AA}$

Max. differential flux at sample

- $5 \cdot 10^5$ neutrons s⁻¹ cm² at 10 Å
- $2 \cdot 10^5$ neutrons s⁻¹ cm² polarised

Analyzer

- 2 cavities
- 2 bender
- ³He-spin filter

Detector

- 20 x 20 cm² 2-D PSD with 1 x 2 mm² resolution
- 1 inch ³He finger detector
- 20 x 20 cm² 2-D PSD, time resolution < 1 ps

MIRA-2

Primary beam

- Neutron guide: NL6-S
- Dimensions: 60 x 120 mm² (width x height)
- Coating: sides m = 2.0, top/bottom m = 2

Monochromator

- Horizontal focussing HOPG $\Delta\lambda/\lambda \approx 2\%$
- $3.5 \text{ \AA} < \lambda < 6 \text{ \AA}$

Max. differential flux at sample

- $1 \cdot 10^7$ neutrons s⁻¹ cm² at 4.7 Å (2015)
- $1 \cdot 10^6$ neutrons s⁻¹ cm² polarised

Analyzer

- 2 cavities
- S-bender, transmission polariser
- ³He-spin filter

Detector

- 20 x 20 cm² 2-D PSD with 1 x 2 mm² resolution
- 1 inch ³He finger detector
- 20 x 20 cm² 2-D PSD, time resolution < 1 ps
- with low background < 0.1 cps

Dr. Robert Georgii

Dr. Klaus Seemann

Email: robert.georgii@frm2.tum.de klaus.seemann@frm2.tum.de
Office: +49.(0)89.289.14986 +49.(0)89.289.14668
Instrument: +49.(0)89.289.14877



www.mlz-garching.de/mira
DOI: 10.17815/jlsrf-1-21



POLI is a versatile two axes single crystal diffractometer mostly dedicated to the investigation of magnetic structures in single crystals using neutron spin polarisation. The using of the separate non-polarising variably focussing monochromators and optimised polarisers allows the structural investigations using polarised and non-polarised neutrons with short wave length and high resolution. Two standard options are currently implemented on POLI:

- Zero-field spherical neutron polarimetry (SNP) using third generation polarimeter CRYOPAD;
- Non-polarised diffraction under special conditions e.g. very low temperatures (< 3 K), magnetic and electric fields, high pressures, high temperatures, etc. using dedicated sample environments and out-of-plane lifting counter.

A third measurement method named flipping-ratio is under development and planned to be operational using a new 8 T magnet in 2016.

SNP allows distinguishing between the depolarisation and polarisation rotation occurring during the scattering in the sample. X, Y, Z components of the scattered polarisation are measured for each orientation of the incoming polarisation and hence a polarisation matrix of 9 elements for an individual Bragg reflection can be found. Determining the relationship between the directions of incident and scattered polarisations gives access to the 16 independent correlation functions involved in the most

general nuclear and magnetic scattering processes. Generally, this leads to the direction's determination of the magnetic interaction vector of magnetic structures. For those structures coinciding nuclear and magnetic reflections in reciprocal space, SNP leads to the amplitude's determination of the magnetic interaction vectors, and hence to the magnetisation distribution. SNP could also be employed to study the magnetic domain distribution.


Actually, polarised ^3He spin filters are used on POLI in order to produce and to analyse neutron polarisation. This method offers best performance in the polarising of hot neutrons. Tuning the pressure of ^3He in the filter cell it is possible to optimise the polariser's performance depending on wave length or experimental needs. The automatic correction of the time-dependent polarisation is applied. Because of the ^3He filters' sensitivity about non-homogeneous magnetic fields, a dedicated solid state supermirror bender will be used in the future for experiments with high magnetic fields.

Typical Applications

- Complex commensurate and incommensurate magnetic structures studied in ground state (zero-field) – very useful especially for superconductors.
- Studies of magnetic or magneto-electric domains using SNP on samples cooled in zero-field as well as in external magnetic field (up to 7.5 T).
- The combination of magnetic and electric fields applied on the sample – important for the studies on multiferroic materials.
- Determination of magnetic form factors (sometimes also in antiferromagnets).

Sample Environment

- Standard closed-cycle cryostat (4 K – 425 K)
- Variox LHe cryostat (1.5 – 300 K)
- Lower temperature inserts on request
- 2.2 T HTS magnet / 7.5 T magnet on request
- 8 T dedicated magnet planned from 2016
- Electric field up to 10 kV

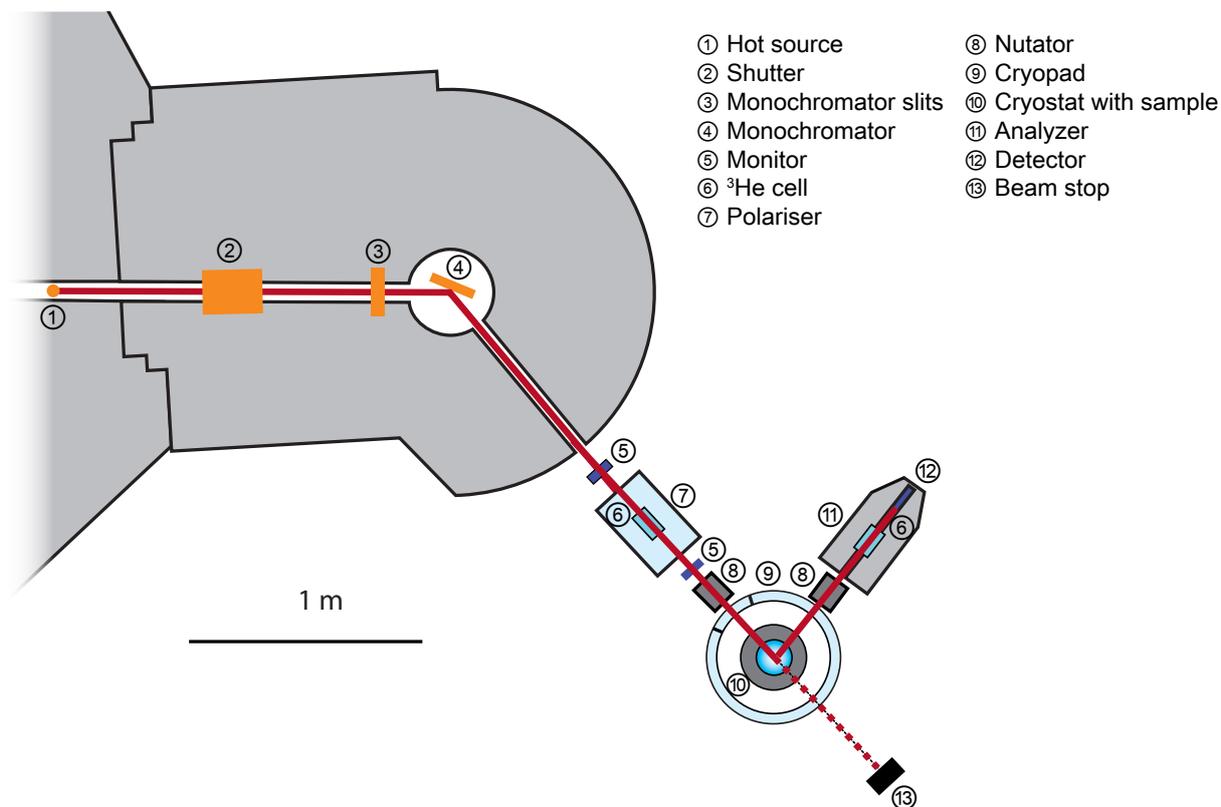

- | | |
|-----------------------|------------------------|
| ① Hot source | ⑧ Nutator |
| ② Shutter | ⑨ Cryopad |
| ③ Monochromator slits | ⑩ Cryostat with sample |
| ④ Monochromator | ⑪ Analyzer |
| ⑤ Monitor | ⑫ Detector |
| ⑥ ^3He cell | ⑬ Beam stop |
| ⑦ Polariser | |

Technical Data

Primary beam: SR-9a on hot source

- Focussing monochromators

Crystal	λ Å	flux n/s/cm ²	λ Å	flux n/s/cm ²
Cu (220)	0.55	$4 \cdot 10^6$	0.9	$2.4 \cdot 10^7$
Si (311)	0.7	$7 \cdot 10^6$	1.15	$2.8 \cdot 10^7$

Neutron polarisation: ^3He spin filter

- ^3He polarisation: 0.75 – 0.7
- Neutron polarisation: 0.93 – 0.8
- Cell holding time: 2 days
- Cells replacement time: 5 – 10 min

Diffractometer angles

with Cryopad	without Cryopad
$-10^\circ < 2\theta < 120^\circ$	$-30^\circ < 2\theta < 130^\circ$
$-180^\circ < \omega < 180^\circ$	$-180^\circ < \omega < 180^\circ$
$-4^\circ < \xi_1 < 4^\circ$	$-5^\circ < \xi_1 < 5^\circ$
$-4^\circ < \xi_2 < 4^\circ$	$-5^\circ < \xi_2 < 5^\circ$
$\nu = 0$	$-4.2^\circ < \nu < 30^\circ$

Cryopad (zero-field polarimeter)

- Precision in polarisation control: better 1°
- Low background / low absorption
- LHe autonomy: 10 days
- LN₂ refill (automatic): daily
- Sample space for closed cycle cryostat or orange type cryostat

Dr. Vladimir Hutanu

Email: vladimir.hutanu@frm2.tum.de
 Office: +49.(0)89.289.12153
 Instrument: +49.(0)89.289.14480



www.mlz-garching.de/poli
 DOI: 10.17815/jlsrf-1-22



The diffractometer RESI is designed for using a maximum of thermal neutron intensity, allowing optimum measurements of weak diffraction phenomena in a large portion of the reciprocal space on single crystalline samples.

Typical Applications

Structure analysis with thermal neutrons ($\lambda = 0.8 \text{ \AA}$ to 2 \AA) is complementary to structure analysis with X-rays. The measurement possibilities provided by this instrument are crucial for many scientific questions:

- **Structure analysis, bonding theory, electron densities:** Due to the interaction with atomic cores and the diffraction angle independence of the atomic form factor, it is possible to measure Bragg scattering up to high diffraction angles.
- **Real crystals** and compounds of interest for material science are often not perfectly ordered. The elucidation of these real structures requires the analysis of the corresponding diffuse scattering. The diffuse scattering - off the Bragg reflections - is normally differentially weak and distributed continually (anisotropic) in the reciprocal space.
- **Partially crystalline** compounds, like **fibre structures**, show a specific scattering, which is highly anisotropic and continuously distributed in the reciprocal space. Therefore, diffractometers with area detectors like RESI are best suited for this kind of problems.

- A new class of **aperiodic crystals** (“quasi crystals”) show dense, but discrete reflex patterns, where more than 90% of the reflexes are very weak. Additionally, due to the fact that quasi crystals often contain two or more transition metals (which are almost isoelectronic), neutrons offer much higher contrast than X-ray methods.
- **Structural phase transitions** can be accompanied by continuous reflection shifting.
- **Modulated structures** show satellite reflections at “incommensurable” positions. Both areas require analysis of large portions of the reciprocal space.
- **Twinned crystals** and **multi-domain/multi-phase crystals** are often difficult to measure on single-counter instruments. The area detector at RESI allows for easy detection and in many cases separation of reflections in such systems.

The advantages of the high-resolution area detector can be utilised best, if the reciprocal space is not too empty. That means, that RESI is optimal for cells of ca. 1000 \AA^3 to ca. 20000 \AA^3 . Typical crystal sizes range from 5 mm^3 to 25 mm^3 .


Sample Environment

Dedicated sample environment of RESI:

- Oxford Cryosystems Cryostream 700
temperature range 100 K - 400 K
consumption $\sim 20 \text{ l L-N}_2/\text{d}$
- Oxford Instruments Helijet
temperature range 15 K - 100 K
consumption $\sim 2 \text{ l L-He / h}$
sample size $1 \times 1 \times 1 \text{ mm}^3$ max.

Standard sample environment usable with RESI

- Closed-cycle cryostat CC, 2.5 K – 300 K
- Closed-cycle cryostat CCR, 3 K – 100 K
using ^3He insert, 500 mK – 4 K
using $^3\text{He}/^4\text{He}$ dilution, 50 mK – 1 K
- Vacuum furnace, 340 K – 2100 K
- Mirror furnace, RT – 1250 K

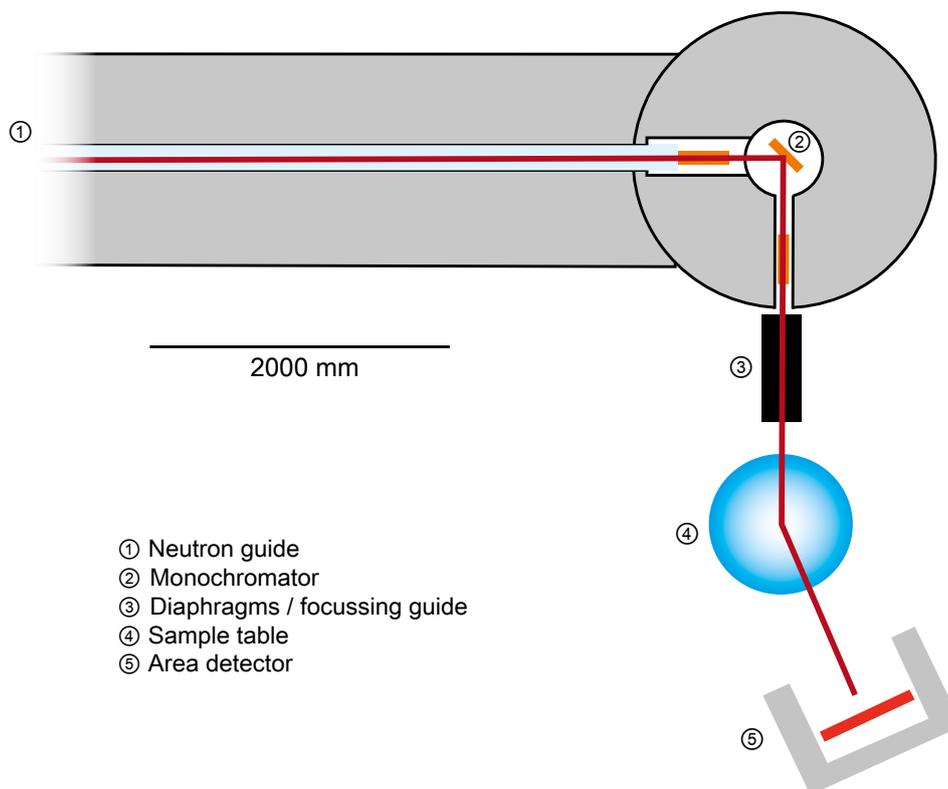

Technical Data

Primary beam

- Beam tube SR-8b
- Neutron guide
Length: 12 m, focussing vertical / horizontal
section: 70 × 40 mm → 60 × 30 mm
- Coatings: m = 3 top/bottom; m = 1 side

Monochromators

- Vertically focussing lamella type, fixed take-off 90°
- Cu-422, 20' mosaic, 1 Å : $2 \cdot 10^6 \text{ n cm}^{-2}\text{s}^{-1}$
 - Ge-511, 25' mosaic (deformed wafer stack)
1.5 Å : $6 \cdot 10^6 \text{ n cm}^{-2}\text{s}^{-1}$

Secondary neutron guide

- Vertically focussing elliptical guide-in-guide
- Length: 1 m
 - Focus 400 mm after guide exit
 - Coating: m = 5

Available goniometers

- Kappa-Goniometer:
Bruker-Nonius Mach3
carrying capacity: max 100 g
- Eulerian cradle Huber 420:
higher carrying capacity, e.g. for closed-cycle
cryostat
- Huber 2-circle goniometer with tilting head
highest carrying capacity, e.g. for CCR with ^3He
insert

Available detectors

- MAR345 image plate detector:
345 mm diameter, N-sensitive image plate
- Single counter ^3He with optional analyzer
for pure elastic scattering

Dr. Bjørn Pedersen


Email: bjoern.pedersen@frm2.tum.de
 Office: +49.(0)89.289.14707
 Instrument: +49.(0)89.289.14827



www.mlz-garching.de/resi
 DOI: 10.17815/jlsrf-1-23

SPODI

high resolution powder diffractometer

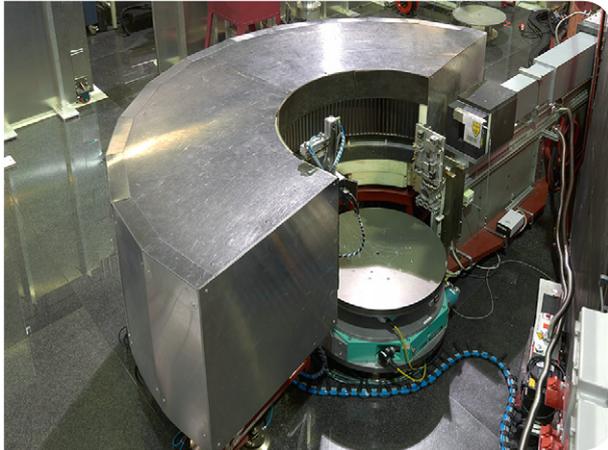

The high resolution powder diffractometer SPODI is designed for structure solution and Rietveld refinement of structural parameters on crystalline powders. The instrument is characterised by a very high monochromator take-off angle of 155° (standard configuration). Optionally, a take-off angle of 135° is available.

The detector array consists of 80 ^3He position sensitive detector tubes (300 mm active height) with fixed Soller collimators of $10'$ horizontal divergence. The multidetector of SPODI spans an angular range of $2\theta = 160^\circ$. Each detector covers 2° corresponding to $160^\circ / 80$ detectors. Therefore the data collection is performed via stepwise positioning of the detector array to obtain a diffraction pattern of the desired step width (typically $2^\circ / 40$ steps resulting in $\Delta(2\theta) = 0.05^\circ$).

The two-dimensional raw data are evaluated to provide diffraction patterns corresponding to different detector heights ranging from 10 mm to 300 mm and variable detector height, accounting for vertical beam divergence effects. Thus, asymmetric broadenings at quite low and high scattering angles are overcome, while the full detector height in the medium 2θ regime can be used [1].

Various sample environmental devices enable the characterisation of materials under special conditions: A rotatable tensile rig allows in-situ studies under tensile stress, compression stress or torsion

while the load axis can be oriented with respect to the scattering plane. A potentiostat for charging/discharging of Lithium ion batteries is available as well as a device to apply high electric fields on ferroelectrics.

Typical Applications

- Determination of complex crystal and magnetic structures
- Structural evolutions and phase transformations under various environmental conditions.
- Static and thermal disorder phenomena

Research Areas

- Ionic conductors
- Materials for lithium ion batteries
- Ferroelectrics, multiferroics
- Hydrogen storage materials
- Shape memory alloys
- Superalloys
- Correlated electron systems
- Superconductors
- Minerals

Sample Environment


Standard sample environment of FRM II

- Closed cycle cryostat 3 – 550 K (with ^3He insert: $T_{\min} = 500$ mK)
- Vacuum high temperature furnace
 $T_{\max} = 1900^\circ\text{C}$
- Cryomagnet
 B_{\max} at SPODI: 5 T
- Sample changer (six samples, ambient temperature)

Special sample environment

- Rotatable tensile rig
 $F_{\max} = 50$ kN, $M_{\max} = 100$ Nm
- Device for electric fields
 $V_{\max} = 35$ kV
- Potentiostat for electrochemical treatment of materials VMP3 and SP240

[1] Hoelzel, M. et al., Nucl. Instr. A 667, 32-37 (2012).

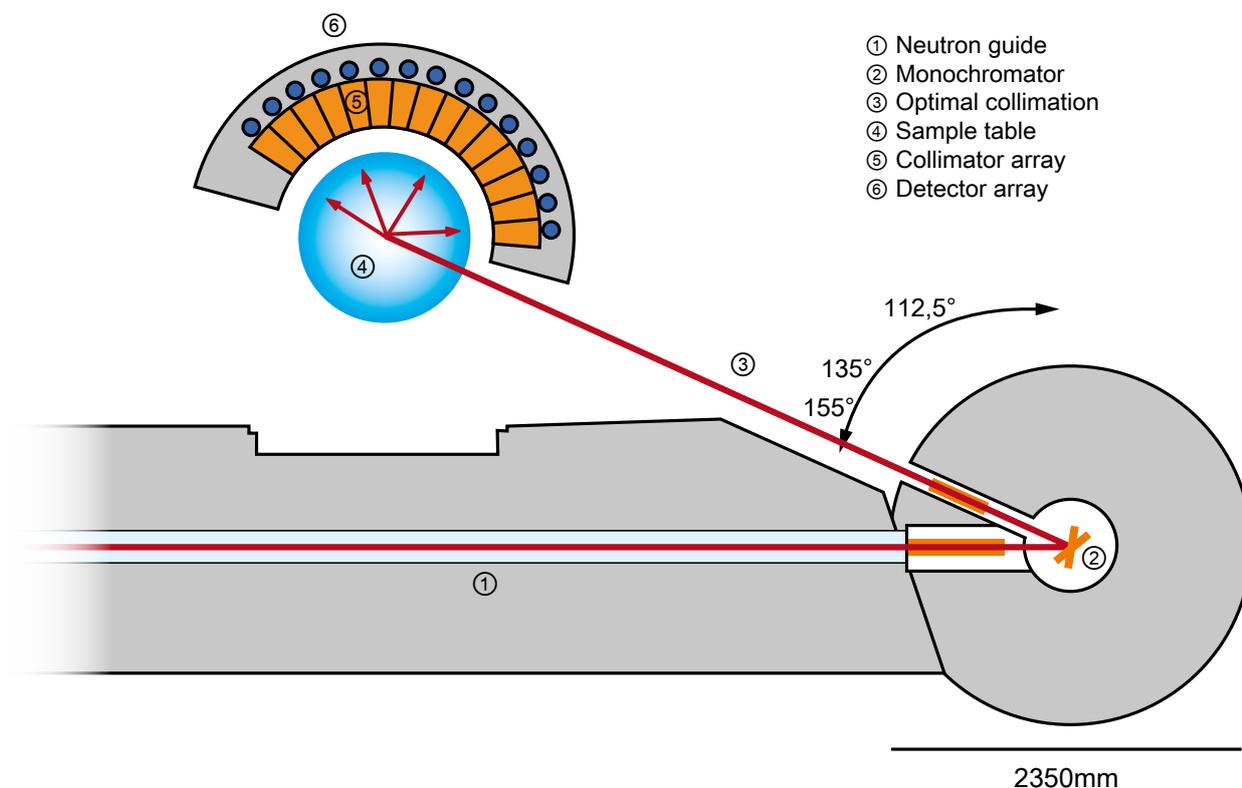

Technical Data

Monochromator

- Ge(551) wafer stack crystals
- standard configuration: take-off angle 155°
 - Ge(551): 1.548 Å
 - Ge(331): 2.436 Å
 - Ge(711): 1.111 Å

Collimation

- $\alpha_1 \approx 20'$ (neutron guide)
- $\alpha_2 = 5', 10', 20', 25'$ nat. (for 155°)
- $\alpha_2 = 10', 20', 40'$ nat. (for 135°)
- $\alpha_3 = 10'$

Detector array

- 80 position-sensitive ^3He tubes,
 - angular range $2\theta = 160^\circ$,
 - effective height: 300 mm

Dr. Markus Hölzel

Email: markus.hoelzel@frm2.tum.de
 Office: +49.(0)89.289.14314
 Instrument: +49.(0)89.289.14826

Dr. Anatoliy Senyshyn

anatoliy.senyshyn@frm2.tum.de
 +49.(0)89.289.14316

Dr. Oleksandr Dolotko


odolotko@frm2.tum.de
 +49.(0)89.289.14373



www.mlz-garching.de/spodi
 DOI: 10.17815/jlsrf-1-24

STRESS-SPEC

materials science diffractometer

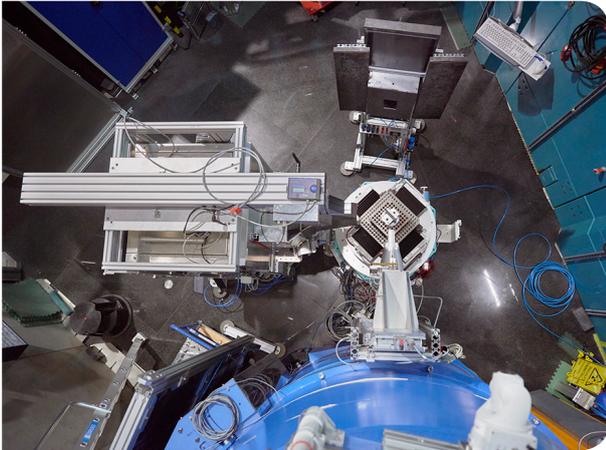

In response to the development of new materials and the application of materials and components in new technologies the direct measurement, calculation and evaluation of textures and residual stresses has gained worldwide significance in recent years.

STRESS-SPEC, the materials science diffractometer is located at the thermal beam port SR-3 of the FRM II and can easily be configured either for texture or stress analysis.

The set-up utilises three different monochromators: Ge (511), bent silicon Si (400) and pyrolytic graphite PG(002). This selection of monochromators and the possibility to vary automatically the take-off angles from $2\theta_M = 35^\circ$ to 110° allows to find a good compromise between resolution and intensity for each measuring problem.

The gauge volume defining optical system of primary and secondary slits is designed with regard to reproducibility of geometrical alignment and sturdiness. Both slit systems are linked to the sample table and the detector in such a way that the center of the beam remains the same under all conditions. Samples can be aligned using theodolites and a camera system. In addition the possibility to scan surfaces of components offline using a CMM laser scanner is available at STRESS-SPEC.

Typical Applications

Residual stress analysis [1]

- Industrial components
- Welds
- Superalloys
- Strain mapping
- Surface measurements from 150 μm possible [2]

Texture determination [3]

- Global textures
- Local textures
- Strain pole figures
- FWHM pole figures

Structural applications

- Phase transformation dynamics
- Spatially resolved phase analysis (e.g. batteries)

Sample Environment


- XYZ-table
capacity 300 kg, Travel $xy = \pm 120$ mm,
 $z = 300$ mm, accuracy ~ 10 μm
- Load frame
 ± 50 kN, heatable to 1000°C
- Full circle Eulerian cradle (max. load 5 kg)
- $\frac{1}{4}$ circle Eulerian cradle for heavy samples
- Standard sample environment (e.g. furnace, cryostat)

A positioning system consisting of a Stäubli-6-axes robotic arm for texture and strain measurements (payload up to 30 kg) can be mounted instead of the standard sample table. It offers more flexibility than an Eulerian cradle and can be also used as automatic sample changer for texture measurements.

[1] Hofmann, M. et al., *Physica B*, 385-386, 1035-1037 (2006).

[2] Saroun, J. et al., *J. Appl. Cryst.*, 46, 628-638 (2013).

[3] Brokmeier, H.-G. et al., *Nucl. Instr. Meth. A*, 642, 1, 87-92 (2011).

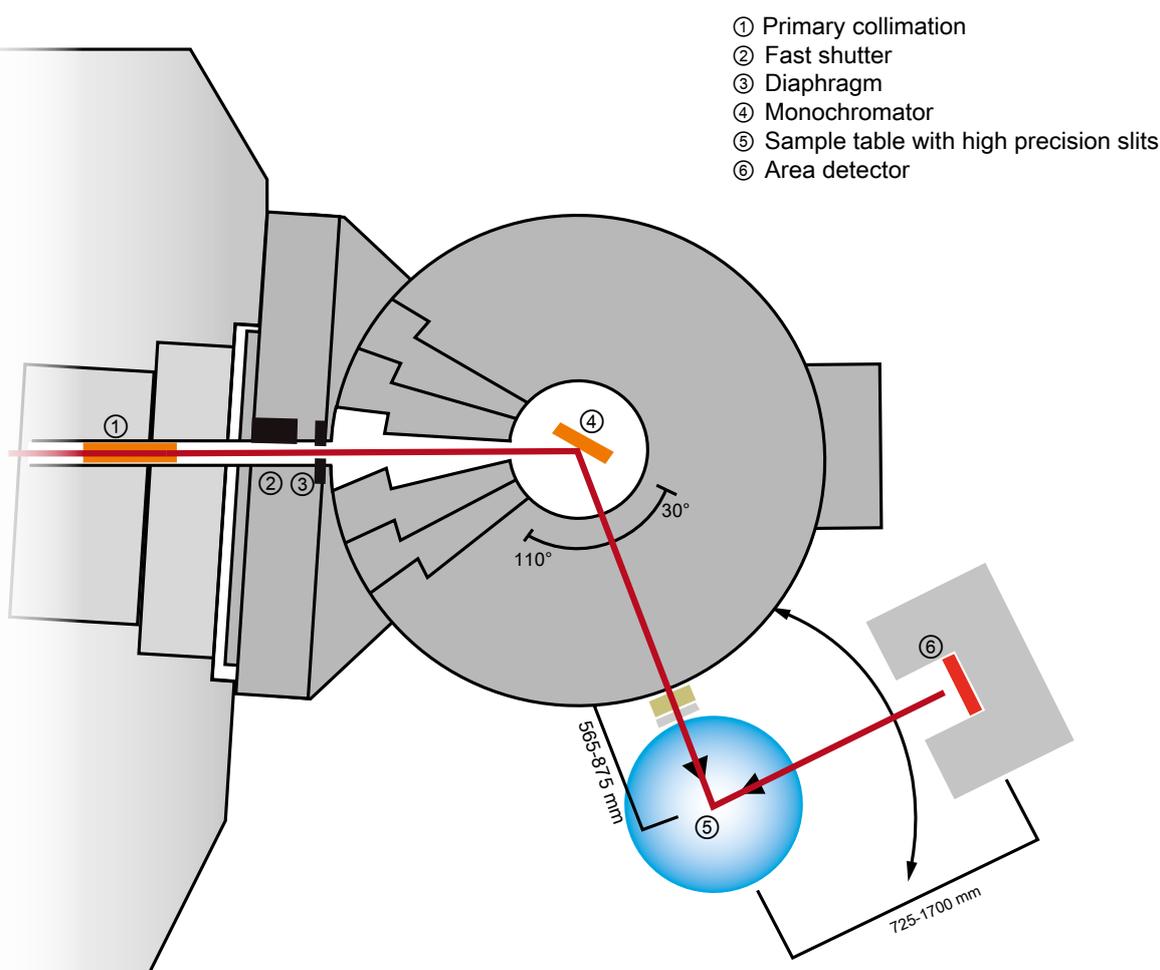

- ① Primary collimation
- ② Fast shutter
- ③ Diaphragm
- ④ Monochromator
- ⑤ Sample table with high precision slits
- ⑥ Area detector

Technical Data

Neutron beam

- SR-3 thermal neutrons
- Collimators ('in-pile') 15', 25', open

Monochromators

- Ge(511), Si(400), PG(002)
- $2\theta_M$ 35° – 110° continuous
- Wavelength 1 Å – 2.4 Å ; ($2.5 \text{ \AA}^{-1} < Q < 10.5 \text{ \AA}^{-1}$)

Possible slit size - Residual Stress

- Primary slit: automatic continuously variable up to 7 x 17 mm² (W x H)
- Secondary slit: continuously variable up to 15 mm
- Radial collimators (FWHM = 1 mm, 2 mm, 5 mm, 10 mm)

Possible slit size – Textures

- Primary slit: max. 30 x 40 mm² (W x H)
- Secondary slit: continuously variable up to 15 mm or open

Detector

- ³He-PSD, 25 x 25 cm²; 256 × 256 pixel

Dr. Michael Hofmann

Dr. Joana Rebelo-Kornmeier

Dr. Weimin Gan

Email: michael.hofmann@frm2.tum.de
 Office: +49.(0)89.289.14744
 Instrument: +49.(0)89.289.14814

joana.kornmeier@frm2.tum.de
 +49.(0)89.289.14710

weimin.gan@hzg.de
 +49.(0)89.289.10766



www.mlz-garching.de/stress-spec
 DOI: 10.17815/jlsrf-1-25



KWS-1
small angle scattering diffractometer



KWS-2
small angle scattering diffractometer



KWS-3
very small angle scattering diffractometer with focusing mirror



MARIA
magnetic reflectometer with high incident angle



NREX
neutron reflectometer with X-ray option



REFSANS
horizontal time-of-flight reflectometer with GISANS option



SANS-1
small angle neutron scattering

Large Scale Structures

KWS-1

small angle scattering diffractometer



The KWS-1 is dedicated to high resolution measurements [1] due to its 10% wavelength selector. This property is interesting for highly ordered or highly monodisperse samples. With the foreseen chopper the wavelength uncertainty can be reduced further to ca. 1%. The scientific background of KWS-1 is placed in magnetic thin films. Magnetic samples will be studied with the full polarisation analysis including incident beam polarisation and polarisation analysis of the scattered neutrons. In front of the collimation, a 3-cavity polariser with V-shaped mirrors is placed. The full bandwidth of 4.5 to 20 Å will be covered with min. 90% (95% typical) polarisation. A radio frequency spin flipper allows for changing the polarisation. The polarisation analysis will be realised with ^3He -cells which will be optimised for the used wavelength and scattering angle. Vertical magnets will be provided to render the magnetic field at the sample position. Thin films can be well studied in the grazing incidence geometry – the method is called **g**razing incidence **s**mall **a**nge **n**eutron **s**cattering (GISANS). A newly installed hexapod will allow for positioning the sample with 0.01 mm and 0.01° precision.

Classical soft-matter systems will be investigated on KWS-1 if the resolution is needed. Biological samples can be handled due to the detector distance of ca. 1 m, which will allow for maximal scattering angles of $Q = 0.5 \text{ \AA}^{-1}$.

The MgF_2 lenses are used for the high flux mode with large sample areas, while the resolution stays in the classical SANS range. These enhanced intensities allow for real time measurements in the 1/10 second region (typical 1 s).

The chopper in parallel allows for studying faster dynamics in the ms range. The so-called TISANE mode interlocks the chopper frequency with the excitation field frequency and with the detection binning. The precise consideration of the flight times allows for higher precision compared to classical stroboscopic illuminations.


Typical Applications

- Grain boundaries
- Alloys
- Magnetic structures
- Flow lines
- Soft matter and biology (as for KWS-2)
- Complex fluids near surfaces
- Polymer films
- Magnetic films
- Nanostructured films

Sample Environment

- Rheometer shear sandwich
- Rheowis-fluid rheometer (max. shear rate 10000 s^{-1})
- Anton-Paar fluid rheometer
- Stopped flow cell
- Sample holders: 9 horizontal x 3 vertical (temperature controlled) for standard Hellma cells 404-QX and 110-QX
- Oil & water thermostats (range -40 – +250) electric thermostat (RT – 200°C)
- 8-positions thermostated (Peltier) sample holder (-40°C – 150°C)
- Magnet (horizontal, vertical)
- Cryostat with sapphire windows
- High temperature furnace
- Pressure cells (500 bar, 2000 bar, 5000 bar)

[1] Feoktystov, A. et al., J. Appl. Cryst. 48 (2015), 61-70.

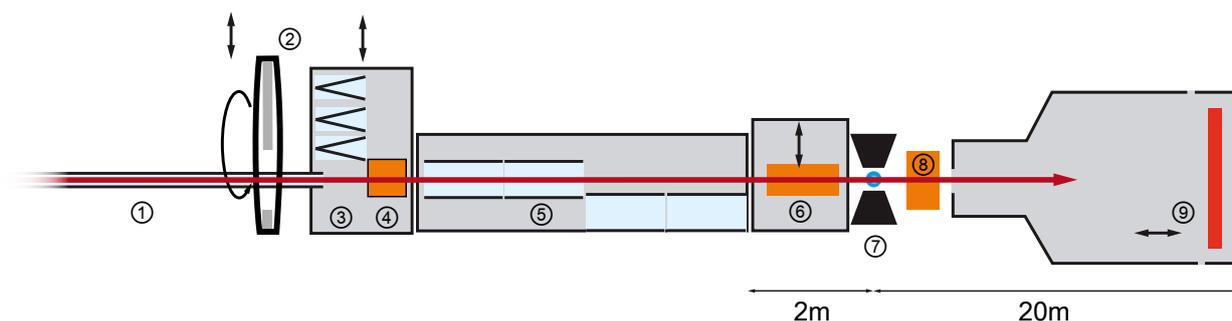

- | | |
|---|--|
| ① Neutron guide NL3 | ⑥ MgF_2 focussing lenses |
| ② High-speed chopper
$\Delta\lambda/\lambda=1\%$ | ⑦ Sample position with magnet |
| ③ Changeable polarisers | ⑧ ^3He spin filter
with reversible polarisation
(to be implemented) |
| ④ Spin flipper | ⑨ Anger-type scintillation detector |
| ⑤ Neutron guide sections 18 x 1 m | |

Technical Data

Overall performance

- $Q = 0.0007 - 0.5 \text{ \AA}^{-1}$
- Maximal flux: $1.5 \cdot 10^8 \text{ n cm}^{-2} \text{ s}^{-1}$
- Typical flux: $8 \cdot 10^6 \text{ n cm}^{-2} \text{ s}^{-1}$ (collimation 8 m, aperture $30 \times 30 \text{ mm}^2$, $\lambda = 7 \text{ \AA}$)

Velocity selector

- Dornier, FWHM 10%, $\lambda = 4.5 \text{ \AA} - 12 \text{ \AA}$, 20 \AA

Chopper

- For TOF-wavelength analysis, FWHM 1%

Polariser

- Cavity with V-shaped supermirror, all wavelengths
- Polarisation better 90%, typical 95%

Spin-flipper

- Radio-Frequency spin flip probability better than 99.8%

Active apertures

- 2 m, 4 m, 8 m, 14 m, 20 m

Aperture sizes

- Rectangular $1 \times 1 \text{ mm}^2 - 50 \times 50 \text{ mm}^2$

Sample aperture

- Rectangular $1 \times 1 \text{ mm}^2 - 50 \times 50 \text{ mm}^2$

Neutron lenses

- MgF_2 , diameter 50 mm, curvature 20 mm
- Packs with 4, 6, 16 lenses

Sample stage

- Hexapod, resolution better than 0.01° , 0.01 mm

Detector 1

- Detection range: continuous 1.5 m – 20 m
- ^6Li -Scintillator 1 mm thickness + photomultiplier
- Efficiency better than 95%
- Spatial resolution $5.3 \times 5.3 \text{ mm}^2$, 128 x 128 channels
- Max. countrate 0.6 MHz ($\tau_{\text{dead}} = 0.64 \text{ \mu s}$)

Dr. Henrich Frielinghaus

h.frielinghaus@fz-juelich.de
+49.(0)89.289.10706
+49.(0)89.289.14324

Dr. Artem Feoktystov

a.feoktystov@fz-juelich.de
+49.(0)89.289.10746

Dr. Ida Berts

i.berts@fz-juelich.de
+49.(0)89.289.10758

Dr. Gaetano Mangiapia

g.mangiapia@fz-juelich.de
+49.(0)89.289.54810



www.mlz-garching.de/kws-1
DOI: 10.17815/jlsrf-1-26

KWS-2

small angle scattering diffractometer



KWS-2 [1] represents a classical pinhole SANS instrument where, combining the pinhole mode using different neutron wavelengths and detection distances with the focusing mode using MgF_2 lenses, a wide Q-range between 1×10^{-4} and 0.5 \AA^{-1} can be explored. It is dedicated to high intensity/ wide-Q investigation of mesoscopic structures and structural changes due to rapid kinetic processes in soft condensed matter, chemistry, and biology. The high neutron flux, comparable with that of the world leading SANS instruments, which is supplied by the neutron delivery system (cold source, selector, guides) [2, 3], and the possibility to use large sample area using focussing lenses, enable high intensity and time-resolved studies. On demand, the instrument resolution can be tuned using the double-disc chopper with adjustable opening slit, which allows the variation of the wavelength spread between 2 and 20%. This offers a high flexibility in optimising the instrument performance towards improved characterisation of structural details and accurate beam characteristics (avoid the gravity and chromatic effects while using the lenses).

Typical Applications

- Colloids, nanocomposites, polymer gels, networks
- Polymer blends, diblock copolymers
- Microemulsions, complex fluids, micelles
- Membranes, films; in-situ adsorption – desorption/ humidifying – drying phenomena
- Kinetics of demixing, formation, aggregation
- Shear induced micelle deformation, rubber network deformation, nanocomposite ordering

- Protein structure and folding/ unfolding
 - Pressure dependence of phase diagrams, fluctuations, molecular interactions
 - In-situ crystallisation semi-crystalline polymer
- Self-assembly of block-copolymers in micellar structures is a widely studied topic at KWS-2. The properties of block-copolymer micelles tuned by changing e.g. solvent quality, temperature, solvent selectivity, block copolymer composition, and molecular weight are investigated thoroughly benefiting from the adjustable instrumental resolution between 2 and 20%.

Another kind of typical application relate to fast structural changes of micellar systems (formation, transformation or chain exchange at equilibrium) or polymer crystallisation which are investigated by time-resolved SANS in the second or sub-second (up to 50 ms) regimes. More recently, the determination and control of the morphological parameters of biocompatible gels and amphiphiles became an important topic of study stimulated by the demands from nanomedicine related to the design of new functional drug delivery vehicles.

Sample Environment

- Anton-Paar fluid rheometer
- Stopped flow cell
- Sample holders: 9 horizontal x 3 vertical (temperature controlled) for standard Hellma cells
- Oil & water thermostats (typical 10 ... 100°C)
- 8-positions thermostated (Peltier) sample holder (-40°C ... 150°C)
- Magnet (1.5 T, vertical)
- Cryostat with sapphire windows
- Pressure cells (500 bar, 5000 bar)
- Humidity chamber, 5% ... 95% for 10°C ... 60°C


Complementary in-situ techniques (optional at sample aperture, see instrument plan)

- FT-IR spectroscopy
- DLS & SLS
- ^3He spin analyzer (SEOP)

[1] Radulescu, A. et al., *J. Phys. Conf. Series* 351, 012026 (2012).

[2] Radulescu, A., Ioffe, A., *Nucl. Inst. Meth. A*, 586, 55 (2008).

[3] Radulescu, A. et al., *Nucl. Inst. Meth. A* 689, 1 (2012)

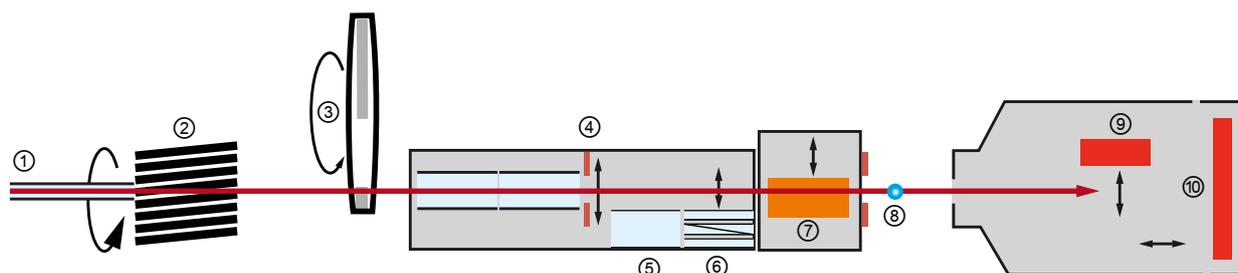

- | | |
|--|---|
| ① Neutron guide | ⑥ Transmission polariser |
| ② Velocity selector $\Delta\lambda/\lambda=20\%$ | ⑦ MgF_2 focussing lenses |
| ③ High-speed chopper $\Delta\lambda/\lambda=1\%$ | ⑧ Sample aperture |
| ④ Entrance aperture | ⑨ High resolution position-sensitive detector |
| ⑤ Neutron guide sections 18 x 1m | ⑩ ^3He tubes array detector |

Technical Data

Overall performance

- $Q = 0.0001 \dots 1 \text{ \AA}^{-1}$
- Maximal flux: $2 \times 10^8 \text{ n cm}^{-2} \text{ s}^{-1}$
- Typical flux: $2.5 \times 10^7 \text{ n cm}^{-2} \text{ s}^{-1}$
(collimation 8 m, $50 \times 50 \text{ mm}^2$, $\lambda = 5 \text{ \AA}$)

Velocity selector

- Astrium, $\Delta\lambda/\lambda = 20\%$, $\lambda = 3 \dots 20 \text{ \AA}$;

Chopper

- Tunable $\Delta\lambda/\lambda$: 20%... 2% (TOF analysis)

Active apertures

- 2 m, 4 m, 8 m, 14 m, 20 m, sample position

Aperture sizes

- Rectangular $1 \times 1 \text{ mm}^2 - 50 \times 50 \text{ mm}^2$

Neutron lenses

- MgF_2 , diameter 50 mm, curvature 20 mm
- Packs with 4, 6, 16 lenses

Polariser

- Transmission, $P > 95\%$ for $\lambda > 4.5 \text{ \AA}$

Sample stage

- XYZ θ translational-rotational stage + cradle
- Accuracy better than 0.01° , 0.01 mm

Detector 1

- Detection range: continuous 1 – 20 m
- ^3He tubes array, active area $\sim 0.9 \text{ m}^2$, count rate for no deadtime $> 2 \text{ MHz}$, resolution = $< 8 \text{ mm}$, stability of pixel response $\sim 0.1\%$, efficiency 85% for 5 \AA

Detector 2 (high res.)

- Spatial resolution $0.45 \times 0.45 \text{ mm}^2$
- Active area: $\varnothing = 8.7 \text{ cm}$
- ^6Li -Scintillator 1 mm thickness
- Fixed position: 17 m after sample position

Dr. Aurel Radulescu

Dr. Noemi Kinga Szekely

Dr. Marie-Sousai Appavou

Email: a.radulescu@fz-juelich.de
Office: +49.(0)89.289.10712
Instrument: +49.(0)89.289.14326/14873

n.szekely@fz-juelich.de
+49.(0)89.289.10739

m.s.appavou@fz-juelich.de
+49.(0)89.289.10747



www.mlz-garching.de/kws-2
DOI: 10.17815/jlsrf-1-27

KWS-3

very small angle scattering diffractometer with focusing mirror



KWS-3 is a very small angle neutron scattering (VSANS) instrument running on the focussing mirror principle. The principle of this instrument is a one-to-one image of an entrance aperture onto a 2D position sensitive detector by neutron reflection from a double-focussing toroidal mirror.

The instrument's standard configuration with a 9.5 m sample-to-detector distance allows performing scattering experiments with a wave vector transfer resolution between $4.0 \cdot 10^{-5}$ and $2.5 \cdot 10^{-3} \text{ \AA}^{-1}$, bridging a gap between Bonse-Hart and pinhole cameras. A second sample position at 1.3 m sample-to-detector distance extends the Q-range of the instrument to $2.0 \cdot 10^{-2} \text{ \AA}^{-1}$ and reaches more than one-decade overlapping with the classical pinhole SANS instruments. Another "mobile" sample position can be installed to adept sophisticated sample environment between 8 and 2 m sample-to-detector distance according to the requested instrumental resolution.

The instrument covers the Q range of small angle light scattering instruments. Especially when samples are turbid due to multiple light scattering, V-SANS gives access to the structural investigation. Thus, the samples do not need to be diluted. The contrast variation method allows for highlighting of particular components.


Small-angle scattering is used for the analysis of structures with sizes just above the atomic scale, between 1 and about 100 nm, which can not be assessed or sufficiently characterised by microscopic techniques. KWS-3 is an important instrument extending the accessible range of scattering angles to very small angles with a superior neutron flux when compared to a conventional instrumental set up with pinhole geometry. Thus, the length scale that can be analysed is extended beyond 10 μm for numerous materials from physics, chemistry, materials science, and life science, such as alloys, diluted chemical solutions, and membrane systems.

Typical Applications

- High-flux bridge between Bonse-Hart and conventional SANS diffractometers
- Colloid science: mixtures of particles, particles of micron size, silicon macropore arrays
- Materials science: filled polymers, cements, microporous media
- Polymer science: constrained systems, emulsion polymerisation
- Bio science: aggregations of bio-molecules, protein complexes, crystallisation of proteins
- Hierarchical structures
- Multilamellar vesicles

Sample Environment

- Anton-Paar fluid rheometer
- Stopped flow cell
- Sample holders:
 - 4 horizontal x 2 vertical (temperature controlled) for standard Hellma cells 404-QX
 - 9 horizontal x 2 vertical (room temperature) for standard Hellma cells 404-QX
- Oil & water thermostats (typical 10 – 100°C)
- Electric thermostat (RT - 200°C)
- 6-positions thermostated (Peltier) sample holder (-40 – 150°C)
- Magnet (2 T, vertical)
- Magnet (5 T, horizontal)
- Cryostat with sapphire windows
- High temperature furnace
- Pressure cells (500 bar, 2000 bar, 5000 bar)

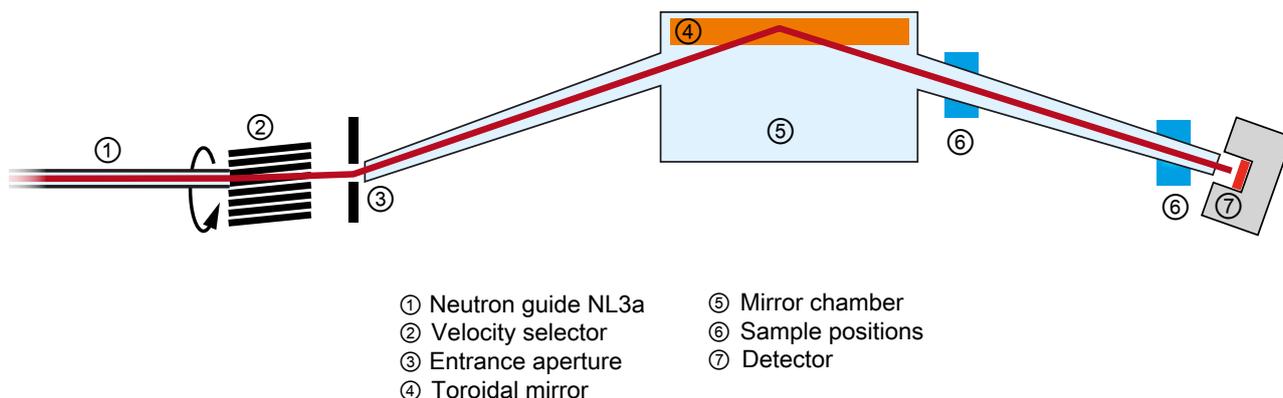

Technical Data

Overall performance

- Resolution:
 $\delta Q = 10^{-4} \text{ \AA}^{-1}$ (extension to $4 \cdot 10^{-5} \text{ \AA}^{-1}$ possible)
- Q-range:
 $1.0 \cdot 10^{-4} - 3 \cdot 10^{-3} \text{ \AA}^{-1}$ at 9.5 m distance
 $1.5 \cdot 10^{-3} - 2 \cdot 10^{-2} \text{ \AA}^{-1}$ at 1.3 m distance
- Neutron flux:
 high-resolution mode: $> 10000 \text{ n s}^{-1}$
 high-intensity mode: $> 60000 \text{ n s}^{-1}$

Monochromator

- MgLi velocity selector
- Wavelength spread $\Delta\lambda/\lambda = 0.2$
- Wavelength range $\lambda = 10 - 30 \text{ \AA}$ (maximal flux at 12.8 \AA)

Aperture size (focus)

- $1 \times 1 \text{ mm}^2 - 5 \times 5 \text{ mm}^2$

Beam size at 9.5 m

- $0 \times 0 \text{ mm}^2 - 100 \times 25 \text{ mm}^2$

Beam size at 1.3 m:

- $0 \times 0 \text{ mm}^2 - 15 \times 10 \text{ mm}^2$

Dr. Vitaliy Pipich

Dr. Zhendong Fu

Email: v.pipich@fz-juelich.de
 Office: +49.(0)89.289.10710
 Instrument: +49.(0)89.289.14873

z.fu@fz-juelich.de
 +49.(0)89.289.10716



www.mlz-garching.de/kws-3
 DOI: 10.17815/jlsrf-1-28

MARIA

magnetic reflectometer with high incident angle



The neutron reflectometer MARIA with polarisation analysis was designed for the investigation of thin magnetic layered structures down to the monolayer scale and lateral structures. The reflection of polarised neutrons allows to determine individually the density and the modulus and the direction of the magnetisation vector of buried layers.

MARIA is optimised for layer thicknesses between 3 – 300 Å and lateral structure sizes from nm to μm sizes. Consequently the instrument is designed for small focused beam and sample sizes of 1 cm² at $\lambda = 4.5 \text{ \AA}$ (available: $4.5 \text{ \AA} < \lambda < 40 \text{ \AA}$) in a vertical orientation with a maximum incident angle of 180° and outgoing angle ranging from -14° to 100°. MARIA provides polarisation analysis in standard operation, where the beam is polarised by a polarising guide (z-geometry; $4.5 \text{ \AA} < \lambda < 10 \text{ \AA}$) and analysed by a wide angle ³He-cell.

Beside the above described reflectometer mode with good resolution in the horizontal scattering plane, MARIA can be used in the GISANS mode with additional resolution in the vertical direction. The latter mode allows for measuring lateral structures down to the nm scale.

At the sample position, a Hexapod with an additional turntable (360°) is installed, which can take a load up to 500 kg. In the standard configuration magnetic fields are provided up to 1.3 T (Bruker electromagnet) and cryogenic temperatures down to 4 K (He

closed cycle cryostat). Beside this standard setup the complete sample environment of the JCNS can be adopted to MARIA so that magnetic fields up to 5 T and temperatures from 50 mK to 500 K are available.

All parts of MARIA are controlled by a computer system according to the “Jülich-Munich” standard based on a Linux workstation. This allows a flexible remote control with automatic scan programs, including the control of sample environment as cryostat and electromagnet.

Typical Applications


With scattering under grazing incidence we investigate depth-resolved the laterally-averaged magnetisations and the correlations between their lateral fluctuations. With an additionally polarised neutron beam we derive a vector information on the laterally-averaged magnetisation (reflectivity) and on the correlations between their lateral fluctuations (off-specular scattering – μm length scale, GISANS – nm length scale).

In general, MARIA can be used for measurements of magnetic roughness, the formation of magnetic domains in thin layered structures, lateral structures, etc. (polarised mode) and density profiles, structures of solid polymer layers, etc. (unpolarised mode with higher intensity).

Furthermore possible without the need for multilayers investigation of:

- Diluted semiconductors
- Influence of the substrate
- Interfaces between oxide materials

Additionally the instrument in non-polarised beam mode can be used for reflectometry and GISANS studies of “soft” layers at the solid/ liquid interface by the use of appropriate liquid cells that are available at the beamline. Candidate systems for such investigations include polymer brushes, polyelectrolyte multilayers, biomimetic supported membranes, adsorbed proteins etc. For typical applications involving deuterated solvents the dynamic range that can be expected covers 7 order of magnitude.

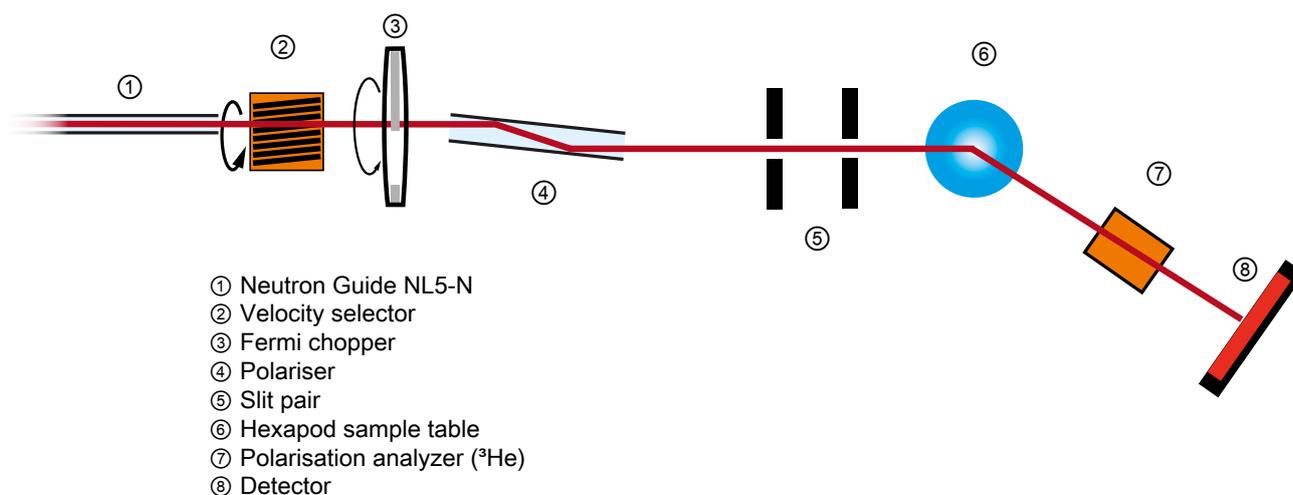

Sample Environment

The optimal sample size for MARIA is $10 \times 10 \text{ mm}^2$ with the following parameters:

- Thin magnetic layers down to sub mono layers
- Polarisation analysis as standard
- Layer thickness of $1 - 300 \text{ \AA}$ optimised, but $- 1000 \text{ \AA}$ (multi layers) should be feasible
- Lateral structures of nm to μm

- Temperature controlled liquid cells for soft matter, accomodating various substrates

Besides the described cryogenic temperatures and magnetic fields MARIA can provide users with a fully equipped Oxid-MBE (**M**olecular **B**eam **E**pitaxy). The typical sample sizes are $10 \times 10 \text{ mm}^2$ and as targets we can provide Al, Cr, Pr, Fe, La, Nb, Ag, Nd, Tb, Sr, Mn, Ti and Co.

Technical Data

Primary beam

- Neutron guide NL5-N: vertically focussing elliptic guide
- Monochromator: Velocity selector
- Wavelength: $4.5 \text{ \AA} - 10 \text{ \AA}$ (polarised)
 $4.5 \text{ \AA} - 40 \text{ \AA}$ (unpolarised)
- Resolution: 10% velocity selector
1%, 3% Fermi chopper
- Double reflection polariser
- Horizontal scattering plane

Flux at sample

- Expected pol. flux $5 \cdot 10^7 \text{ n cm}^{-2} \text{ s}^{-1}$ for 3 mrad collimation

Distances and angles

- 4100 mm distance S1 – S2 (collimation)
- 400 mm distance S2 – sample
- $50 \text{ mm} \times 40 \text{ mm}$ (w \times h) max. opening S2
- 1910 mm distance sample – detector
- 120° maximum detector angle
- GISANS option: 4 m collimation length

Accessible Q-range

- Reflectometry:
 Q_z - range $0.002 \text{ \AA}^{-1} - 3.2 \text{ \AA}^{-1}$
 Q_x - range $6 \cdot 10^{-5} \text{ \AA}^{-1} - 0.001 \text{ \AA}^{-1}$
 α_r - $14^\circ - 100^\circ$
- GISANS option:
 Q_y - range $0.002 \text{ \AA}^{-1} - 0.2 \text{ \AA}^{-1}$

Polarisation analysis

- ^3He cell

Dr. Stefan Mattauch

Dr. Sabine Pütter

Dr. Alexandros Koutsoumpas

Email: s.mattauch@fz-juelich.de

s.puetter@fz-juelich.de

a.koutsoumpas@fz-juelich.de

Office: +49.(0)89.289.10709

+49.(0)89.289.10742

+49.(0)89.289.11674

Instrument: +49.(0)89.289.10799



www.mlz-garching.de/maria

DOI: 10.17815/jlsrf-1-29

NREX

neutron reflectometer with X-ray option



The neutron/ X-ray contrast reflectometer NREX, operated by the Max Planck Institute for Solid State Research, is designed for the determination of structural and magnetic properties of surfaces, interfaces, and thin film systems.

The instrument is an angle-dispersive fixed-wavelength machine with a default wavelength of 4.3 Å. A horizontal focussing monochromator gives the possibility to switch between modes “high intensity/ relaxed resolution” and “high resolution/ reduced intensity” and provides a beam especially for small samples (down to 5 x 5 mm²). A Beryllium filter attenuates higher order reflections. Transmittance supermirrors $m = 3.5$ with a polarising efficiency of $P = 99\%$ and high efficiency gradient RF field spin flippers are used for a full 4 spin channel polarisation analysis.

The sample is aligned horizontally. By tilting the sample the incident angle is varied. The detector arm can move for GISANS horizontally as well as vertically for specular and diffuse scattering measurements. Neutrons are detected with a 20 x 20 cm² position sensitive or a pencil detector. An X-ray reflectometer can be mounted on the sample table orthogonal to the neutron beam. It allows for the in-situ characterisation of sensitive soft matter samples and neutron/ X-ray contrast variation experiments.


Typical Applications

The instrument provides specular and off-specular reflectometry as well as grazing incidence small angle diffraction both in polarised and non-polarised modes. While the specular reflectivity allows determining the scattering length density profiles (20 – 1500 Å) with nm precision along the surface normal, the off specular reflectivity is sensitive to in-plane-inhomogeneity like roughness, (magnetic) domains, vortices in superconductors- and clusters- in the μm -range. To probe lateral (in-plane) structures in the order of atom distance (down to few Å) at the surface, grazing incidence diffraction is provided.

Sample Environment

A closed cycle cryostat (down to 3.5 K) and an electromagnet for fields up to 3.5 T applicable in all three space-directions are provided. Additionally the standard sample environment (magnets up to 7.5 T and ³He inserts for the cryostat down to 50 mK) are available. To the instrument pool belong a cell for investigations at the solid/ liquid interface and a gastight chamber for experiments under defined environmental conditions (arbitrary atmospheres: for example defined relative humidity) at the solid/ air interface.



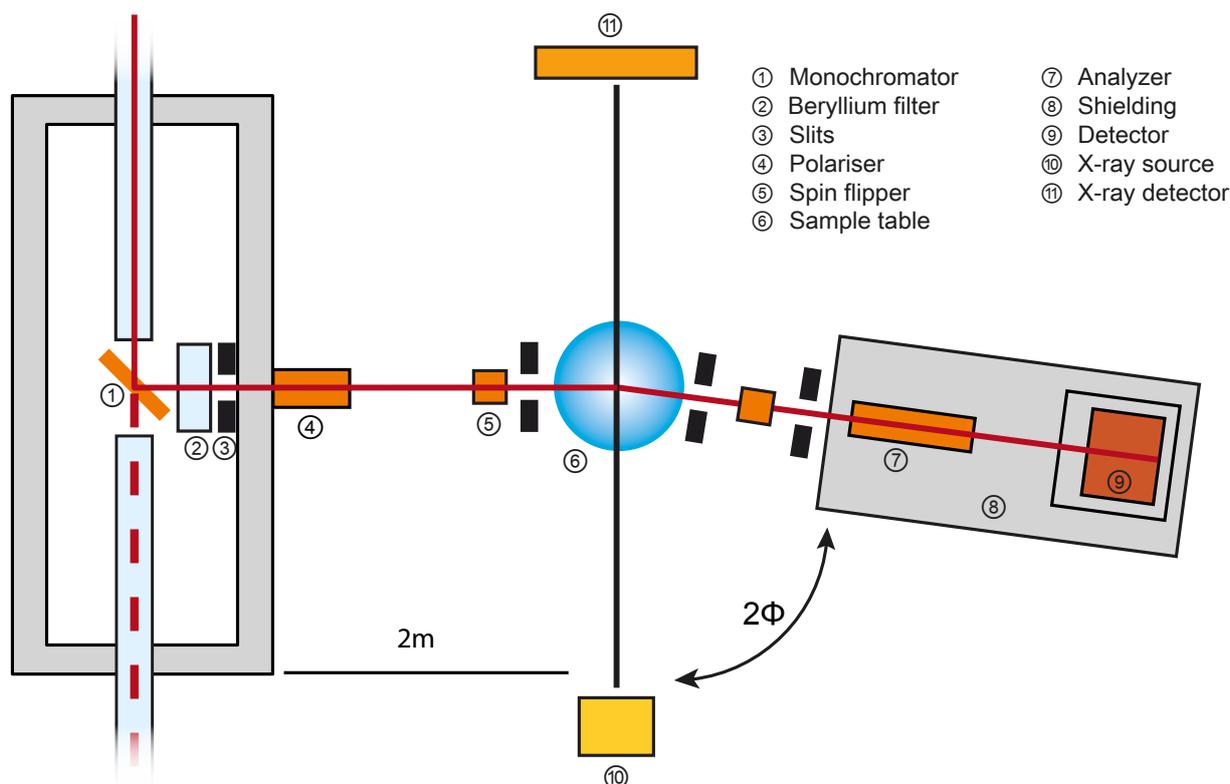

Technical Data

Monochromator

- Type 7 × 5 HOPG crystals
Horizontal focussing
- Wavelength 4.3 Å
- Wavelength resolution 1...2%
- Distance to sample 2500 mm
- Higher order filter cooled Be

Collimation

- Vertical
 - Slit sizes 0.2 – 6 mm
 - Divergence 0.05 – 1.4 mrad
- Horizontal
 - Slit 0.2 – 100 mm

Polarisation

- Beam polarisation > 99%
- Flipper efficiency > 99%

Detector

- 2 pencil detectors ³He
- 2D area detector ³He wire chamber
 - Active area 200 × 200 mm²
 - Lateral resolution 3 mm
- Distance to sample 2465 mm

Dynamical- and Q-Range

- Specular reflectivity 1 : 1 × 10⁻⁶ (@ 5 × 5 mm² sample and full polarisation)
- Q_z 0.005 - 0.5 Å⁻¹
- δQ_z (0 < Q_z < 0.2 Å⁻¹) < 0.002 Å⁻¹

X-rays

- Source: Cu-Kα fixed anode (1.54 Å)
- Monochromator: Goebel Mirror & Double Ge-Crystal
- Detector: 0-dimensional NaI scintillation counter

Dr. Yury Khaydukov

Dr. Olaf Soltwedel

Dr. Thomas Keller

Email: Y.Khaydukov@fkf.mpg.de
Office: +49.(0)89.289.14769
Instrument: +49.(0)89.289.14878

O.Soltwedel@fkf.mpg.de
+49.(0)89.289.14769


T.Keller@fkf.mpg.de
+49.(0)89.289.12164



www.mlz-garching.de/nrex
DOI: 10.17815/jlsrf-1-30

REFSANS

horizontal time-of-flight reflectometer with GISANS option

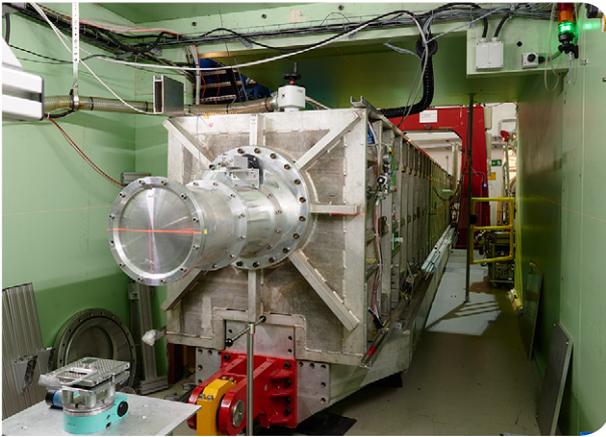

The horizontal reflectometer REFSANS was designed to enable specular reflectometry as well as grazing incidence neutron scattering studies of both solid samples and liquid-air interfaces.

By using a polychromatic incident neutron beam and time-of-flight (TOF) wavelength resolution, REFSANS gives simple access to a large Q range. Typical reflectometry curves are recorded using three incident angles to cover the $0 - 2 \text{ nm}^{-1} Q_z$ domain. In the case of GISANS, the TOF mode provides direct information about the full penetration curve from a single incident angle.

The instrument versatility relies on the one hand on the fact that the wavelength resolution can be tuned between 0.2 and 10%, on the other hand on the possibility to independently control the horizontal and vertical divergence by means of a complex optic. These two characteristics make it possible to optimally perform reflectometry and GISANS. One can easily switch between these two configurations for a given sample and thereby fully investigate its structure without having to alter externally applied fields or constraints (temperature, chemical environment).

For reflectometry, a horizontally smeared out beam of up to 80 mm width is used in order to maximise intensity. For GISANS, up to 13 point beams are impinging on the sample and point focused on the 2D position sensitive detector placed at a distance

of 9 m. This setup allows to resolve lateral structures with dimensions up to several micrometer. In all other cases the detector can be placed at any distance between 1.5 m and 12 m from the sample, thereby making it easy to control the explored angular range and optimise the resolution/ background intensity trade-off.

Typical Applications

The TOF reflectometry and GISANS techniques can be used to characterise thin films in general. Reflectometry provides information about the structure along the sample's normal, while GISANS gathers information about the in-plane correlations. Typical reflectometry experiments include:

- Characterisation of polymer thin film structure and their swelling behaviour in presence of various vapors
- Biological systems such as solid or liquid supported membranes (e.g determination of the morphology and localisation of proteins at interfaces)
- Metallic multilayers (e.g magnetically active films)
- Coatings


GISANS complements these measurements and has been successfully applied to polymer thin films (lateral correlations e.g in dewetted systems, detection and identification of polymer lamellae in immiscible blends or semicrystalline systems), composites, nanopatterned metallic surfaces for which Bragg truncation rods have been reconstructed.

Sample Environment

The optimal sample size is $70 \times 70 \text{ mm}^2$. Various environments are available:

- Simple sample changer for three substrates
- Vibration controlled Langmuir trough for liquid-air interfaces studies
- Magnetic fields up to 7 Tesla
- Cryostats

A heavy load Huber goniometer (max. load 200 kg) is normally used to carry the experimental set-up, but it can easily be removed and replaced by custom equipments.

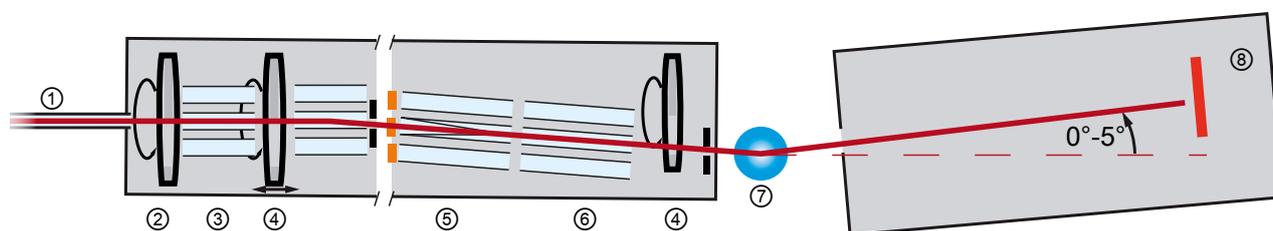

- | | |
|--------------------------|--------------------------|
| ① Neutron guide NL2b | ⑤ Changeable polariser |
| ② Master chopper | ⑥ Neutron guide elements |
| ③ Neutron guide elements | ⑦ Sample position |
| ④ Slave chopper 1+2 | ⑧ Detector |

Technical Data

Primary beam

- Neutron guide NL2b
- Astrium choppers with wavelength resolution to be chosen in the range 0.2 – 10% for wavelengths in the range 2 – 20 Å. Rotation speed up to 6000 rpm
- Collimation: 2 vertical adjustable slits (0 – 12 mm) separated by 8.68 m
- For reflectometry, the horizontal divergence is maximized by use of supermirrors ($m = 2 - 3$)

Flux at sample

Typical values ($\Delta Q / Q = 3\%$):

- $1 \cdot 10^4 \text{ n s}^{-1}$ (incident angle 0.2°)
- $3 \cdot 10^6 \text{ n s}^{-1}$ (at 2.5°)
in the wavelength range 2 to 6 Å for a $60 \times 60 \text{ mm}^2$ sample.

Accessible Q-range

- Reflectometry:
 Q_z up to 0.3 \AA^{-1} for reflectivities down to the 10^{-7} range.
- GISANS:
 $Q_y = 9.5 \cdot 10^{-5} \text{ \AA}^{-1}$ to 0.18 \AA^{-1} (corresponding to distances from 6 \mu m down to 3.5 nm)

Detector

- High performance 2D $500 \times 500 \text{ mm}^2$ multiwire ^3He detector (pixel size 2.7 mm , efficiency 80% at 7 Å, gamma sensitivity $< 10^{-6}$) positioned between 1.5 m and 12 m from the sample. The detector is installed in a liftable vacuum tube in order to reach exit angles up to 6 degrees at the maximum distance.

TOF analysis

- The data are acquired in list mode, each neutron arrival time and impact position being stored for later analysis. This makes it possible to perform various rebinnings in order to tune the resolution/ intensity trade-off.

Dr. Jean-François Moulin Martin Haese

Email: jean-francois.moulin@hzg.de martin.haese@hzg.de
Office: +49.(0)89.289.10762 +49.(0)89.289.10763
Instrument: +49.(0)89.289.14880



www.mlz-garching.de/refsans
DOI: 10.17815/jlsrf-1-31

SANS-1

small angle neutron scattering



The new small angle scattering instrument SANS-1, a joint project of TUM and HZG, has completed commissioning and is in regular user service [1]. SANS-1 is located at the end of neutron guide NL4a in the Neutron Guide Hall West.

SANS-1 is a standard pinhole SANS instrument with both 20 m collimation distance and 20 m sample detector distance, respectively. SANS-1 has been optimised by Monte-Carlo simulations to fit the restrictions in both available space and optimal usage of the provided neutron beam [2]. A vertical S-shaped neutron guide with extreme suppression of fast background neutrons is optimised for complementary wavelength packages, followed by the selector tower with two selectors for high and low resolution, respectively. Adjacent to the selector tower, a collimation system with four parallel horizontal tracks provides vast flexibility: The first track is occupied by a neutron guide system, the second track carries the collimation system with additional background apertures on track three. One track remains empty for various future applications such as focussing lenses or a longitudinal spin echo option. Two Fe/Si transmission polarisers have been optimised to cover the whole wavelength band from 4.5 – 30 Å.

The acentric mounting of the detector tube with around 2.4 m inner diameter allows to use a primary detector of $1 \times 1 \text{ m}^2$ with lateral movement of more than 0.5 m, significantly expanding the accessible Q-range to around $Q_{\text{max}} \approx 1 \text{ \AA}^{-1}$. The primary detector is made up of an array of 128 position sensitive tubes to provide $8 \text{ mm} \times 8 \text{ mm}$ spatial resolution. A second high resolution (3 mm) detector, installed downstream of the primary detector is foreseen for 2016.

A TISANE chopper disk set-up will be available in 2015 which allows to perform kinetic neutron scattering experiments in the μs regime and simultaneously sets the stage for a later upgrade to a complete time-of-flight option for SANS-1.


Typical Applications

The instrument SANS-1 is dedicated to study the structure of materials on length scales of 10 to 3000 Å. With its polarised beam option, the flexible sample goniometer, the wide non-magnetic sample space and the specialised set of sample environment, SANS-1 is particularly adapted for the needs of materials research and magnetism. The precise sample goniometer carries various loads up to 750 kg and fulfills the rising demand on diffraction experiments at low scattering angles, for instance for studies of superconducting vortex lattices and other large magnetically ordered systems.

- Precipitates and segregation in alloys
- Chemical aggregation
- Defects in materials, surfactants, colloids
- Ferromagnetic correlations in magnetism
- Magnetic domains
- Polymers, proteins, biological membranes, viruses, ribosomes and macromolecules
- Superconducting vortex lattices
- Large magnetic structures such as helical magnets and skyrmion lattices

[1] Gilles, R. et al., *Physica B*, 385 386, 1174-1176 (2006).

[2] Gilles, R. et al., *J. Appl. Cryst.*, 40, s428-s432 (2007).

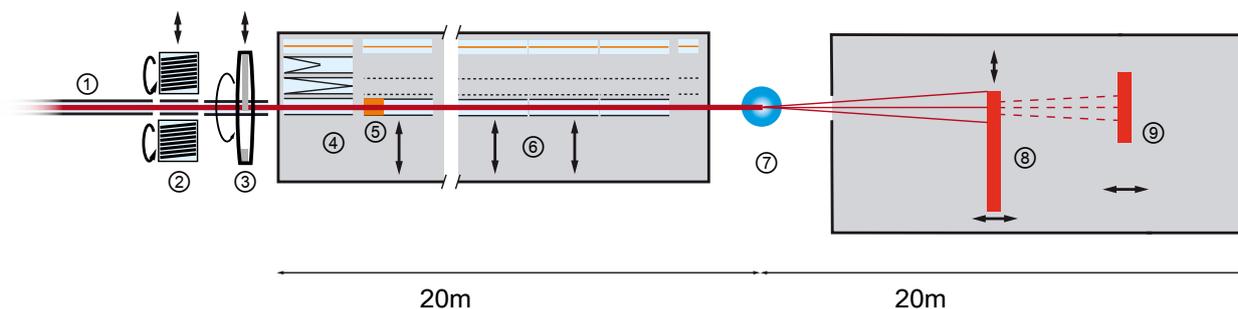

- | | |
|--|--|
| ① Neutron guide NL4a | ⑦ Sample position |
| ② Velocity selector 1+2 | ⑧ Position sensitive area detector, 1 x 1 m ² |
| ③ TISANE Chopper | ⑨ High resolution position-sensitive area detector, 0.5 x 0.5 m ² (installation 2016) |
| ④ Changeable polarisers | |
| ⑤ Spin flipper | |
| ⑥ 4 collimation sections 19 m (neutron guide, collimation slits) | |

Sample Environment

- Standard sample changer with 22 positions
- Different types of high temperature furnaces up to 1900°C
- Deformation-rig with heating
- Set of magnets (5 T horizontal, parallel and perpendicular access, 7.5 T vertical)
- Sample changer with thermostat (-20 +200°C), 11 positions
- Different cryostats with optional ³He insert (460 mK base temp. with 5 T magnet, 50 mK with 7.5 T magnet)
- Polarisation analysis with ³He cell

Technical Data

Primary beam

- S-shaped neutron guide (NL 4a), 50 × 50 mm²
- Mechanical velocity selectors with variable speed
 - 1) $\Delta\lambda/\lambda = 10\%$ medium resolution
 - 2) $\Delta\lambda/\lambda = 6\%$ high resolution
- Wavelength range: 4.5 Å – 30 Å
- TISANE chopper setup with μ s time resolution

Polarisation

- Two V-shaped polarisers

Collimation system (source-to-sample distance)

- 1 m, 2 m, 4 m, 8 m, 12 m, 16 m to 20 m in steps via insertion of neutron guide sections

Sample size

- 0 – 50 mm diameter

Q-range

- $0.0005 \text{ \AA}^{-1} < Q < 1 \text{ \AA}^{-1}$ with primary detector
- $Q_{\min} = 0.0001 \text{ \AA}^{-1}$ with secondary high resolution detector

Detectors

- Primary detector: Array of 128 ³He position-sensitive tubes with an active area of 1000 × 1020 mm² and 8 mm resolution. Lateral detector movement up to 0.5 m, counting rate capability up 1 MHz.
- Secondary high resolution detector (3 mm) and an active area of 500 × 500 mm² to be installed 2016.

Dr. André Heinemann

Dr. Sebastian Mühlbauer

Email: andre.heinemann@hzg.de
 Office: +49.(0)89.289.14534
 Instrument: +49.(0)89.289.12818/14992

sebastian.muehlbauer@frm2.tum.de
 +49.(0)89.289.10784



www.mlz-garching.de/sans-1
 DOI: 10.17815/jlsrf-1-32



DNS
diffuse scattering time-of-flight spectrometer



J-NSE
neutron spin echo spectrometer



PANDA
cold three axes spectrometer



PUMA
thermal three axes spectrometer



RESEDA
resonance spin echo spectrometer



SPHERES
backscattering spectrometer



TOFTOF
cold neutron time-of-flight spectrometer




TRISP
three axes spin echo spectrometer

Spectroscopy

DNS

diffuse scattering neutron time-of-flight spectrometer

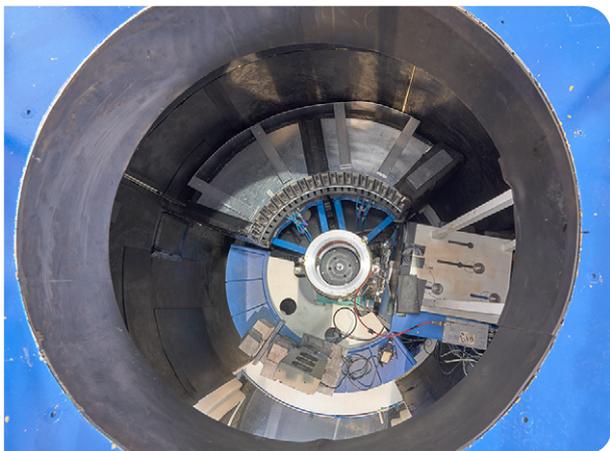

DNS is a versatile diffuse scattering cold neutron time-of-flight spectrometer with polarisation analysis. It allows the unambiguous separation of nuclear coherent, spin incoherent, and magnetic scattering contributions simultaneously over a large range of scattering vector Q and energy transfer E . With its compact size DNS is optimised as a high intensity instrument with medium Q - and E - resolution.

New chopper, neutron velocity selector, and position sensitive detector systems have recently been installed at DNS. This is expected to largely improve possibilities for single-crystal time-of-flight spectroscopy with efficient measurements in all four dimensions of $S(Q,E)$. With its unique combination of single-crystal time-of-flight spectroscopy and polarisation analysis, DNS is also complimentary to many modern polarised cold neutron three axes spectrometers.

Typical Applications


With the increased flux and efficiency delivered by the FRM II, DNS becomes ideal for the studies of complex spin correlations, such as in highly frustrated magnets and strongly correlated electrons, as well as of the structures of soft condensed matter systems, such as the nanoscale confined polymers and proteins, via polarisation analysis. The exploration of unusual magnetic properties can also be efficiently undertaken on single-crystal samples by reciprocal space mapping. In addition to the separation of magnetic cross section from nuclear and spin-incoherent ones, polarisation analysis also al-

lows to distinguish in detail the anisotropy of spin correlations. It has also been well demonstrated that polarised powder diffraction on DNS is complementary to standard neutron powder diffraction and may be extremely useful for magnetic structure refinements, particularly in case of small moments by improving the signal to background ratio. DNS also represents a powerful instrument for the soft condensed matter community for the separation of nuclear coherent scattering from often dominating spin incoherent scattering background. The main applications can be summarised:

- **Application of polarisation analysis:** uniaxial-, longitudinal-, and vector-PA
- **Magnetic, lattice, and polaronic correlations:** geometrically frustrated magnets, strongly correlated electrons, emergent materials
- **Single-crystal and powder time-of-flight spectroscopy:** single-particle excitations, magnons and phonons
- **Soft condensed matters:** separation of coherent scattering from hydrogenous materials, polymer, liquids and glasses

Sample Environment

- Top-loading CCR
- Closed-cycle cold head
- Orange-type cryostat
- Cryo-furnace
- Dilution $^3\text{He}/^4\text{He}$ cryostat insert (~ 20 mK)
- Cryomagnet (self-shielding, vertical field up to 5 T)

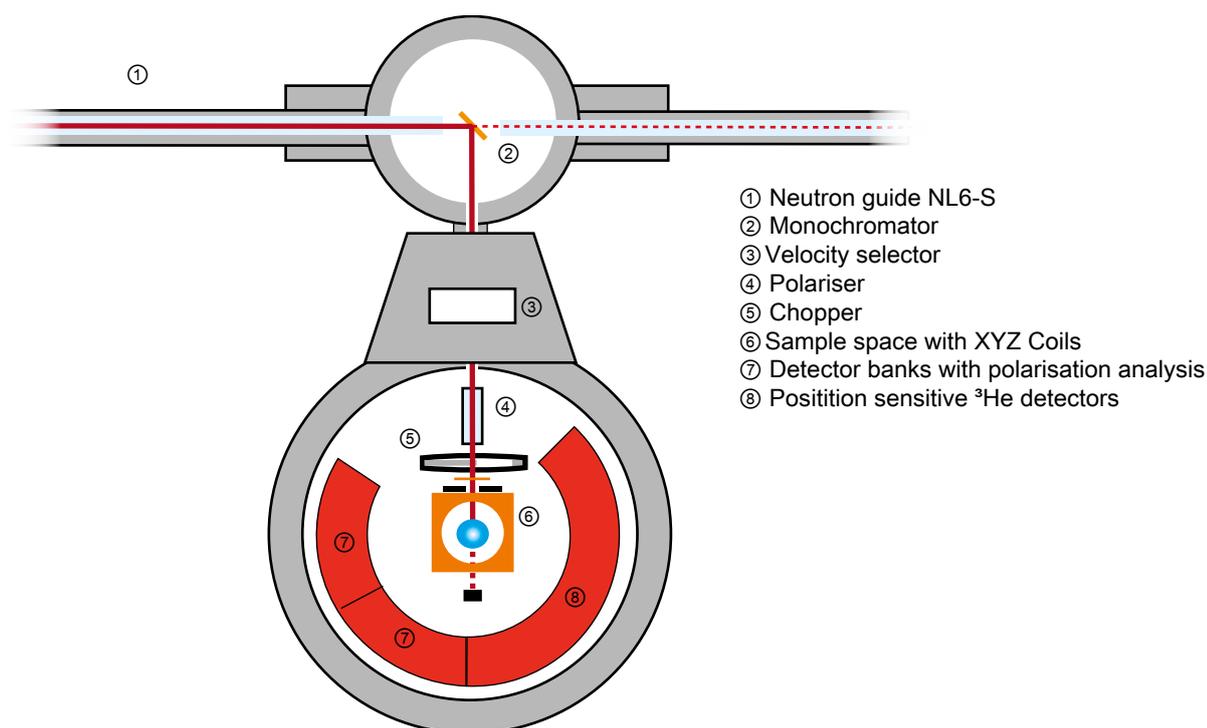

Technical Data

Monochromator

- Neutron guide NL6-S
Horizontally and vertically adjustable, double-focusing
- PG(002), $d = 3.355 \text{ \AA}$
- Crystal dimensions: $2.5 \times 2.5 \text{ cm}^2$ (5 \times 7 crystals)
- Wavelength range: $2.4 \text{ \AA} < \lambda < 6 \text{ \AA}$

Neutron velocity selector

- Minimum wavelength: 1.5 \AA @ 22000 rpm
- Resolution $\Delta\lambda/\lambda$: 30 – 40%

Chopper

- Chopper frequency $\leq 300 \text{ Hz}$
- Repetition rate $\leq 900 \text{ Hz}$
- Chopper disc: Titanium, 3 slits, $\varnothing = 420 \text{ mm}$

Flux at sample

- Non-polarised: $\sim 10^8 \text{ n cm}^{-2} \text{ s}^{-1}$
- Polarised: $\sim 5 \cdot 10^6 - 10^7 \text{ n cm}^{-2} \text{ s}^{-1}$
(polariser: $m = 3$ supermirror benders)

Detector banks for non-polarised neutrons

- 128 position sensitive ^3He tubes
 $\varnothing = 1.27 \text{ cm}$, height $\sim 100 \text{ cm}$
- Total solid angle covered: 1.9 sr
- Covered scattering angles in the horizontal plane: $0^\circ < 2\theta \leq 135^\circ$

Detector banks for polarised neutrons

- 24 detection units:
Polarisation analysis by $m = 3$ supermirror benders
 ^3He detector tubes, $\varnothing = 2.54 \text{ cm}$, height 15 cm
- Covered scattering angle in the horizontal plane: $0^\circ < 2\theta \leq 150^\circ$
- Q_{max}
 $\lambda_i = 2.4 \text{ \AA}$ ($E_i = 14.2 \text{ meV}$): 4.84 \AA^{-1}
 $\lambda_i = 6 \text{ \AA}$ ($E_i = 2.28 \text{ meV}$): 1.93 \AA^{-1}

Energy resolution

- $\lambda_i = 2.4 \text{ \AA}$ ($E_i = 14.2 \text{ meV}$): 1 meV
- $\lambda_i = 6 \text{ \AA}$ ($E_i = 2.28 \text{ meV}$): 0.1 meV

Dr. Yixi Su

Email: y.su@fz-juelich.de
 Office: +49.(0)89.289.10740
 Instrument: +49.(0)89.289.14339 / .14876

Dr. Kirill Nemkovskiy


k.nemkovskiy@fz-juelich.de
 +49.(0)89.289.10779

Dr. Sultan Demirdis

s.demirdis@fz-juelich.de
 +49.(0)89.289.10717



www.mlz-garching.de/dns
 DOI: 10.17815/jlsrf-1-33

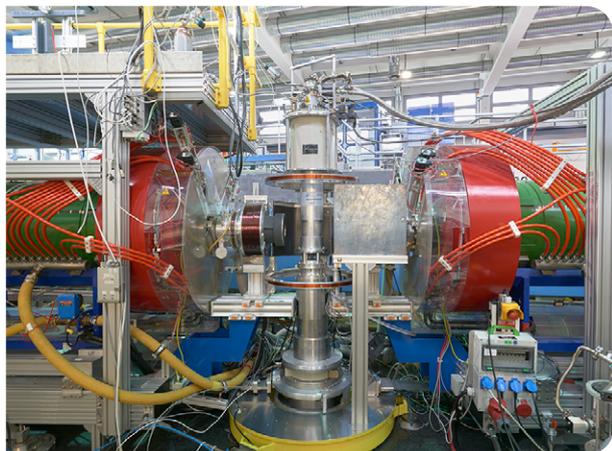

The neutron spin echo technique NSE uses the neutron spin as an indicator of the individual velocity change the neutron suffered when scattered by the sample. Due to this trick, the instrument accepts a broad wavelength band and at the same time is sensitive to velocity changes down to 10^{-5} . However the information carried by the spins can only be retrieved as the modulo of any integer number of spin precessions as intensity modulation proportional to the cosine of a precession angle difference. The measured signal is the cosine transform $S(Q, \tau)$ of the scattering function $S(Q, \omega)$. All spin manipulations only serve to establish this special type of velocity analysis. For details see "Neutron Spin Echo", ed. F. Mezei, Lecture Notes in Physics, Vol. 128, Springer Verlag, Heidelberg, 1980.

Due to the intrinsic Fourier transform property of the NSE instrument it is especially suited for the investigation of relaxation-type motions that contribute at least several percent to the entire scattering intensity at the momentum transfer of interest. In those cases the Fourier transform property yields the desired relaxation function directly without numerical transformation and tedious resolution deconvolution. The resolution of the NSE may be corrected by a simple division.

For a given wavelength the Fourier time range is limited to short times (about 2 ps for J-NSE set-up) by spin depolarisation due to vanishing guide field and to long times by the maximum achievable field integral J . The time is proportional to $J \times \lambda^3$. The J-NSE may achieve a $J = 0.5 \text{ Tm}$ corresponding to $\tau = 48 \text{ ns}$ at $\lambda = 8 \text{ \AA}$.

The instrument itself consists mainly of two large water-cooled copper solenoids that generate the precession field. The precession tracks are limited by the $\pi/2$ -flippers and the π -flipper near the sample position. The embedding fields for the flippers are generated by Helmholtz-type coil pairs around the flipper locations. After leaving the last flipper the neutrons enter an analyzer containing 60 CoTi supermirrors located in a solenoid set. These mirrors reflect only neutrons of one spin direction into the multidetector. By the addition of compensating loops the main coils and the analyzer coil are designed such that the mutual influence of the different spectrometer components is minimised.

Typical Applications


The spin echo spectrometer J-NSE is especially suited for the investigation of slow (~ 1 to 100 ns) relaxation processes. Typical problems from the fields of "soft matter" and glass transition are:

- Thermal fluctuations of surfactant membranes in microemulsions
- Polymer chain dynamics in melts
- Thermally activated domain motion in proteins, which is an important key for understanding the protein function

Sample environment

- Circulation thermostat furnace (260 – 360 K)
- Cryofour (3 – 650 K)
- Furnace (300 – 510 K)
- CO_2 - pressure cell (500 Bar)

Other specialised sample environments on request.

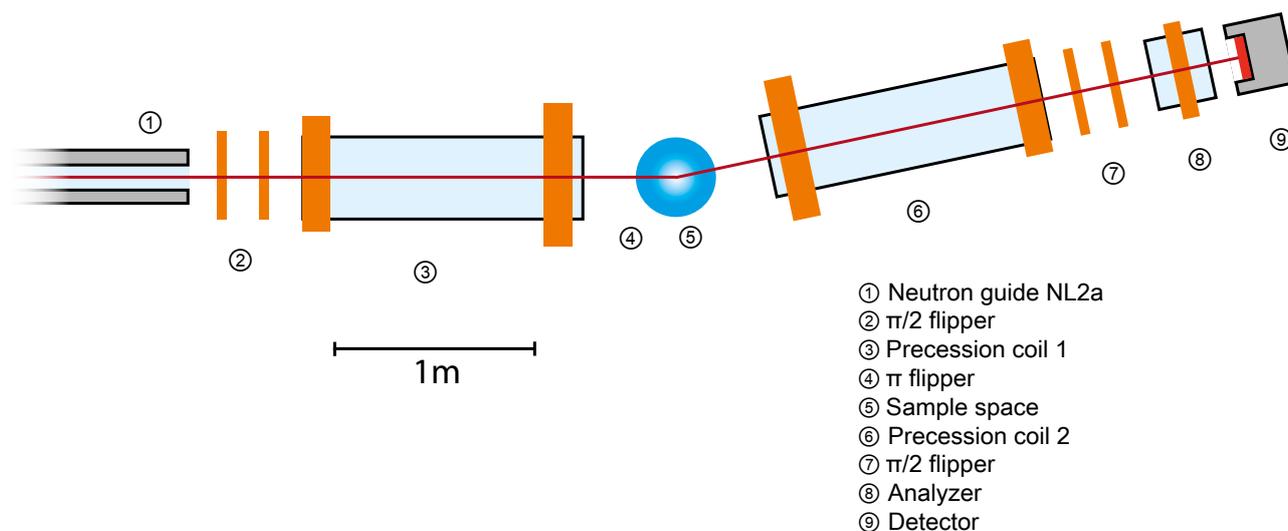

Technical Data

Main parameters

- Polarised neutron flux at sample position
 7 Å: $1 \cdot 10^7 \text{ n cm}^{-2} \text{ s}^{-1}$
 12 Å: $6.8 \cdot 10^5 \text{ n cm}^{-2} \text{ s}^{-1}$
- Momentum transfer range:
 $0.02 - 1.5 \text{ \AA}^{-1}$
- Fourier time range:
 $2 \text{ ps (4.5 \AA)} < \tau < 350 \text{ ns (16 \AA)}$
- Max. field integral: 0.5 Tm

Primary beam

- Neutron guide NL2a
- Polarisation:
 Short wavelength by bent section with FeSi
 $m = 3$ remanent supermirror coating
 Long wavelength by FeSi polariser at entrance
 of the spectrometer
- Cross section of guide: 6 cm × 6 cm
- Max. sample size: 3 cm × 3 cm
- Collimation:
 By source and sample size or wire collimators
 $0.5^\circ \times 0.5^\circ$

Analyzer

- 30 × 30 cm² CoTi supermirror Venetian blind

Detector

- 32 × 32 1 cm² cells ³He multidetector

Dr. Olaf Holderer


Email: o.holderer@fz-juelich.de
 Office: +49.(0)89.289.10707
 Instrument: +49.(0)89.289.14903

Dr. Oxana Ivanova

o.ivanova@fz-juelich.de
 +49.(0)89.289.10730



www.mlz-garching.de/j-nse
 DOI: 10.17815/jlsrf-1-34

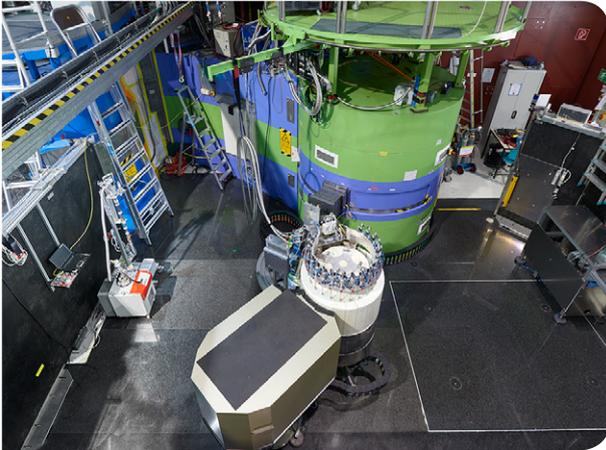

The cold three axes spectrometer PANDA offers high neutron flux over a large dynamic range keeping the instrumental background comparably low. PANDA is situated on the cold neutron beam-tube SR-2 in the Experimental Hall. The high flux is achieved by neutron guide elements in the beam tube, a short source-to-monochromator distance, and the double-focussing monochromator and analyzer crystals. Options for high energy and high q-resolution are available. With dedicated sample environments for very strong magnetic fields and very low temperatures, PANDA is ideally suited for the studies of magnetism and superconductivity on single crystals. Lattice dynamics and magnetic structures are investigated successfully, too.

A polarised neutron set-up using both Heusler monochromator and analyzer as well as a sample-space Helmholtz-coil set for longitudinal polarisation analysis is available.

Typical Applications

Magnetic properties

- Spin-waves
- Crystal field excitations
- Excitations in low dimensional systems
- Magnetic vs nuclear scattering

Lattice dynamics

- Phonon dispersion
- Anharmonic effects
- Polarisations vectors

General

- Critical scattering at phase transitions
- Magnon - phonon interaction
- Soft mode
- Central peak
- Diffraction:
 - without analyser: integral E method
 - with analyser:
 - dE close to 0, high E & Q resolution

Sample Environment

The sample table of PANDA allows for a variety of sample environments, and may easily be adapted to user specific devices. Among other, PANDA disposes routinely operated sample environment for:

Low temperature


- Closed cycle cryostat (3 K < T < 300 K)
 - Variox cryostat (1.5 K < T < 100 K)
 - ³He insert (0.4 K < T < 300 K)
 - Dilution insert (50 mK < T < 6 K)
- typical dimensions for sample space
Ø 50 mm, h = 70 mm

Vertical magnetic field:

- Cryomagnet V15T with optional ³He-⁴He-dilution insert; H_{max} = 13.5 T
(50 mK) 1.5 K < T < 100 K
max. sample diameter: (12 mm) 19 mm
split of coils 20 mm
- Closed-cycle magnet V7.5T
H_{max} = 7.5 T
field at low and high temperatures available

High temperature

- High temperature furnace
300 K < T < 2100 K
sample space: Ø 50 mm, h = 50 mm

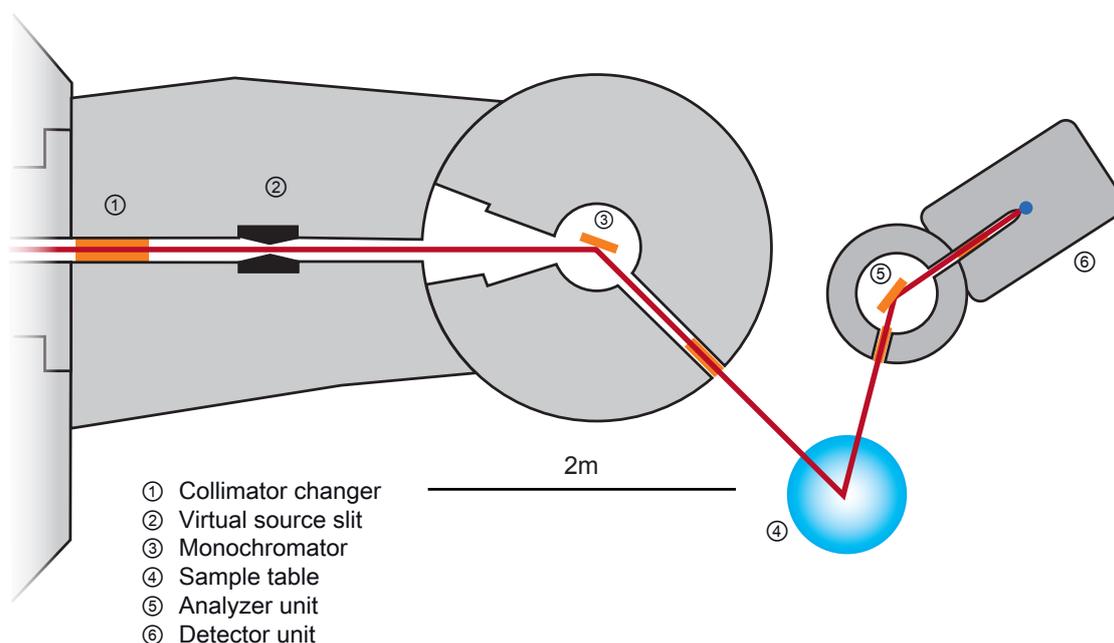

Technical Data

Monochromators

- PG(002) ($d = 3.355 \text{ \AA}$)
 $20^\circ < 2\Theta_M < 132^\circ$
 $1.05 \text{ \AA}^{-1} < k_i < 4.0 \text{ \AA}^{-1}$
 variable horizontal and vertical focussing
- Heusler ($d = 3.35 \text{ \AA}$, polarised neutrons)
 $20^\circ < 2\Theta_M < 120^\circ$
 $1.1 \text{ \AA}^{-1} < k_i < 4.0 \text{ \AA}^{-1}$
 variable vertical focussing

Analysers

- PG(002)
 $-130^\circ < 2\Theta_A < 100^\circ$
 $1.05 \text{ \AA}^{-1} < k_f$
 variable horizontal focussing
- Heusler (polarised neutrons)
 $-130^\circ < 2\Theta_A < 100^\circ$
 $1.05 \text{ \AA}^{-1} < k_f$
 variable horizontal focussing

Detectors

- 1st ^3He tube (focussing mode)
- 2nd ^3He tube (collimated mode))

Flux at sample

PG monochromator vertically focussed, horizontal flat, no collimation:

- $1.9 \cdot 10^7 \text{ n cm}^{-2} \text{ s}^{-1}$ for $k_i = 1.55 \text{ \AA}^{-1}$ Be Filter
- $5.5 \cdot 10^7 \text{ n cm}^{-2} \text{ s}^{-1}$ for $k_i = 2.662 \text{ \AA}^{-1}$ PG Filter

Main parameters

- Scattering angle at the sample:
 $5^\circ < 2\Theta_s < 125^\circ$ (moveable beam-stop)
- Energy transfer
 up to 20 meV
- Momentum transfers
 up to $Q = 6 \text{ \AA}^{-1}$ (depending on k_i)

Filters for higher order suppression:

- PG ($l = 60 \text{ mm}$); $k_f = 2.57 \text{ \AA}^{-1}$ or 2.662 \AA^{-1}
- Be (closed-cycle cryostat, $T \leq 45 \text{ K}$);
 $k_f = 1.55 \text{ \AA}^{-1}$
- BeO (liq.- N_2 cooled); $k_f = 1.33 \text{ \AA}^{-1}$

Dr. Astrid Schneidewind


Dr. Petr Čermák

Email: a.schneidewind@fz-juelich.de
 Office: +49.(0)89.289.14749
 Instrument: +49.(0)89.289.14869 / .14813

p.cermak@fz-juelich.de
 +49.(0)89.289.11773



www.mlz-garching.de/panda
 DOI: 10.17815/jlsrf-1-35

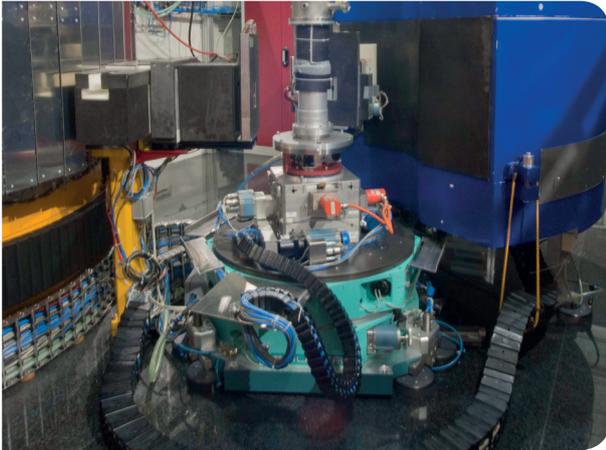

Three axes spectrometers allow the direct measurement of the scattering function $S(Q, \omega)$ in single crystals at well defined points of the reciprocal lattice vector Q and frequency ω and thus represent the most general instrument type.

PUMA is characterised by a very high neutron flux as a result of the efficient use of focussing techniques. Three different vertical openings and a horizontal slit with a maximum opening of 40 mm define the virtual source, which is two meters before the monochromator. To reduce the primary beam's contamination by epithermal neutrons, a sapphire filter can be placed in front of the monochromator. PUMA has a remote controlled monochromator changing unit which allows to place one out of four different monochromators inside the drum. All of them are equipped with double focussing devices that allow for optimum focussing conditions over a wide range of incident wavevectors k_i . The horizontal divergence of the beam can be defined using a series of four Soller collimators. The two inside the drum, before and after the monochromator, can be remotely changed, whereas the two in the analyzer housing can be changed manually. An Eulerian cradle can optionally be used to access the four dimensional Q - ω -space.

An innovative option of the spectrometer is the multianalyzer/ detector system. It allows a unique and flexible type of multiplexing. Using this option a scattering angle range of 16° can be measured simultaneously and flexible Q - ω paths can be realised without repositioning the instrument. Mapping

of excitations is equally well possible as kinetic single shot experiments on time scales that have not been accessible so far.

A unique feature of the instrument is the possibility to perform stroboscopic, time resolved measurements of both elastic and inelastic signals on time scales down to the microsecond regime. Using this technique, the sample is periodically perturbed by an external variable such as temperature, electric field, etc. The signal is then recorded not only as a function of momentum and energy transfer, but also given a time stamp, relative to the periodic perturbation.


Typical Applications

- Phonons
 - Electron-phonon interaction
 - Phonon anharmonicities
 - Soft mode phase transitions
- Magnons
 - Spin waves in (anti)ferromagnets
 - Electron-magnon interaction
 - Unconventional superconductors
 - Crystal fields
- Time resolved/ stroboscopic measurements
 - Temperature cycling (excitations during demixing processes)
 - Electrical field cycling (polarisation processes in ferroelectrics)
 - Temperature/ pressure cycling
- Diffraction; purely elastic signals
 - Superstructures/ satellites
 - Diffuse scattering

Sample Environment

Besides standard sample environment, we provide:

- Closed-cycle cryostates 3.5 – 300 K; – 650 K with adaptable heating device
 - Cryofurnace 5 K – 750 K
 - Paris-Edinburgh type pressure cell $p < 10$ GPa
- Along with the detector electronics required for time resolved measurements, special sample environment for the rapid cycling is available:
- Furnace for fast temperature jumps (~ 5 K/s cooling rate; < 620 K; ambient atmosphere)
 - Switchable HV power supply (< 500 Hz; ± 10 KV)

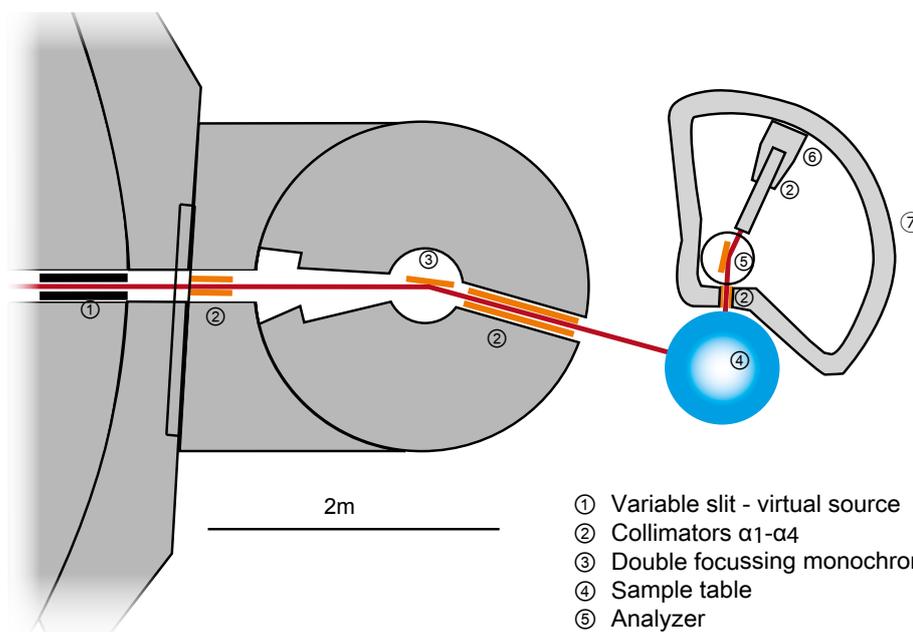

- ① Variable slit - virtual source
- ② Collimators α_1 - α_4
- ③ Double focussing monochromators (PG, Si, Cu, Ge)
- ④ Sample table
- ⑤ Analyzer
- ⑥ Detector
- ⑦ Common analyzer and detector shielding

Technical Data

Primary beam

- Beam tube SR-7 (thermal)
- Beam tube entrance $140 \times 90 \text{ mm}^2$
- Virtual source dimensions:
 - horizontal: 0 – 40 mm
 - vertical: (90, 110, 130 mm)

Distances

- Beam tube entrance - monochromator: 5.5 m
- Virtual source – monochromator: 2.0 m
- Monochromator – sample: $2.0 (\pm 0.1) \text{ m}$
- Sample – analyzer: $1.0 (\pm 0.1) \text{ m}$
- Analyzer – detector: 0.9 m

Collimation

Remote controlled:

- α_1 : 20', 40', 60'
- α_2 : 14', 20', 24', 30', 45', 60'

Manually changeable:

- α_3 : 10', 20', 30', 45', 60'
- α_4 : 10', 30', 45', 60'

Monochromators

- Crystals: PG(002), Cu(220), Cu(111), Ge(311)
size: $260 \times 162 \text{ mm}^2$;
- Focus vertically and horizontally adaptable to incident energy

Analyzer

- Crystals : PG(002), Ge(311); $210 \times 150 \text{ mm}^2$
vertical fixed focus
horizontally adaptable to incident energy

Sample table

- Diameter 800 mm
- Max. load 900 kg
- Amagnetic goniometer ($\pm 15^\circ$)
- Z translation ($\pm 20 \text{ mm}$)
- Optional Eulerian cradle

Main parameters

- Monochromator take-off angle
 $-15^\circ < 2\Theta < -115^\circ$
- Scattering angle sample
 $-70^\circ < 2\Theta < 120^\circ$
(dependent on monochromator take-off angle)
- Analyzer scattering angle
 $-120^\circ < 2\Theta < 120^\circ$
- Incident energy range 5 meV – 160 meV
- Momentum transfer range < 12 \AA^{-1}
- Energy transfer < 100 meV


Dr. Oleg Sobolev

Dr. Jitae Park

Email: oleg.sobolev@frm2.tum.de jitae.park@frm2.tum.de
Office: +49.(0)89.289.14754 +49.(0)89.289.13983
Instrument: +49.(0)89.289.14914



www.mlz-garching.de/puma
DOI: 10.17815/jlsrf-1-36

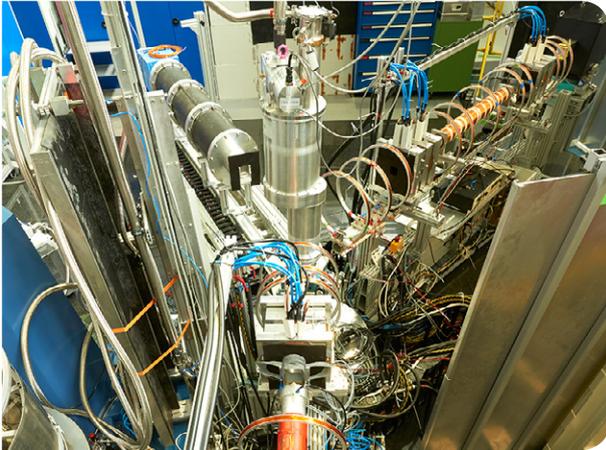

RESEDA (resonance spin echo for diverse applications) is a high-resolution resonance spin-echo spectrometer installed at the cold neutron guide NL5-S in the Neutron Guide Hall West. The instrument gives access to a large time and scattering vector range for quasi-elastic measurements.

RESEDA supports longitudinal neutron resonance spin echo (LNRSE, time range from 0.001 to 20 ns for $\lambda = 8 \text{ \AA}$) and modulation of intensity with zero-effort (MIEZE, time range from 0.001 to 20 ns for $\lambda = 8 \text{ \AA}$) experiments. At RESEDA, the analysis of $S(Q, \tau)$ provides characteristic parameters, e.g. relaxation time and amplitude of the dynamic processes in the sample investigated. The determination of $S(Q, \tau)$ is feasible for different Q-values and/ or different temperatures and pressures.

NRSE experiments require non-depolarising sample environment conditions. For MIEZE experiments (and in contrast to NRSE) the spin manipulation and analysis is realised solely before the sample. Therefore, the MIEZE method enables high-resolution study of depolarising samples, under magnetic field and/ or within depolarising sample environments. However, as a consequence of the polarisation analysis before the sample, MIEZE experiments are limited to a smaller Q-range than NRSE measurements.

Next to ^3He detectors, a 2D CASCADE detector with an active area of $20 \text{ cm} \times 20 \text{ cm}$ characterised by a spatial resolution of 2.6 mm^2 and a time dynamics of the order of a few MHz is available [1, 2].

Hence, RESEDA is in addition suited to (polarised) small angle neutron scattering (SANS) applications.

Typical Applications

- Quasi-elastic measurements:
 - e.g. to determine the dynamics of water in porous media, polymer melts, diffusion processes in ionic liquids as well as magnetic fluctuations in single crystals, powder samples and thin films
- (Polarised) SANS:
 - e.g. to investigate the diffraction pattern of magnetic structures and vortex lattices to choose suited reflections for a line-width determination
- Spherical polarisation analysis

Sample Environment

At RESEDA the whole sample environment of the MLZ is applicable. Depolarising conditions are limited to MIEZE experiments.


- Available temperature range: 50 mK (dilution insert, see below) up to more than 1300 K (high temperature furnace, non-depolarising)
- Maximal pressure: 7 GPa
- Maximal magnetic field: 7.5 T

Available cryostats:

- Closed cycle cryostat: ($3 \text{ K} < T < 300 \text{ K}$)
- ^3He insert: ($450 \text{ mK} < T < 300 \text{ K}$)
- Dilution insert: ($50 \text{ mK} < T < 6 \text{ K}$)

[1] Schmidt, C.J. et al., *Journal of Physics: Conf. Ser.* 251, 012067 (2010).

[2] Hauler, W. et al., *Rev. Sci. Instr.* 82, 045101 (2011).

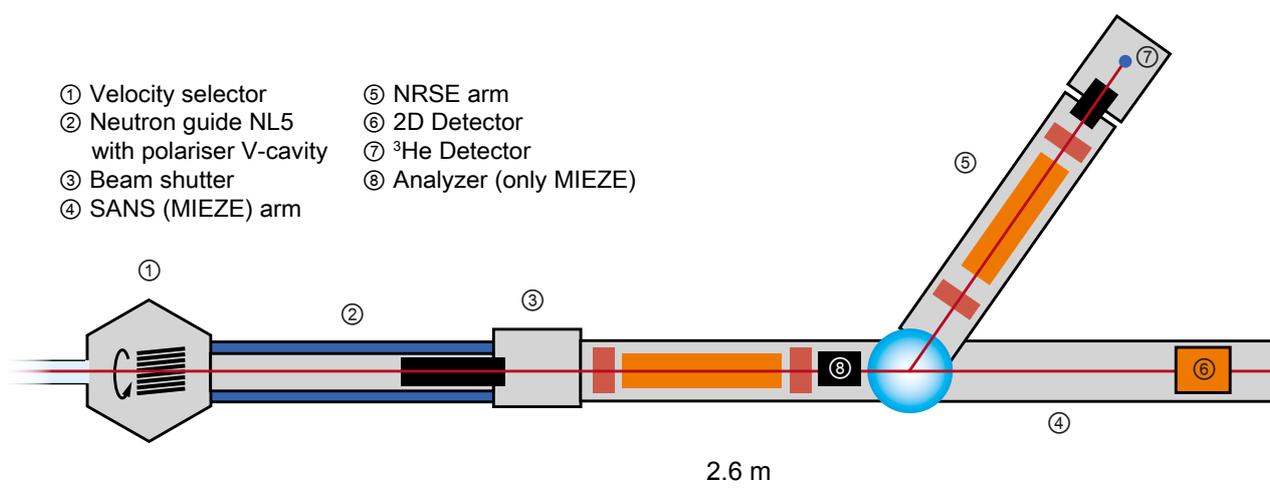

Technical Data

Primary beam

- Neutron guide: NL5-S
- Guide cross section: 29 x 34 mm²
- Wavelength selection: Velocity selector (max. 28000 rpm)
- Wavelength range: $\lambda = 3 - 12 \text{ \AA}$
- Wavelength bandwidth at sample position: $\Delta\lambda/\lambda = 9 - 20\%$
- Polariser: V-cavity (length: 2 m, coating: $m = 3$)

Spectrometer

- Optional polariser before sample: V-cavity (length: 30 cm, coating: $m = 4$)
- Length of the spectrometer arms: 2.6 m
- Two secondary spectrometer arms: SANS (MIEZE) arm and LNRSE arm
- For polarisation analysis available: V-cavity, Bender
- Detectors: ³He counter or 2D detector (CASCADE)

Characteristic parameters

- Flux at sample position: $\phi \geq 10^6 \text{ n cm}^{-2} \text{ s}^{-1}$ (at $\lambda = 5.3 \text{ \AA}$)
- Frequency range of RF coils: 35 kHz – 1.7 MHz
- Maximum scattering angle: $2\theta = 93^\circ$
- Maximum scattering vector: $Q = 2.5 \text{ \AA}^{-1}$ (at $\lambda = 3 \text{ \AA}$)
- Spin echo time range: $\tau = 0.001 - 20 \text{ ns}$ for $\lambda = 8 \text{ \AA}$
- Energy resolution: 0.03 $\mu\text{eV} - 0.1 \text{ meV}$

Dr. Christian Franz

Email: christian.franz@frm2.tum.de
 Office: +49.(0)89.289.14760
 Instrument: +49.(0)89.289.14874

Dr. Thorsten Schröder

thorsten.schroeder@frm2.tum.de
 +49.(0)89.289.14760



www.mlz-garching.de/reseda
 DOI: 10.17815/jlsrf-1-37

SPHERES

backscattering spectrometer



SPHERES (**s**pectrometer for high energy **r**esolution) is a third-generation backscattering instrument with focussing optics and a phase-space-transform chopper. It is a versatile spectrometer for the investigation of atomic and molecular dynamics on a GHz scale.

The necessary filtering of neutron energies is achieved by Bragg reflection from perfect monochromator and analyzer crystals under angles close to 180° . The backscattering geometry makes it unavoidable to use a primary beam deflector and a duty-cycle chopper. At SPHERES, these two functions are both realised by a chopper that bears deflector crystals on its circumference. This leads to a particularly compact spectrometer layout so that full use can be made of the focussing neutron guide. As an additional advantage, the fast motion of the deflector crystals achieves a phase-space transform of the primary spectrum, thereby enhancing the flux at the sample.

The principal figures of merit qualify SPHERES as one of the best of its class [1]. Count rates and signal-to-noise ratio have been improved by filling the instrument housing with argon, thereby avoiding air scattering in the secondary spectrometer. Another gain in flux will be achieved by a more efficient phase-space transform chopper which is in the commissioning phase. The new designed chopper will be more efficient due to optimised rotation

speed and higher reflectivity and mosaicity of the graphite crystals. The resolution of the small angle detectors have been improved by reducing the azimuth angle range of the analyzers [2].

As a multi-detector instrument with relaxed angular resolution, SPHERES is particularly suited for studying tagged-particle motion by incoherent scattering. Typical applications include for example dynamical processes in polymers and biological systems [3]. The high resolution of the spectrometer allows to investigate the dynamics of water in confined geometry. The unprecedented sensitivity of SPHERES helps us to detect the onset of quasielastic scattering deep in the supercooled state [4]. Other important applications are hyperfine splitting in magnetic materials [5] and rotational tunneling [6]. The high count rates allow inelastic temperature scans [7] and real-time kinetic experiments [8].

Raw histograms are accumulated on an equidistant ω grid. A script driven program, SLAW [9], is provided to normalise the raw counts, to perform optional binning, and to deliver $S(q, \omega)$ in a variety of output formats so that users are not bound to any specific data analysis program. In data fitting, it is critically important to convolute theoretical models with the measured resolution function in an efficient and numerically stable way. We strive to support best practice through our FRIDA package [10].

[1] Wuttke, J. et al., *Rev. Sci. Instrum.* 83, 075109 (2012).

[2] Wuttke, J., Zamponi, M., *Rev. Sci. Instrum.* 84, 115108 (2013).

[3] Gallat, F.-X. et al., *J. Am. Chem. Soc.* 134, 13168 (2012).

[4] Doster, W. et al., *Phys. Rev. Lett.* 104, 098101 (2010).

[5] Chatterji, T. et al., *Phys. Rev. B* 78, 012411 (2008).


[6] Bator, G. et al., *Chem. Phys.* 410, 55 (2013).

[7] Häußler, W. et al., *Neutron News* 22, 24 (2011).

[8] Leon, A., Wuttke, J., *J. Phys.: Condens Matter* 23, 254214 (2011).

[9] Wuttke, J.: SLAW - a neutron histogram to scattering law converter, <http://apps.jcns.fz-juelich.de/sl原因>

[10] Wuttke, J.: FRIDA - fast reliable interactive data analysis, <http://apps.jcns.fz-juelich.de/doku/frida/start>

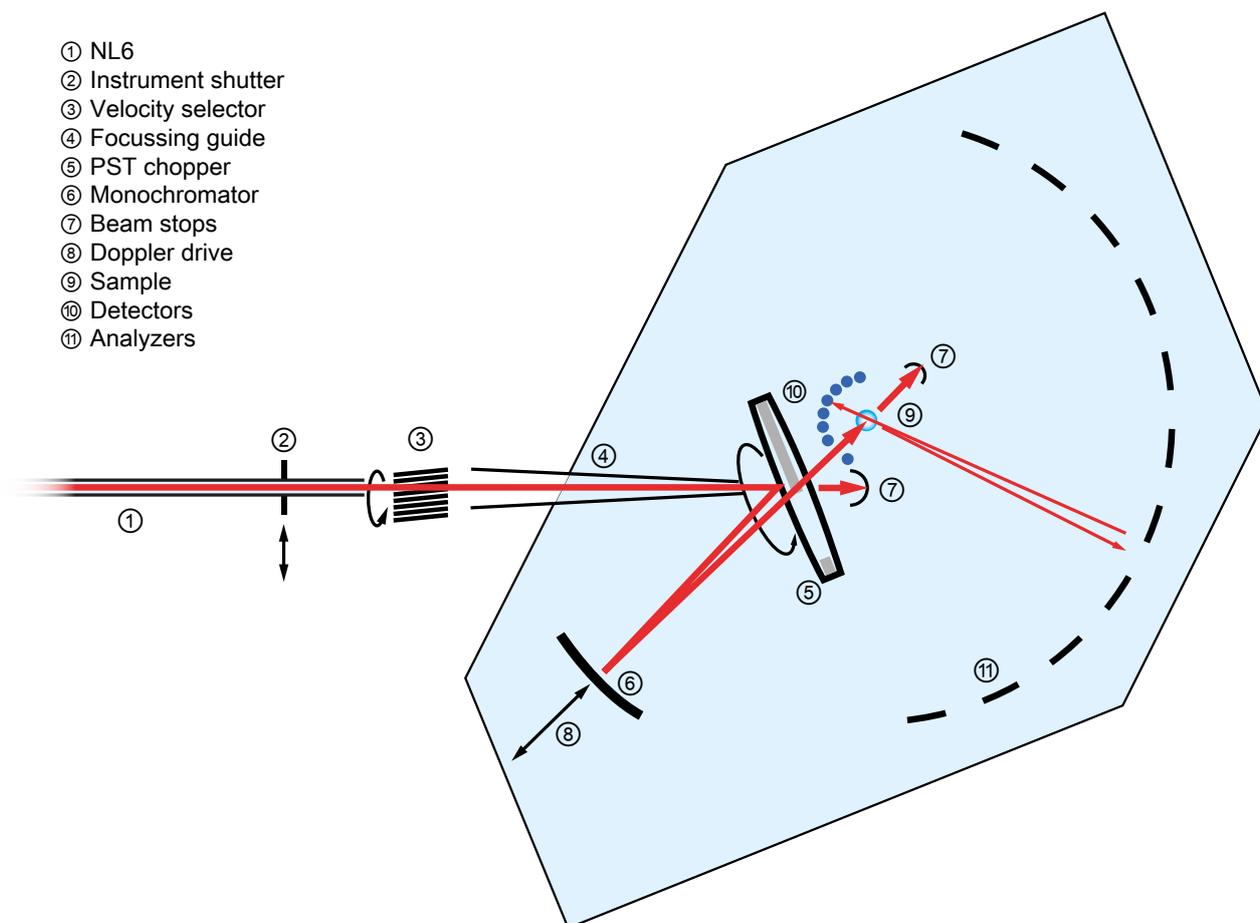

Typical Applications

- Hyperfine splitting
- Molecular reorientations and rotational tunneling
- Dynamic signature of phase transitions
- Hydrogen diffusion
- Liquid dynamics
- Polymer relaxation
- Protein aggregation

Sample Environment

- Cryofurnace 3...700 K
- Dilution inset 20 mK
- Furnace

Technical Data

Primary beam

- | | |
|----------------------|----------|
| • Neutron guide | NL6-S |
| • Neutron wavelength | 6.27 Å |
| • Neutron energy | 2.08 meV |

Main parameters

- | | |
|-----------------------|---------------------------------|
| • Resolution FWHM | 0.62 – 0.65 μeV |
| • Dynamic range | $\pm 31 \mu\text{eV}$ |
| • Q range | 0.2 – 1.8 \AA^{-1} |
| • Flux after selector | 10^{10} s^{-1} |
| • Flux at sample | $1.8 \cdot 10^6 \text{ s}^{-1}$ |
| • Illuminated area | $40 \times 30 \text{ mm}^2$ |

Dr. Michaela Zamponi

Email: m.zamponi@fz-juelich.de
 Office: +49.(0)89.289.10793
 Instrument: +49.(0)89.289.14875

Dr. Marina Khanefit

m.khanefit@fz-juelich.de
 +49.(0)89.289.11676



www.mlz-garching.de/spheres

DOI: 10.17815/jlsrf-1-38

TOFTOF

cold neutron time-of-flight spectrometer



TOFTOF is a direct geometry disc-chopper time-of-flight spectrometer located in the Neutron Guide Hall West. It is suitable for both inelastic and quasi-elastic neutron scattering and the scientific questions addressed range from the dynamics in disordered materials in hard and soft condensed matter systems (such as polymer melts, glasses, molecular liquids, or liquid metal alloys), properties of new hydrogen storage materials to low-energy magnetic excitations in multiferroic compounds, and molecular magnets.

A cascade of seven fast rotating disc choppers which are housed in four chopper vessels is used to prepare a monochromatic pulsed beam which is focussed onto the sample by a converging super-mirror section. The scattered neutrons are detected by 1000 ^3He detector tubes with a time resolution up to 50 ns. The detectors are mounted at a distance of 4 m and cover 12 m² (or 0.75 sr). The high rotation speed of the chopper system (up to 22 000 rpm) together with a high neutron flux in the wavelength range of 1.4 -14 Å allows free tuning of the energy resolution between 3 meV and 2 μeV .

The 60 m primary neutron guide has an s-shape which efficiently suppresses fast neutron background. This enables the investigation of weak signals. The prototype of a new focussing neutron guide has been installed recently, as alternative option in the last section of the guide system. The ex-

isting linearly tapered neutron guide yields a beam spot size of 23 x 47 mm². Using the focussing guide, an intensity gain up to a factor of 3 (wavelength dependent) is observed on a sample area of 10 x 10 mm².

Typical Applications

TOFTOF represents a versatile instrument combining high energy resolution, high neutron flux (also at short wavelengths), and an excellent signal-to-background ratio. It is perfectly suited for both inelastic and quasielastic neutron scattering and scientific topics include e.g.:

- Diffusion in liquid metals and alloys
- Hydrogen dynamics in soft matter systems such as molecular liquids, polymer melts or colloids
- Molecular magnetism, quantum criticality in heavy fermion compounds, low energy excitations in multiferroic materials and novel magnetic phases
- Dynamic properties of energy storage materials, such as solid state hydrogen storage materials, electrolytes for batteries and fuel cells, or gas storage materials
- Energy-resolved quasi-elastic neutron scattering on proteins, vesicles, and biological materials
- Kinetic studies of hydrogen binding, e.g. in concrete
- Aging effects in disordered media and low frequency dynamics in glasses
- Biological activity and functionality of proteins and cells under pressure

Sample Environment

Standard sample environment:


- CCR Cryostat (4 - 600 K)
- ^3He insertion device (down to 0.5 K)
- Circulation thermostat furnace (255 - 450 K)
- High temperature furnace (300 – 2100 K)
- 2.5 T magnet

Sample environment provided by collaborators:

- Electromagnetic levitator
- Electrostatic levitator
- Hydraulic pressure cells (up to 3.5 kbar)
- Clamp pressure cells (few GPa)

- ① Neutron guide NL2a-u
- ② PCR chopper-pair
- ③ 1st higher order removal chopper
- ④ 2nd higher order removal and frame overlap chopper
- ⑤ MCR chopper-pair

- ⑥ Sample position
- ⑦ Radial collimator
- ⑧ Beam stop
- ⑨ Shielding
- ⑩ Detector bench

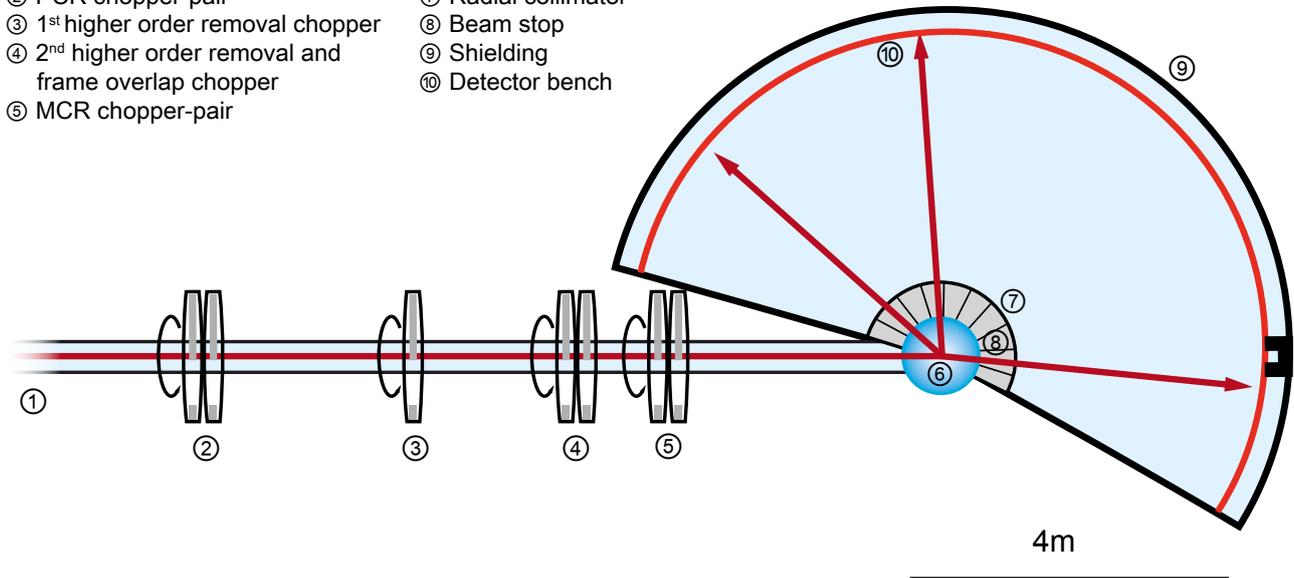

Technical Data

Primary Beam

- | | |
|---|---|
| • Neutron guide | NL2a-u |
| • Number of chopper discs | 7 |
| • Chopper frequency range | 400 min ⁻¹ – 22000 min ⁻¹ |
| • Diameter of chopper disc | 600 mm |
| • Cross section of neutron guide at the entrance | 44 × 100 mm ² |
| • Cross section of neutron guide, 20 cm in front of sample position | 23 × 47 mm ² |
| • Cross section of focussing guide | minimal 12 × 25 mm ² |

Main Parameters

- | | |
|--|---|
| • Adjustable range of incident neutrons | 1.4 – 16 Å |
| • Elastic energy resolution | 2 μeV – 3 meV |
| • Range of energy transfers | -30 meV – 50 meV |
| • Integral neutron flux of the white beam at sample position | 10 ¹⁰ n cm ⁻² s ⁻¹ |
| • Angular range of the detector bank | -15° to -7° and 7° to 140° |

Dr. Wiebke Lohstroh

Email: wiebke.lohstroh@frm2.tum.de
 Office: +49.(0)89.289.14735
 Instrument: +49.(0)89.289.14881

Dr. Zachary Evenson


zachary.evenson@frm2.tum.de
 +49.(0)89.289.14975



www.mlz-garching.de/toftof
 DOI: 10.17815/jlsrf-1-40

TRISP

three axes spin echo spectrometer

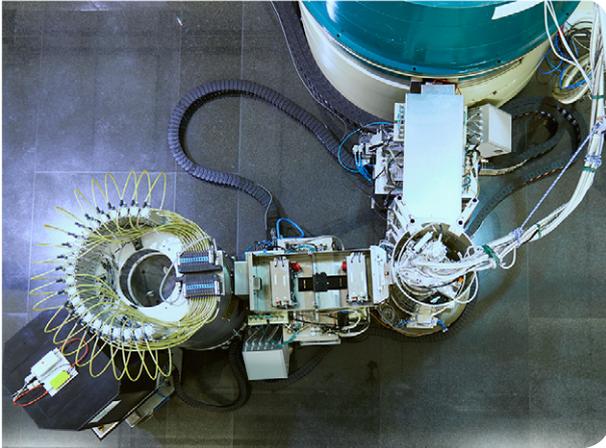

TRISP is a high-resolution neutron spectrometer combining the three axes and neutron resonance spin echo (NRSE) techniques. The design of TRISP is optimised for the study of intrinsic linewidths of elementary excitations (phonons, magnons) with an energy resolution in the μeV region over a broad range of momentum and energy transfers. Compared to conventional three axes spectrometers (TAS), this corresponds to an improvement of the energy resolution of one to two orders of magnitude.

TRISP also incorporates the Larmor diffraction (LD) technique, which allows to measure lattice spacings with a relative resolution $\Delta d/d = 1.5 \cdot 10^{-6}$, i.e. one to two orders of magnitude better than conventional neutron or X-ray diffraction. Absolute d -values can be determined by calibrating the instrument against an Si standard. The main applications of LD include thermal expansion under pressure and low or high temperature, and distributions of lattice constants (second order stresses). LD thus is unique in a parameter region, where standard methods such as dilatometry fail.


Typical Applications

- Measurement of the intrinsic linewidths of phonons
- Measurement of the intrinsic linewidths spin excitations
- Larmor diffraction is used to determine thermal expansion and second order stresses under pressure and at low or high temperature

Sample Environment

Besides the standard sample environment a dedicated dilution cryostat with a base temperature of 6 mK is available.



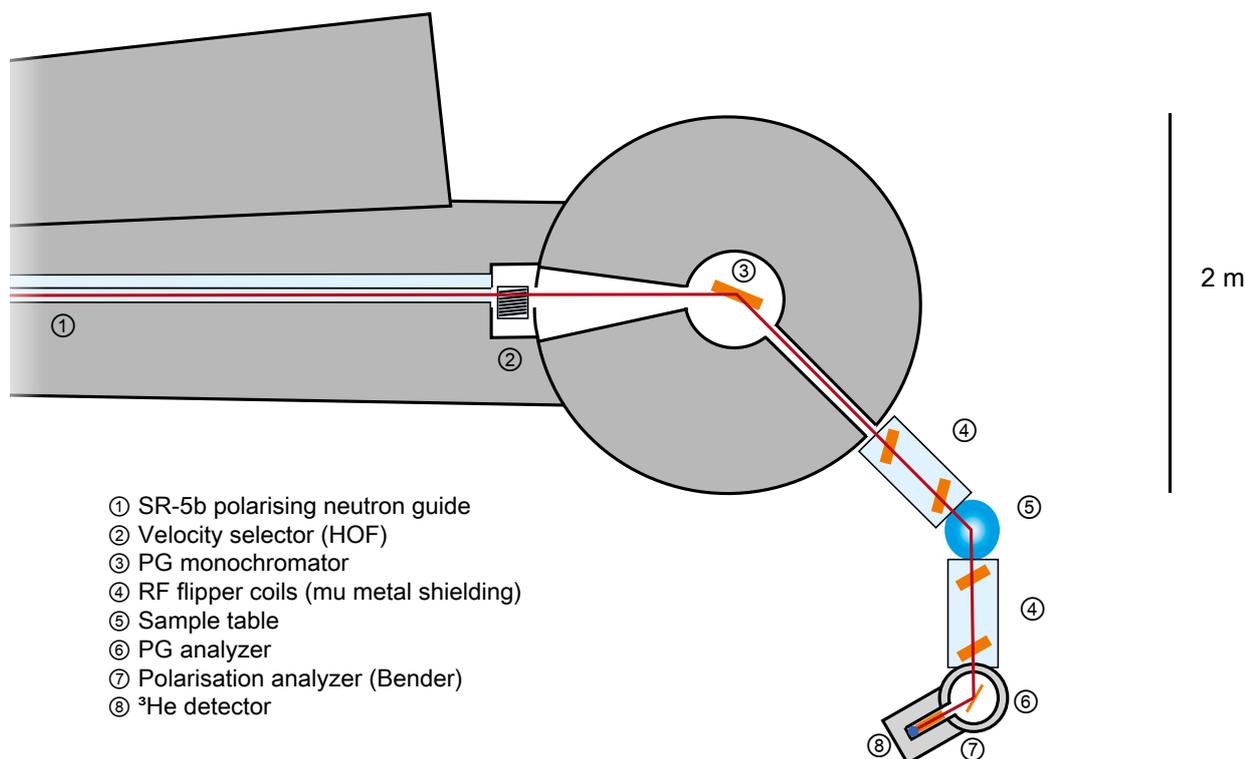

Technical Data

Primary beam

- thermal beam tube SR-5b polarising supermirror bender
 $1.3 \text{ \AA}^{-1} < k_i < 7.0 \text{ \AA}^{-1}$
- Velocity selector
 Astrium type, as higher order wavelengths filter

Monochromator

- PG(002) or (004)
 variable focussing horizontal and vertical

Analyzer

- PG(002)
 variable horizontal focussing
- Heusler (111) (polarised neutrons)
 variable horizontal focussing

Spin echo

- Resonance spin echo, enclosed by mu-metal magnetic screen.

Dr. Thomas Keller

Prof. Dr. Bernhard Keimer

Email: thomas.keller@fkf.mpg.de
 Office: +49.(0)89.289.12164
 Instrument: +49.(0)89.289.14816

b.keimer@fkf.mpg.de
 +49.(0)711.689.1650



www.mlz-garching.de/trisp
 DOI: 10.17815/jlsrf-1-41



ANTARES

cold neutron radiography and tomography station



MEDAPP

fission neutron beam for science, medicine, and industry



NECTAR

radiography and tomography station using fission neutrons



PGAA

prompt gamma and in-beam neutron activation analysis facility

Imaging & Analysis

ANTARES

cold neutron radiography and tomography station



The neutron imaging facility ANTARES is located at the cold neutron beam port SR-4a. Based on a pinhole camera principle with a variable collimator located close to the beam port, the facility provides the possibility for flexible use in high resolution and high flux imaging. ANTARES offers two different detector positions in chamber 2 and 3, which may be chosen according to the requirements for beam size, neutron flux and spatial resolution. Both chambers offer abundant space for user-provided experimental systems or sample environment.

Chamber 1 is separately accessible for the optional installation of beam and spectrum shaping devices provided by the user. At this position, ANTARES also offers built-in options such as a velocity selector, double crystal monochromator, interference gratings, and a Be-filter which are readily available for standard user operation.

Additionally we can provide access to a 300 kV microfocus X-ray CT setup for complementary investigations with a spatial resolution as good as 1 μm .

Typical Applications

The ANTARES neutron imaging facility is designed to deliver radiographs and computed tomography of samples, similar to an X-ray machine. The resulting information is often complementary to X-ray measurements with its most important features,


the high penetration of metals (Fe \sim 4 – 5 cm, Al \sim 20 – 30 cm, Pb \sim 10 – 20 cm) and the high sensitivity for hydrogen. These allow to visualise metal machine parts as well as liquids, sealants, and plastics inside of metal parts. Liquid contrast agents can be employed for crack and void detection. Examples of different techniques and their typical applications:

- **Standard neutron radiography:** Moisture in sandstone, O-rings in machine parts, aerospace pyrotechnical components, fuel cells
- **Computed tomography:** Geological samples, mineral phases, voids in carbon fiber structures (using contrast agents), machine parts, biological samples like e.g. lung tissue
- **Continuous radioscopy:** Video speed radiography of dynamic processes like boiling in refrigerators or water boilers
- **Stroboscopic imaging:** Visualisation of repetitive processes with high time resolution: Oil distribution in running combustion engines
- **Phase contrast:** Edge enhancement, aluminium foams, interfaces of similar alloys
- **Energy / wavelength scan:** Scanning for Bragg edges, phase or material identification, examination of welds
- **Polarised neutron imaging:** Metallurgical homogeneity of ferromagnetic materials, fundamental research on ferromagnetic phase transitions, visualisation of magnetic field profiles
- **Neutron Grating Interferometry:** Measurement of the spatially resolved SANS or USANS signal of the sample. Detection of microstructures on length scales of 500 nm – 10 μm , porous materials, magnetic and superconducting vortex lattice domains

Sample Environment

Standard sample environment can be used at ANTARES:

- Closed-cycle cryostats CC, CCR:
T = 50 mK – 300 K
- Electro magnet: 0 – 300 mT
- Cooling water and pressurised air

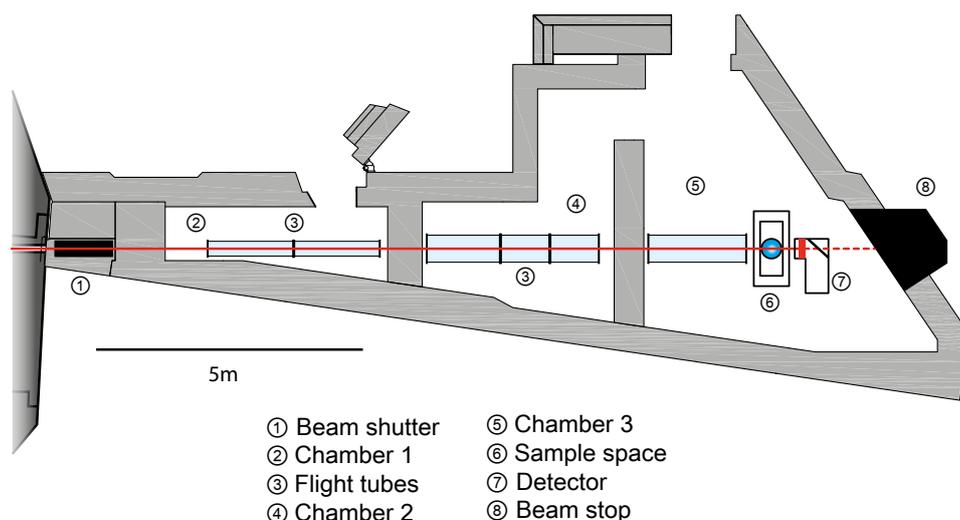

Technical Data

Collimation and flux at the sample position

- L/D = 200, $4 \cdot 10^8 \text{ n cm}^{-2} \text{ s}^{-1}$
- L/D = 400, $1 \cdot 10^8 \text{ n cm}^{-2} \text{ s}^{-1}$
- L/D = 800, $2.6 \cdot 10^7 \text{ n cm}^{-2} \text{ s}^{-1}$
- L/D = 8000, $2.6 \cdot 10^5 \text{ n cm}^{-2} \text{ s}^{-1}$
- Beam size up to $35 \times 35 \text{ cm}^2$

Neutron beam optics (optional)

- Double crystal monochromator:
 $1.4 \text{ \AA} \leq \lambda \leq 6.0 \text{ \AA}$ ($1\% < \Delta\lambda/\lambda < 3\%$)
- Neutron velocity selector:
 $3.0 \text{ \AA} \leq \lambda \leq 8 \text{ \AA}$ ($\Delta\lambda/\lambda = 10\%$)
- Neutron grating interferometer:
Sensitive to length scales 500nm – 10 μm
- Beam Filters:
Cd filter for epithermal imaging
Be filter to suppress wavelengths $\lambda < 4 \text{ \AA}$
Sapphire filter to suppress fast neutrons
- ^3He neutron spin filter polariser
- polarising supermirror V-cavity

Sample table

XY-Phi-table:

- Capacity: 500 kg
- Travel: x = 800 mm, y = 600 mm
- Rotation table: 360° rotation
- additional high precision 5-axes HUBER table for small samples (< 10 kg)

Detection systems

- various detection systems with spatial resolutions as good as 30 μm
- Camera box with mirror and scintillation screens of different sizes from $6 \times 6 \text{ cm}^2$ to $40 \times 40 \text{ cm}^2$, screen thickness from 10 μm to 200 μm , plus X-ray screens
- Standard detector: ANDOR cooled CCD camera, 2048 x 2048 pixels, 16 bit
- Fast cooled scientific CMOS camera: ANDOR Neo 2560 x 2160 pixels, 16 bit, up to 50 fps full frame
- Intensified triggerable iStar ANDOR cooled CCD camera, 1024 x 1024 pixels, 16 bit
- Intensified NTSC video camera (30 fps) with analog frame grabber, MPEG-2 and DivX recording
- DürrDental Image Plate scanner for arbitrary imaging plates, focus size 12.5 – 100 μm
- Fuji BAS 2500 Image Plate scanner, focus size 25 – 100 μm
- X-ray and neutron imaging plates
- MAR345 image plate detector, 345 mm diameter, N-sensitive image plate

Dr. Michael Schulz

Dr. Burkhard Schillinger

Email: michael.schulz@frm2.tum.de burkhard.schillinger@frm2.tum.de
 Office: +49.(0)89.289.14718 +49.(0)89.289.12185
 Instrument: +49.(0)89.289.14815




www.mlz-garching.de/antares

DOI: 10.17815/jlsrf-1-42

MEDAPP and NECTAR

fission neutron facility for scientific, medical, and industrial applications

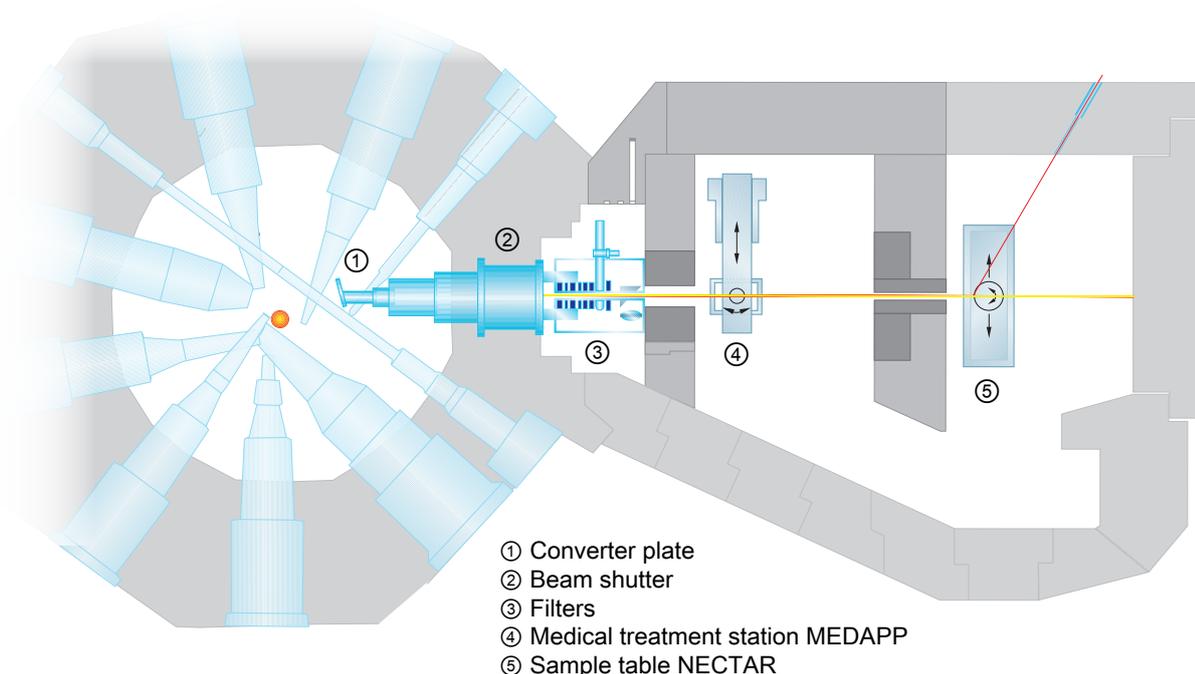

Both the instruments MEDAPP (**m**edical **a**pplications) and NECTAR (**n**eutron **c**omputed **t**omography and **r**adiography) are located at beam tube SR-10 to which a uranium converter is attached. Both instruments are operated with fission neutrons and can be used for a broad variety of different applications. For selected tasks, an alternative use with thermal neutrons is possible.

MEDAPP

MEDAPP is an instrument for the medical treatment of malignant tumours, as well as for biological research, and technical irradiations. Due to their energy spectrum, fast neutrons have the highest biological effectiveness of clinical neutron beams used in cancer treatment, comparable only to the effectiveness of heavy ions. This advantage comes at the expense of penetration depth, which - due to the relatively low energy of 2 MeV - limits their application to near-surface tumours, typically recurrent breast tumours and melanomas.

The particularly large beam cross-section of SR-10 allows the irradiation of rather large objects, such as cell culture flasks. In addition, the FaNGas (**f**ast **n**eutron **g**amma **s**pectrometry) instrument, consisting of a movable shielded HPGe detector system, can be installed within the MEDAPP irradiation chamber to directly measure gamma radiation emitted, e.g., in (n,n') , $(n,2n)$, (n,p) , and (n,α) reactions, and for non-destructive qualitative elemental analysis.

NECTAR

NECTAR is a versatile facility for the non-destructive inspection of various objects by means of fission neutron radiography and tomography, respectively. The obtained images often show complementary or additional information compared to other radiography sources such as X-rays or gamma-radiation. Especially for large objects consisting of dense materials, the deep penetration of fission neutrons allows for their non-destructive investigation maintaining the sensitivity for hydrogen containing materials. The acquired radiographs are available as *.tiff and *.fits files and can be processed with any image processing software. On demand, reconstruction and visualisation software is available for data analysis.

Christoph Genreith

Email: christoph.genreith@frm2.tum.de
Office: +49.(0)89.289.14756
Instrument: +49.(0)89.289.14831



Technische Universität München



www.mlz-garching.de/medapp-nectar

DOI: 10.17815/lsrf-1-43

Technical Data

Neutron source	
Converter facility, consisting of 2 plates of uranium-silicide. (93% ^{235}U , total 540 g), P = 80 kW	
MEDAPP	NECTAR
Neutron spectrum	
<ul style="list-style-type: none"> Fission spectrum <ul style="list-style-type: none"> Mean energy: 1.9 MeV Flux: up to $7 \cdot 10^8 \text{ cm}^{-2} \text{ s}^{-1}$ (depends on filter used) Thermal spectrum of the D_2O moderator <ul style="list-style-type: none"> Mean energy: 28 meV Flux: ca. $2 \cdot 10^9 \text{ cm}^{-2} \text{ s}^{-1}$ 	<ul style="list-style-type: none"> Fission spectrum <ul style="list-style-type: none"> Mean energy: 1.8 MeV Flux: $8.7 \cdot 10^5 \text{ cm}^{-2} \text{ s}^{-1} - 4.7 \cdot 10^7 \text{ cm}^{-2} \text{ s}^{-1}$ (depends on filter used) Thermal spectrum of the D_2O moderator <ul style="list-style-type: none"> Mean energy: 28 meV Flux: up to $1 \cdot 10^7 \text{ cm}^{-2} \text{ s}^{-1}$
Collimation	
<ul style="list-style-type: none"> Multi leaf collimator, individually adjustable 	<ul style="list-style-type: none"> L/D: $\leq 233 \pm 16$ (depending on collimator)
Sample space	
<ul style="list-style-type: none"> Max. 40 cm x 30 cm Max. about 100 kg 	<ul style="list-style-type: none"> Max. 80 cm x 80 cm, thickness dependent on material Max. 500 kg Any standard sample environment Custom environment can be easily attached (e.g. hydrogen supply)
Detection systems	
<ul style="list-style-type: none"> Ionisation chambers for dosimetry in custom-made phantoms 50%-HPGe detection system shielded with PE, B_4C, and Pb Custom systems can temporarily be installed by users 	<ul style="list-style-type: none"> CCD-based (ANDOR DV434-BV, pco. 1600) detection systems with different converters, e.g., pp-converter with 30% ZnS and 30 cm x 30 cm x 0.24 cm available
Typical applications	
<ul style="list-style-type: none"> Neutron medical treatment of malign tumours in the MEDAPP-room Biological dosimetry, e.g., irradiations of cell cultures Irradiations of electronic components (also on-line tests) Fundamental physics In-beam gamma spectrometry 	<p>Cultural Heritage</p> <ul style="list-style-type: none"> Restoration and conservation of objects Inner structure of large archaeological objects <p>Technology</p> <ul style="list-style-type: none"> Hydrogen storage Degradation of glue in timber Water or oil in large metallic objects (e.g. gearboxes) <p>Biology</p> <ul style="list-style-type: none"> Water uptake in large wooden samples

Dr. Thomas Bücherl

Dr. Stefan Söllradl

Email: thomas.buecherl@tum.de
 Office: +49.(0)89.289.14328
 Instrument: +49.(0)89.289.11777

stefan.soellradl@frm2.tum.de
 +49.(0)89.289.14768



Technische Universität München



www.mlz-garching.de/medapp-nectar

DOI: 10.17815/jlsrf-1-45

prompt gamma and in-beam neutron activation analysis facility




Prompt gamma-ray activation analysis (PGAA) is typically used for the determination of elemental composition and concentration of solid samples (ca. down to ppm range). Liquids and gaseous samples can also be measured.

The PGAA method is based on the neutron capture in nuclei of the sample material and the subsequent detection of prompt gamma-rays emitted during de-excitation of the compound nuclei: ${}^A_Z(n, \gamma){}^{A+1}_Z$.

PGAA is a non-destructive tool for the analysis of major and minor components, especially advantageous for the assay of light elements (unique for H and B) and certain trace elements (Cd, Hf, rare earths). In the strong neutron beam at FRM II, however, neutron activation can also be performed, and thus many more trace elements can be detected (elements in the 4 – 6 period).

Typical Applications

- Archaeology and cultural heritage objects (ceramics, coins, metals, conditionally bronze objects)
- Cosmochemistry (meteorites)
- Geology, petrology (macerals, sediments)
- Environmental research (air pollution, river pollution)
- Medicine (B, Li, Cd in tissues, nano-particles for cancer therapy, radiation damage of DNA)
- Semiconductor or superconductor research and industry
- Analysis of new chemical materials (catalysts, clathrates, crystals)
- Reactor physics (shielding material, new fuel element), radiation hardness testing with cold neutrons (chips, scintillators)
- Fundamental research (nuclear data, low-spin excited states in nuclei, partial and total neutron capture cross-section measurements)
- Conditionally NAA after the PGAA irradiation

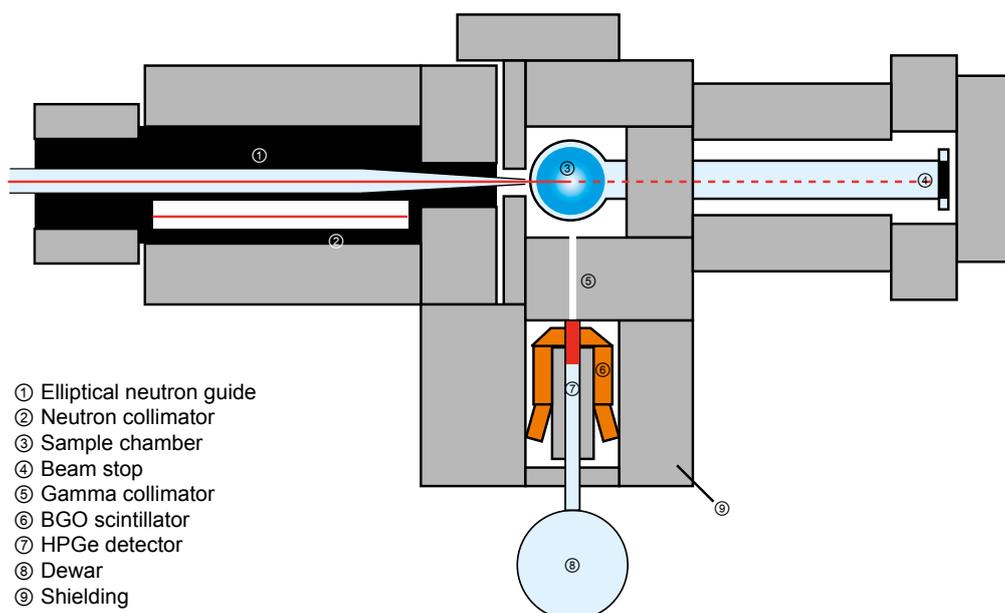

Technical Data

Neutron beam

- Cold neutron spectrum from NL4b (last section of 5.8 m elliptical focussing) with an average energy of 1.83 meV (6.7 Å)
- Two measuring conditions:
 - for large samples with collimation:
 Beam size: 20 x 30 mm²
 Neutron flux max.: $2 \cdot 10^9$ n cm⁻²s⁻¹ thermal n. eq.
 - for small samples with 1.1 m elliptical guide:
 Beam size: 11 x 16 mm²
 Neutron flux max.: $5 \cdot 10^{10}$ n cm⁻²s⁻¹ thermal n. eq.

Detection system

- For the standard PGAA, one Compton-suppressed spectrometer is used (60% HPGe detector surrounded by a BGO scintillator and connected in anticoincidence mode). The signal is processed using a DSPEC-50 digital spectrometer manufactured by Ortec.
- A new low-background counting chamber has also been installed next to the PGAA instrument for the acquisition of decay gamma spectra after activating the samples in the beam. A 30% HPGe detector is used with a DSPEC-50 unit.
- Energy range is from 30 to 12 000 keV.

Measuring conditions

- Low vacuum (0.3 mbar) possible
- Sample weight: 0.1 mg – 10 g
- Max. sample dimensions:
 ca. 40 × 40 × 40 mm³
- Automated measurement for max. six samples in a batch (vertical sample holder with six positions)
- Solid samples are usually sealed into thin FEP bags or other suitable material

Data acquisition and analysis

- An in-house software for the automated measurement of up to six samples in a batch run.
- Evaluation of the spectra and the calibration of the spectrometer (efficiency curve and non-linearity) using the software Hypermet PC developed in Budapest
- Determination of the elemental composition of samples using the Excel macro and Excel sheet package ProSpeRo
- Automated data acquisition using DSPEC-50 is currently under development

Dr. Zsolt Revay

Email: zsolt.revay@frm2.tum.de
 Office: +49.(0)89.289.12694
 Instrument: +49.(0)89.289.14906



www.mlz-garching.de/pgaa
 DOI: 10.17815/jlsrf-1-46




EDM
neutron electric dipole moment measurement



MEPHISTO
facility for particle physics with cold neutrons

Particle Physics

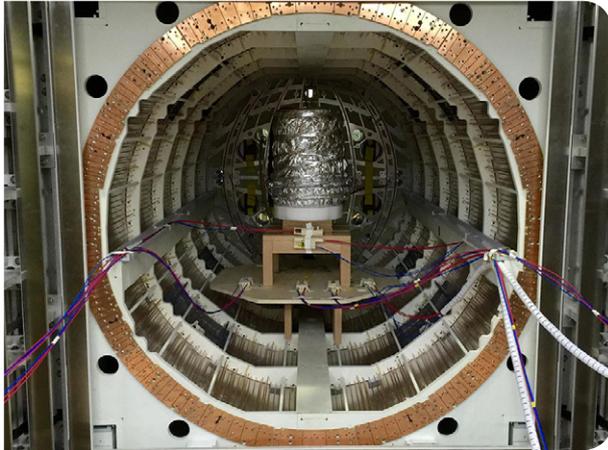

An **electric dipole moment (EDM)** of the neutron would be a clear sign of new physics beyond known particle physics. The search for this phenomenon is considered one of the most important experiments in fundamental physics and could provide key information on the excess of matter versus antimatter in the Universe. With high measurement precision, this experiment aims to ultimately achieve a sensitivity of 10^{-28} ecm for a charge distribution within the neutron. This can be interpreted as a probe of the early Universe, less than 10^{-11} s after the big bang, complementary to experiments at CERN and also beyond the reach of accelerators. In the Neutron Guide Hall East, a new instrument to search for the neutron EDM is currently being commissioned.

At a later stage, **ultra cold neutrons (UCN)** will be provided by a source in beam tube SR-6. In the instrument's centre, the UCN are trapped in a box, where Ramsey's method of separated oscillatory fields is used to precisely measure Larmor precession of UCN over hundreds of seconds as an interferometer in time in a small (about $1 \mu\text{T}$) magnetic field. A deviation of the Larmor precession in the simultaneous presence of the magnetic field and a strong electrostatic field would be an indication for an EDM. While the UCN source will provide the required neutron density of 1000 per cm^3 in the instrument, it itself is optimised to control any systematic effects.

A dominating issue to achieve this precision is the quality of the magnetic environment. Currently, it

provides an extended volume ($\sim 0.5 \text{ m}^3$) of 100 pT or less absolute field and a temporal stability of $< 10 \text{ fT}$ in 1000 s without any active measures. This corresponds to about the magnetic field in between stars in the milky way.

The shielding factor against external distortions of 6 Mio at 1 MHz and $> 10^8$ at 10 Hz. Inside this volume, magnetic fields of up to few μT can be applied vertically and horizontally with a homogeneity of better than up to 10^{-4} relative precision. During the commissioning phase, the volume is used for optical magnetometry and precision spin-clock measurements, as well as for the comparison of precision magnetometers. To apply electrostatic fields, a bipolar 200 kV source with stabilisation on Pockels effect is implemented. The instrument is a joint activity of several institutions, see nedm.ph.tum.de


Other optional measurements inside the instrument include:

- Ramsey-measurements with trapped ultra cold neutrons
- Probing the isotropy of the Universe and its invariance under Lorentz transformations via clock comparison measurements with polarised gases
- Measuring new forces at short distances: search for dark matter candidates (axions etc.)
- Precision magnetometry and magnetometer development
- Ultra low field NMR
- Magnetic material screening on pT level

Technical data

The EDM facility consists of a magnetically shielded room, surrounded by 180 fluxgate magnetometers and 24 coils for external field compensation. The magnetically shielded room consists of an outer and an inner room, which can be used independently or slided into each other. Experiments with $2 \times 2 \times 2 \text{ m}$ dimension can be placed inside the shields in pre-assembled state through a door of $2 \times 2 \text{ m}$. The furniture inside is non-magnetic and either 3D-printed or made from wood, external equipment can be degaussed on-site. Stable current sources with ppm stability over $> 1000 \text{ s}$ can provide DC magnetic fields with up to $5 \mu\text{T}$.

Surrounding field compensation

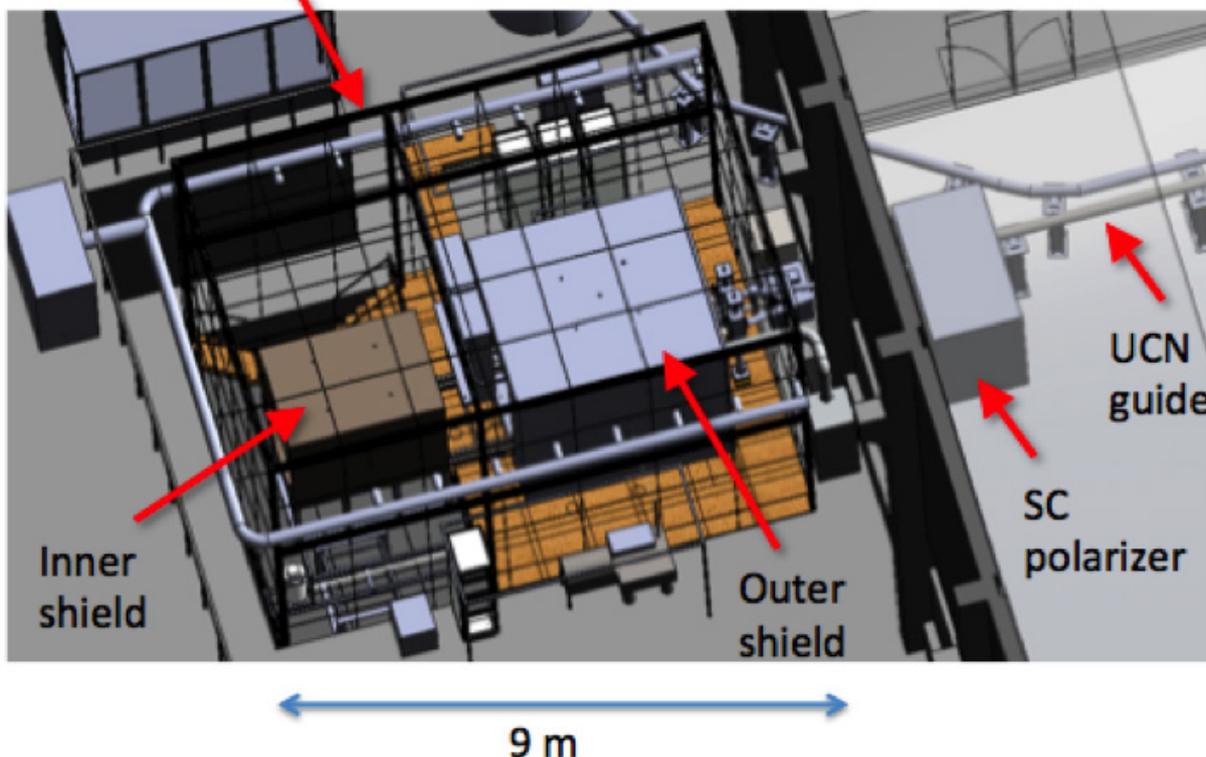

External compensation

- Size 9 x 6 x 6 m
- Number of coils: 24
- Number of magnetometers: 180
- Field < 5 μT (DC)
- Active feedback frequency: 10 Hz (max.)

Outer shield

- Inner dimension:
h = 2.3 m, w = 2.5 m, l = 2.8 m
- Residual field: < 0.5 nT in 1 m³
- Damping factor:
279 at 0.01 Hz, > 10000 at 10 Hz

Inner shield

- Inner dimension:
h = 1.5 m, w = 1.5 m, l = 2.1 m
- Residual field: < 0.1 nT in 0.5 m³
- Damping factor:
6.000.000 at 0.001 Hz, > 10⁸ at 10 Hz

Available magnetometry

- Fully optical Cs magnetometers (< 30 fT in 40 s resolution)
- ¹⁹⁹Hg magnetometers (< 20 fT in 100 s)
- SQUID magnetometers and gradiometers
- SEOP polarisation of ³He and ¹²⁹Xe at the experiment

Prof. Dr. Peter Fierlinger


Email: peter.fierlinger@tum.de
Office: +49.(0)89.35831.7131



www.mlz-garching.de/edm
DOI: 10.17815/jlsrf-1-47

MEPHISTO

facility for particle physics with cold neutrons

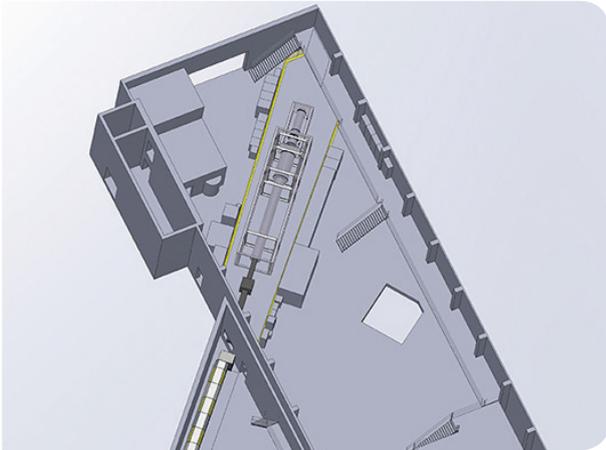

Since the start of reactor the FRM II provides a cold white neutron beam for long term user dedicated experimental setups.

Such an experiment is normally planned and built up by external groups but additional help during the commissioning of the experiment at the reactor is necessary. Therefore, this work must be organised in close contact with the local instrument scientist. The desired precision is reached inter alia by good statistics which means long term experiments over several reactor cycles.

The experimental area MEPHISTO, the measurement facility for particle physics with cold neutrons, is dedicated to those experiments in the field of nuclear and particle physics. Currently, the experimental area moves from the Neutron Guide Hall West to the Neutron Guide Hall East. The solely used neutron guide SR-4b will deliver a white cold spectrum for experiments. A removeable 11% velocity selector at the end of the guide will complete the beam line.

The MC-simulation for this beam with a dimension of $60 \times 106 \text{ mm}^2$ proposes a mean wavelength of 4.5 \AA and a gold capture flux of $2 \cdot 10^{10} \text{ n cm}^{-1}\text{s}^{-1}$. The experimental area is $5 \times 25 \text{ m}^2$, diagonally built-in in the Neutron Guide Hall East. The spectrum shows a shoulder to smaller wavelengths, the maximum of the spectrum is located at 3.3 \AA .

It is planned to install the instrument PERC [1] at the MEPHISTO beam line during the first years of operation in the Neutron Guide Hall East. This instrument is a precise, bright and intense source of protons and electrons from the neutron decay. The instrument PERC itself is open for external user groups with spectrometers to measure the protons and electrons.

Typical Applications

The experiments at MEPHISTO concentrate on induced nuclear reactions of the neutron with atoms or on the free neutron decay with its products.

Some of the experiment types performed at MEPHISTO:


- Free neutron decay and spectroscopy of the decay products
- Spectroscopy of neutron induced fission
- Production of ultra cold neutrons with liquid helium
- Production of ultra cold neutrons with solid gases

Infrastructure

A removeable neutron velocity selector is placed at the end of the neutron guide. The minimal wavelength is 4.5 \AA . The resolution of the passing wavelength is 11%. The selector can be rotated to tune the resolution.

A data system based on VME (ADC, peak ADC, QDC, TDC) is available. For signal forming purpose several NIM inserts exist, a list can be requested from the local instrument scientist. Also available are spectroscopic amplifiers and high voltage sources for detectors.

[1] Dubbers, D. et al., Nucl.Instr.Meth. A, 596, 238-247, (2008).

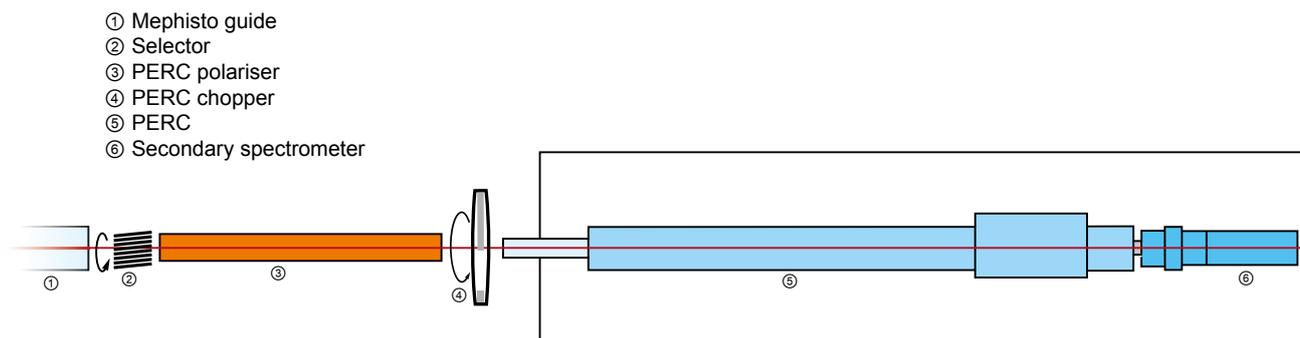

Technical Data

Neutron beam

- End of the cold neutron guide SR-4b (m = 2.5)
- Cross section of the guide: 60 x 106 mm²
- Thermal capture flux (simulated):
 $2 \cdot 10^{10} \text{ n cm}^{-1}\text{s}^{-1}$
- Mean wavelength (simulated): 4.5 Å
- Beam height from floor: ~ 1300 mm
- Experimental area: 5 x 25 m²
- Maximum at 3.3 Å
- Standard neutron spectrum with shoulder to smaller wavelengths

Beam attenuators

- By geometrical attenuation, the beam intensity can be reduced to 20%, 4% and 2%

Polarisation

- A bender (vertical direction) is available to polarise the complete cross section of the beam

Dr. Jens Klenke

Email: jens.klenke@frm2.tum.de
 Office: +49.(0)89.289.14771



www.mlz-garching.de/mephisto
 DOI: 10.17815/jlsrf-1-48



NEPOMUC
neutron induced positron source Munich



CDBS
coincident Doppler-broadening spectrometer



PAES
positron annihilation induced Auger-electron spectrometer



PLEPS
pulsed low energy positron system



SPM
scanning positron microscope

Positrons

NEPOMUC

neutron induced positron source Munich

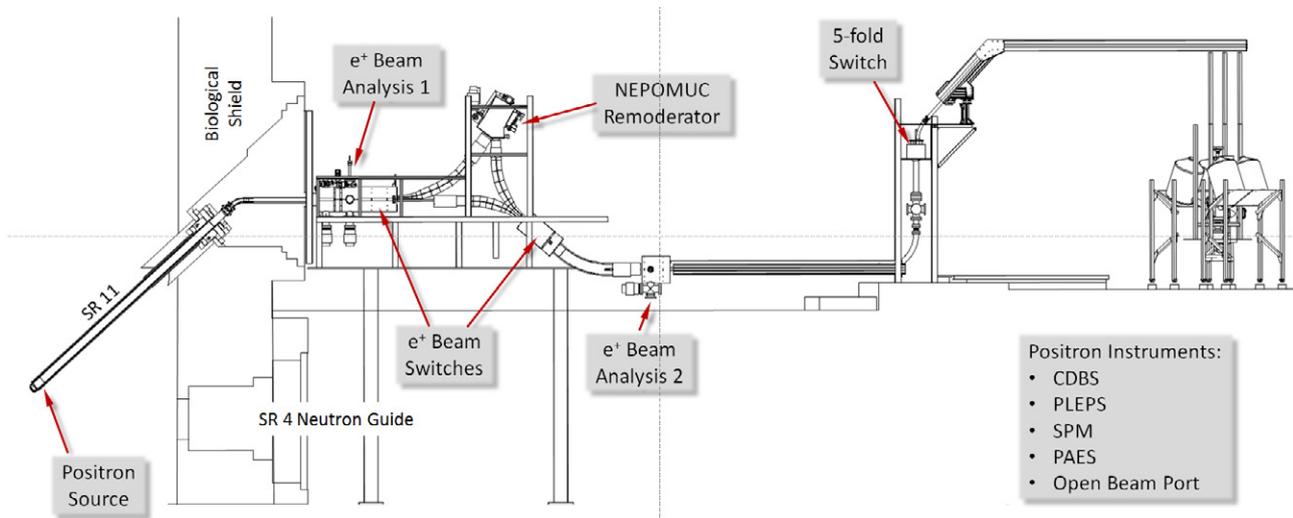


Figure 1: Cross-sectional view of the inclined beam tube SR11: the in-pile positron source is mounted inside the tip. After acceleration, the positron beam is magnetically guided to the remoderation unit outside the biological shield of the reactor.

NEPOMUC provides a high-intensity low-energy positron beam for applications in solid state and surface physics as well as for fundamental research in nuclear and atomic physics.

At NEPOMUC, the positrons are generated by pair production from absorption of high-energy prompt gamma-rays after thermal neutron capture in Cd. A cadmium cap is mounted inside the tip of the inclined beamtube SR11. The released high-energy gamma-radiation is converted into positron – electron pairs in a structure of platinum foils which is mounted inside the cadmium cap. Positive high voltage is applied in order to extract the moderated positrons. The positron beam is magnetically guided in a solenoid field of typically 7 mT.


Technical Data

Key values of the primary positron beam

- E = 1 keV
- Intensity:
> $1.1 \cdot 10^9$ e⁺/s moderated positrons per second
- Diameter of beam spot:
≈ 9 mm (FWHM) in 7 mT beam guiding field.

Key values of the remoderated positron beam

- E = 10 ... 200 eV
- Intensity:
 $5 \cdot 10^7$ remoderated positrons per second
- Diameter of beam spot:
1.85 mm (FWHM) in ≈ 4 mT beam guiding field

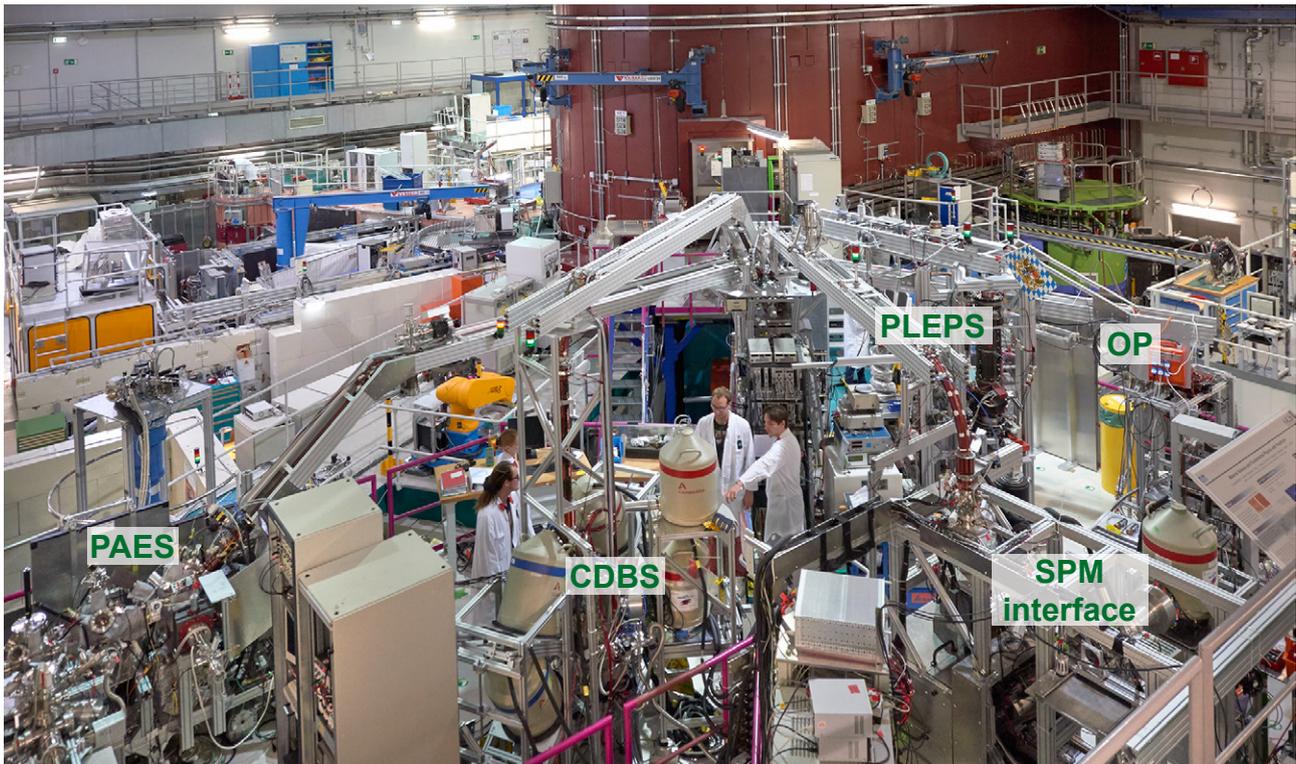

The positron beam facility

The remoderation device of NEPOMUC enhances the brightness of the positron beam and hence enables positron experiments which are high resolved in space or/ and in the time domain. The remoderator is based on the stochastic cooling of the positrons in a W(110) single crystal and the positron reemission of the thermalised positrons into the vacuum with discrete energy. For most of the measurements the brightness enhanced positron beam is used. However, there are also experiments which are not dependent on a very bright beam but need the full intensity of the primary beam. Therefore the primary beam can be used also unaltered via two beam switches e.g. for experiments at the open beam port.

By the central, fivefold beam switch the positron beam is delivered to one of the five experiment beam lines. Currently, at these beam lines four instruments are permanently installed:

- Pulsed low-energy positron system (PLEPS)
- Coincident Doppler-broadening spectrometer (CDBS)
- Positron annihilation induced Auger-electron spectrometer (PAES)
- Scanning Positron Microscope Interface (SPM Interface)

The fifth beam line is for the multi-purpose open beam port (OP) which is used for transportable short-term experimental set-ups.

Dr. Christoph Hugenschmidt


Email: christoph.hugenschmidt@frm2.tum.de
Office: +49.(0)89.289.14609
Instrument: +49.(0)89.289.14774

Dr. Christian Piochacz

christian.piochacz@frm2.tum.de
+49.(0)89.289.12179



www.mlz-garching.de/nepomuc
DOI: 10.17815/jlsrf-1-49

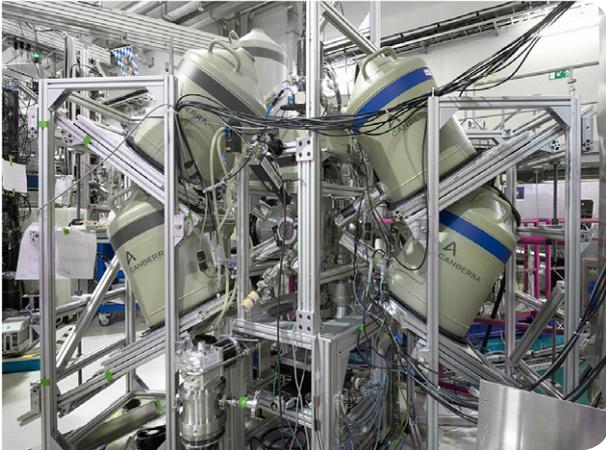

The Doppler broadening of the 511 keV annihilation line contains information of the electron momentum distribution at the positron annihilation site in the sample. Since the probability of core electron annihilation decreases in open volume defects a narrowing of the annihilation line is observed.

For this reason, DBS with the monoenergetic positron beam allows to determine defect profiles, energy dependent 2D imaging of defects, and defect annealing as a function of temperature. In addition CDBS is applied in order to gain elemental information about the positron annihilation site and hence about the chemical surrounding of defects.

Technical Data

Beam properties

- Positron implantation energy: $E = 0.2 - 30$ keV
- Mean positron implantation depth: up to several μm (material dependent)
- Beam size: adjustable between $0.3 - 3$ mm \varnothing

2D x-y-scans

- Scan area: 20×20 mm²
- Step size adjustable between 0.1 and 10 mm

High-purity Ge detectors

- 30% efficiency
- Energy resolution: 1.4 keV at 477.6 keV

Sample

- Size up to $20 \times 20 \times 3$ mm³
- Optimum 4 samples at one sample holder: $< 10 \times 10$ mm²
- Temperature: 100 K – 900 K

Typical measurement times

- ~ 1 – 2 min / spectrum
- ~ 8 h full 2D overview scan (with $\Delta x = \Delta y = 1$ mm)
- ~1 h depth profile ($t = 2$ min, 30 energy values)
- ~ 4 – 6 h/spectrum CDBS

Dr. Christoph Hugenschmidt

Email: christoph.hugenschmidt@frm2.tum.de
 Office: +49.(0)89.289.14609
 Instrument: +49.(0)89.289.14774



Technische Universität München




www.mlz-garching.de/cdbS

DOI: 10.17815/Jsrf-1-50

PAES

positron annihilation induced Auger-electron spectrometer

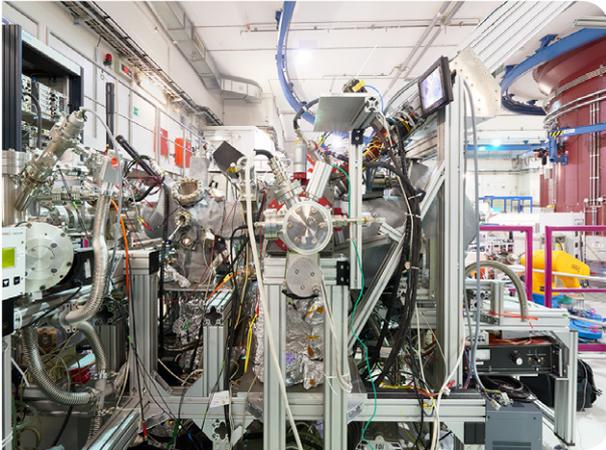

PAES is a newly developed application for surface studies with high elemental selectivity and exceptional surface sensitivity. In PAES, the emission of Auger electrons is initiated by positron-electron annihilation that leads to several major advantages, e.g. topmost layer sensitivity, compared with conventional electron induced AES.

PAES is part of the **surface spectrometer** (SuSpect) which also enables sample preparation in UHV conditions, conventional AES and XPS.

Examples are surfaces with sub-monolayers of foreign atoms, high resolution determination of Auger line shapes, element selective surface studies.

Technical Data

Beam properties

- Positron implantation energy: 20 eV
- Electron energy resolution: $\Delta E/E < 1\%$

Sample

- Sample size: \varnothing 10 mm
- Sample thickness: max. 3 mm

Typical measurement times

- Measurement time (typically): 10 – 15 min

Complementary techniques

- Electron or X-ray induced AES
- XPS
- STM

Dr. Christoph Hugenschmidt


Email: christoph.hugenschmidt@frm2.tum.de
Office: +49.(0)89.289.14609
Instrument: +49.(0)89.289.14774



Technische Universität München



www.mlz-garching.de/paes
DOI: 10.17815/jlsrf-1-51

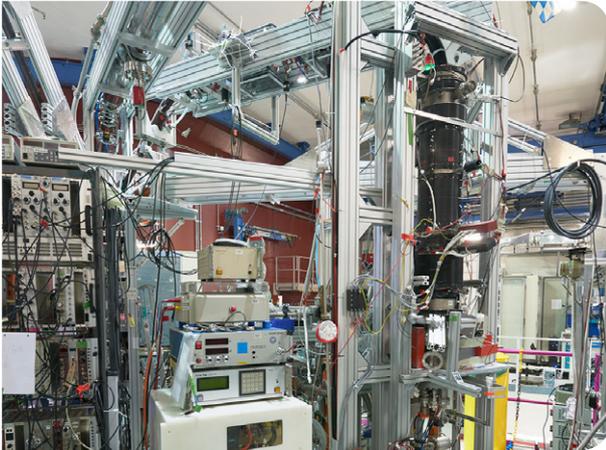

PLEPS is a unique tool for depth profiling of defects with positron annihilation lifetime spectroscopy using a pulsed positron beam of variable energy.

Positron lifetime measurements allow to determine type and size of open volume defects (such as vacancies, vacancy-clusters, dislocations, grain boundaries etc., and free volumes in polymers) in a wide variety of materials and provide information on defect-concentration. In combination with a mono-energetic positron beam of variable energy depth-resolved defect analysis becomes possible.

Typical applications

- Defect identification in thin layers and layered structures of semiconductors and insulators
- Radiation induced defects in materials for fusion and fission reactors
- Characterisation of free volumes in polymers and glasses

Technical Data

Beam properties

- Positron implantation energy: 0.5 – 20 keV
- Beam spot $\varnothing \sim 1$ mm
- Count rate: $\sim 5000 - 10000$ cps

Sample

- Limited to 5×5 mm² – 9×9 mm²

Typical measurement times

- < 10 min per spectrum
($> 3 \cdot 10^6$ counts in the spectrum)
- Depth-profile: 4 – 5 h
(15 – 20 implantation energies, $> 3 \cdot 10^6$ counts in the spectrum)
- Time-window: 20 ns or 40 ns
- Time-resolution: 260 – 280 ps
- Peak/ background $> 50000 : 1$

Dr. Werner Egger

Email: werner.egger@unibw.de
Office TUM: +49.(0)89.289.14609
Office UniBW: +49.(0)89.6004.3007

der Bundeswehr
Universität  München



www.mlz-garching.de/pleps

DOI: 10.17815/ilsrf-1-52



The Munich scanning positron microscope permits positron lifetime measurements with a lateral resolution in the μm range and within an energy range of 1 – 20 keV. Thus, this instrument enables the measurement of high resolved 3D defect maps. Until today, the SPM was operated only in the laboratory at the Universität der Bundeswehr in Munich and was therefore limited by the long measurement times of several days per 2D-scan due to the low intensity of the positron beam produced by a standard ^{22}Na source. This disadvantage will be overcome by installing the SPM at the high intensity positron beam at NEPOMUC.

Therefore, the SPM interface was designed and tested successfully. This device converts the continuous beam of NEPOMUC to a high-brightness, pulsed positron beam, which matches the demands of the SPM. Recently, a sample chamber was connected to the SPM interface which enables spatially resolved positron lifetime measurements with a lateral resolution in the range of 0.1 mm.

Technical Data

Beam properties SPM / SPM Interface

- Positron implantation energy: < 20 keV / < 10 keV
- Beam-Spot < 1 μm / \approx 0.1 mm
- Count rate: > 2000 cps / > 4000 cps
- Time-Window: 20 ns
- Time-Resolution: < 250 ps
- Peak/ Background: > 5000 : 1 / > 2000 : 1

Typical measurement times

- SPM: \approx 1 day for one 2D-Scan (12 x 12 μm^2)
- SPM interface: \approx 0.5 day for one 2D-Scan (1 x 1 mm^2)

M.Sc. Marcel Dickmann

Email: marcel.dickmann@unibw.de
Office TUM: +49(0)89.289.11770
Office UniBW: +49(0)89.6004.3508

Universität  München



Technische Universität München



www.mlz-garching.de/spm
DOI: 10.17815/jlsrf-1-53



Silicon doping facility (SDA)



Standard rabbit irradiation facility (RPA)




Capsule irradiation facility (KBA)



Mechanical irradiation system



Irradiation position in control rod

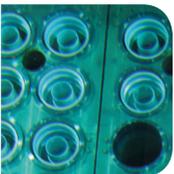

Gamma-ray irradiation facility



Irradiation possibilities at beam tube SR-10 using fission neutrons

Irradiation facilities

Irradiation facilities

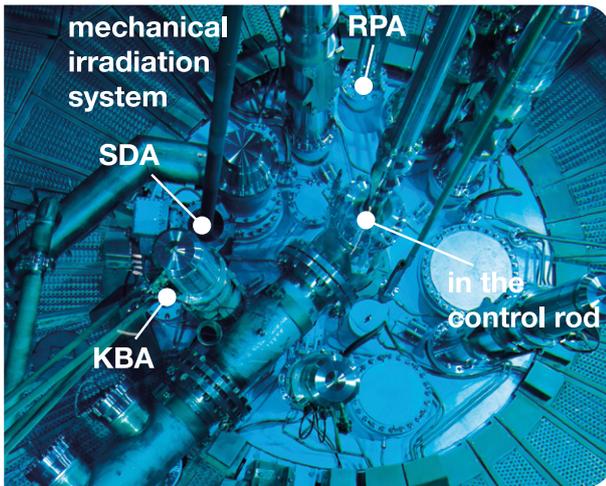


Figure 1: Positions of the different irradiation facilities within the reactor pool: SDA = silicon doping facility, RPA = rabbit irradiation system, KBA = capsule irradiation facility, mechanical irradiation system, and irradiation position in the control rod.



Figure 2: A freshly irradiated Si batch is taken out of the reactor pool.

The irradiation of materials serves many purposes including the doping of silicon and the production of radioisotopes for industry and medicine, as well as the examination of samples with neutron activation analysis.

The neutron source FRM II is equipped with a number of irradiation facilities, which cover a wide range of applications both with regard to the available sample volumes and the achievable neutron fluences.

For proposals and experiments dealing with the irradiation facilities please contact directly the head of the irradiation services Dr. Heiko Gerstenberg.

Silicon doping facility (SDA)

Pure silicon is a poor conductor of electricity. In order to gain the properties which make it interesting for components in the electrical engineering, it needs to be doped with small amounts of host atoms e.g. phosphorous or boron. For high-performance electronic components as thyristors or insulated-gate bipolar transistors (IGBTs) the silicon (Si) needs to have a defined content of phosphorous (P) atoms distributed extremely homogeneously within the Si matrix. At the FRM II, silicon ingots up to a diameter of 200 mm and a height up to 500 mm are irradiated in a position within the moderator tank (see fig. 2). The doping is achieved by neutron capture and the

resulting conversion of individual ^{30}Si atoms into ^{31}P . Due to the neutron moderation by heavy water the facility is particularly useful for the production of high resistivity (up to 1200 Ωcm) neutron transmutation doped (NTD) Si. The silicon doping facility is automated and operated in two working shifts. The typical yearly output sums up to about 15 tons. The customers of the doped ingots are semiconductor producers from Europe and Asia.

Standard rabbit irradiation system (RPA)

Six independent irradiation channels are available within the standard rabbit irradiation system. The positions are vertically staggered in the moderator tank, allowing for the selection of a thermal neutron flux density adapted to the sample. The neutron flux densities range from $5 \cdot 10^{12}$ to $7 \cdot 10^{13} \text{ cm}^{-2}\text{s}^{-1}$. Sample sizes should be less than 8 cm^3 . The samples are packed in polyethylene capsules and conveyed pneumatically by CO_2 into the irradiation position. The available thermal neutron fluence varies between $2 \cdot 10^{14} \text{ cm}^{-2}$ and $3 \cdot 10^{17} \text{ cm}^{-2}$. For the various irradiation channels the ratio of thermal/ fast neutron flux density is as high as 15000 – 60000.

These favourable irradiation conditions avoid the production of undesired radioisotopes due to threshold reactions by fast neutrons and make the pneumatic rabbit system a valuable tool to be used



Figure 3: The standard rabbit irradiation system has six independent irradiation channels and is operated by carbon-dioxide.

in the framework of neutron activation analysis or for the production of particular radioisotopes with well defined low activity values, as required for instance in material sciences. In order to facilitate the cooperation with the neighbouring Institute for Radiochemistry (RCM) the pneumatic rabbit system is connected directly to their laboratories by an additional pneumatic dispatch.

Capsule irradiation facility (KBA)

High dose irradiation spanning periods from several hours to weeks, is carried out in the capsule irradiation facility. It is a pool water-driven hydraulic rabbit system with two mainly identical irradiation channels exhibiting thermal neutron flux densities of up to $1.3 \cdot 10^{14} \text{ cm}^{-2}\text{s}^{-1}$. The samples are packed in aluminium capsules with a volume of up to 20 cm^3 . If required, the samples are additionally packed water tight in an inner capsule made from high purity aluminium or quartz in order to avoid their direct contact to pool water. At maximum three capsules may be irradiated simultaneously in each channel. The ratio of thermal/ fast neutron flux density ranges between 330 and 770. The main application of KBA is the production of high activity radioisotopes to be used in radiopharmaceuticals or for industrial purposes such as density or level sensors. Further activities deal with the study of the deterioration of moving parts in motors.

inside moderator tank

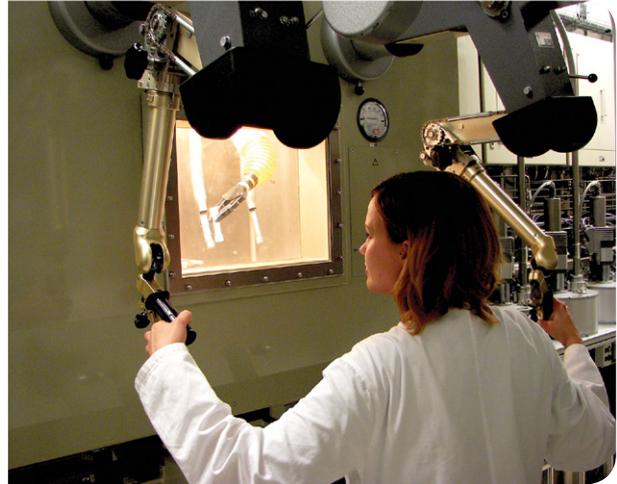


Figure 4: After irradiation in the rabbit system, the samples are prepared via manipulators to be sent to the customers.

Mechanical irradiation system

Samples sized up to 120 cm^3 can be irradiated in a facility in the moderator tank. In spite of the available sample size of 2.5 l, the facility is typically used for short term irradiations of smaller samples. The thermal neutron flux density is $1.1 \cdot 10^{13} \text{ cm}^{-2} \text{ s}^{-1}$. The maximum licensed irradiation time is 2 h corresponding to a thermal neutron fluence of $8 \cdot 10^{16} \text{ cm}^{-2}$. Due to the low radioactivity resulting from short term irradiations at moderate neutron flux density the facility is particularly useful for scientists who are less interested in radioactivity but in other irradiation induced effects on sample properties. The mechanical irradiation system turned out to be perfectly suitable for geochronological studies using the fission track technique.

Irradiation position in control rod

The highest possible fluence of $1.1 \cdot 10^{21} \text{ cm}^{-2}$ can be reached at the irradiation position within the control rod. It is only possible, however, to load and unload it when the reactor is shut down, after the completion of its 60-day cycle.

Dr. Heiko Gerstenberg

Email: heiko.gerstenberg@frm2.tum.de
Office: +49.(0)89.289.12145



Irradiation facilities

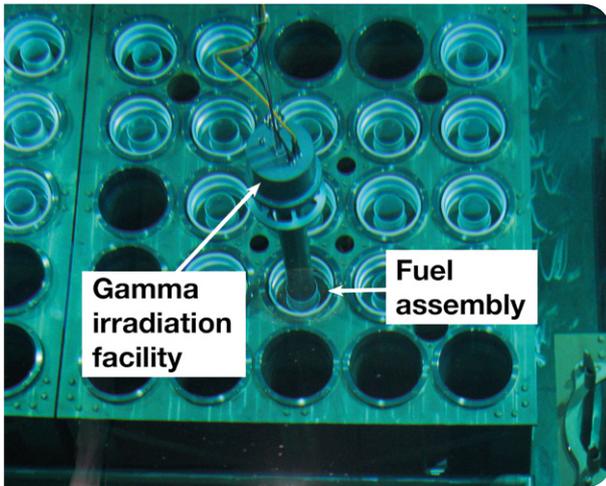


Figure 5: Gamma irradiation facility in spent fuel assembly.

Gamma-ray irradiation facility

The newest facility in operation at the FRM II is gamma irradiation. It uses the intense gamma-radiation emitted by spent fuel assemblies (fig. 5). Its installation was supported by the German Research Association (Deutsche Forschungsgemeinschaft, DFG). It is located in the storage rack of spent fuel elements. The typical gamma dose rate ranges from 1 kGy/h up to 100 kGy/h. The sample container has a diameter of 76 mm and a height of approximately 1 m. Additionally, the sample position can be heated up to 140°C, if required. The irradiation time can vary between several minutes to several weeks. Up to now the application is restricted to inorganic solids. If required the application range could be extended to organic compounds, e.g. cable insulation or organic sealing materials.



Figure 6: ^{177}Lu produced by the company ITG at the FRM II at the capsule irradiation facility. The radioisotope is filled in little vials with a volume of 2 ml and then packed in cans with a diameter of 45 mm and a height of 70 mm.



Figure 7: Loading of the capsule irradiation facility.

Typical applications

Neutron Activation Analysis

Neutron activation analysis (NAA) is used to analyse element composition in a material. Up to 30 or 40 elements can be determined simultaneously down to the ppt and sub-ppt range.

Industrial applications:

- Trace elements in pure silicon
- Environmental monitoring, e.g. retained substrates in the filters of an exhaust of a chemical production

Archeological and geological applications:

- Fingerprint of materials gives clues about the origin of the findings
- Determination of the age and composition of rock by methods requiring reactor irradiation

Radioisotopes for medical and technical applications

Production of radioisotopes in the different facilities (e.g. KBA, RPA)

- ^{60}Co for industrial purposes (e.g. non-destructive material testing)
- ^{177}Lu , which is used for the therapy of neuroendocrine tumours (see fig. 6)
- ^{161}Tb , ^{166}Ho for medical purposes (primarily for tumour therapy)
- The production of ^{99}Mo at the FRM II is planned to start in 2018. Its daughter isotope $^{99\text{m}}\text{Tc}$ is used in more than 80% of all nuclear medical diagnoses.



Figure 8: The tumour irradiation facility (MEDAPP) with a face mask for the accurate positioning of the patient. The Klinikum rechts der Isar (MRI), Department of Radiotherapy and Oncology, is responsible for conducting the therapy.

Irradiation possibilities at beam tube SR-10 using fission neutrons

The beam tube SR-10 can deliver very different beam qualities: A pair of uranium converter plates generates an unmoderated fission spectrum with mean energy of 1.9 MeV, accompanied by prompt and delayed photons. A large, homogeneous, horizontal beam of either fission neutrons (flux up to $7 \cdot 10^8 \text{ cm}^{-2} \text{ s}^{-1}$) or thermal neutrons (flux up to $3.9 \cdot 10^9 \text{ cm}^{-2} \text{ s}^{-1}$) is available for scientific, medical, and commercial testing, available with an adjustable exposure area of up to $40 \times 30 \text{ cm}^2$. More details about the irradiation possibilities at SR-10 can be found in the instrument section about the instruments MEDAPP and NECTAR.

outside moderator tank

Technical Data (all facilities)

Silicon doping facility (SDA)

- $1.7 \cdot 10^{13}$ thermal neutrons $\text{cm}^{-2} \text{ s}^{-1}$
- Si dimensions: height $\leq 500 \text{ mm}$, $\varnothing = 200, 150, 125 \text{ mm}$
- ρ_{target} : $25 \Omega \text{ cm} - 1200 \Omega \text{ cm}$

Standard rabbit irradiation system (RPA)

- $5 \cdot 10^{12}$ to $7 \cdot 10^{13}$ thermal neutrons $\text{cm}^{-2} \text{ s}^{-1}$
- $\phi_{\text{th}} / \phi_{\text{f}}$: $15000 - 60000$

Capsule irradiation facility (KBA)

- Up to $1.3 \cdot 10^{14}$ thermal neutrons $\text{cm}^{-2} \text{ s}^{-1}$
- $\phi_{\text{th}} / \phi_{\text{f}}$: $330 - 770$

Mechanical irradiation system

- $1.1 \cdot 10^{13}$ thermal neutrons $\text{cm}^{-2} \text{ s}^{-1}$
- $\phi_{\text{th}} / \phi_{\text{f}}$: ~ 770

Irradiation position in the control rod

- $2 \cdot 10^{14}$ thermal neutrons $\text{cm}^{-2} \text{ s}^{-1}$
- $2 \cdot 10^{14}$ fast neutrons $\text{cm}^{-2} \text{ s}^{-1}$

Gamma ray irradiation facility

- 1 kGy/h up to 100 kGy/h

Neutron activation analysis

- Irradiation time: seconds – hours
- Sample weight: mg – g
- Rel. efficiency at 1.3 MeV: $> 40\%$
- Energy resolution: $< 0.9 \text{ keV}$ at 122 keV , $< 1.8 \text{ keV}$ at 1.3 MeV

Irradiation with fast neutrons at SR10

- up to $7 \cdot 10^8$ fast neutrons $\text{cm}^{-2} \text{ s}^{-1}$

Dr. Heiko Gerstenberg

Email: heiko.gerstenberg@frm2.tum.de
Office: +49.(0)89.289.12145





SE
magnetic field



LAB
Bio Laboratory



SE
high pressure



LAB
Chemistry Laboratory



SE
specialised equipment



LAB
Materials Science Laboratory



SE
low temperatures



LAB
Sample Preparation Laboratory



SE
high temperatures



UF
TEM




UF
MBE



UF
X-ray Tomography

Sample Environment (SE), Laboratories (LAB), and User Facilities (UF)

MAG-2.2T-HTS

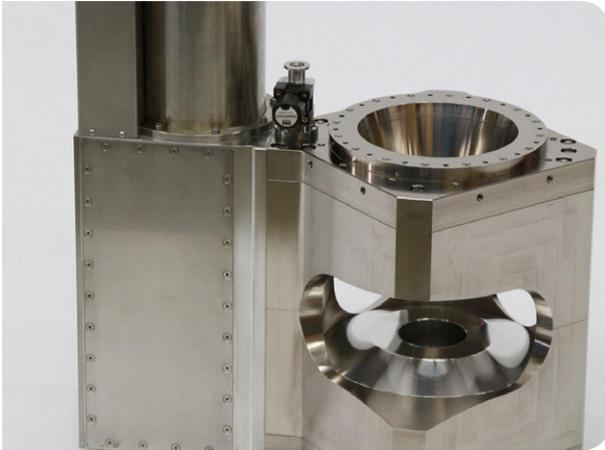

Specifications

- Maximum magnetic field:
- ± 2.2 T, orientation arbitrary
- Homogeneity of the magnetic field (15mm DSV): 1.45%
- Three room-temperature bores: minimum \varnothing 80 mm
- Beam window dimensions: 4 conical windows, opening angle $\pm 20^\circ$, 1 window opening angle 150°
- Cooling system: dry
- Cool down: 22 hours
- Additional sample environment available (CC, CCR)

MAG-H-5.0T



Specifications

- Maximum magnetic field: ± 5.0 T (symm.), ± 2.5 T (asymm.)
- Homogeneity of the magnetic field (25 mm x 25 mm cylinder): 2.0%
- Three room-temperature bores: \varnothing 80 mm
- Beam window dimensions: horizontal and vertical angle: 30°
- Cooling system: dry
- Cool down: 10 days
- Additional sample environment available (CCR)

Dr. Jürgen Peters

Email: juergen.peters@frm2.tum.de
Office: +49.(0)89.289.14700


www.mlz-garching.de/se

Dr. Jürgen Peters

juergen.peters@frm2.tum.de
+49.(0)89.289.14700


www.mlz-garching.de/se

MAG-V-7.5T



Specifications

- Maximum magnetic field: ± 7.5 T
- Homogeneity of the magnetic field ($\varnothing 15$ mm sphere): 0.2%
- Room-temperature bore: $\varnothing 100$ mm
- Beam window dimensions:
vertical gap: 30 mm
vertical open angle: 3°
in plane open angle: 320°
- Cooling system: dry
- Cool down: 3 days
- Total thickness of Al in the beam: 30 mm
- Additional sample environment available (CC, CCR, HTF)

MAG-V-15T



Specifications

- Maximum magnetic field: ± 13.2 (14.5) T
- Homogeneity of the magnetic field (20 mm x 12 mm cylinder): 0.85%
- Low-temperature bore: $\varnothing 20$ mm
- Beam window dimensions: 20 x 20 mm²
- Scattering angle: 320° horizontal
- Cooling system: LHE
- Additional sample environment available: dilution insert
- MAG-V-15T is available at instrument PANDA; for the use at other instruments ask the local contact

Dr. Jürgen Peters

juergen.peters@frm2.tum.de
+49.(0)89.289.14700



www.mlz-garching.de/se

Dr. Astrid Schneidewind

a.schneidewind@fz-juelich.de
+49.(0)89.289.14749


www.mlz-garching.de/se

Electromagnet with CCR

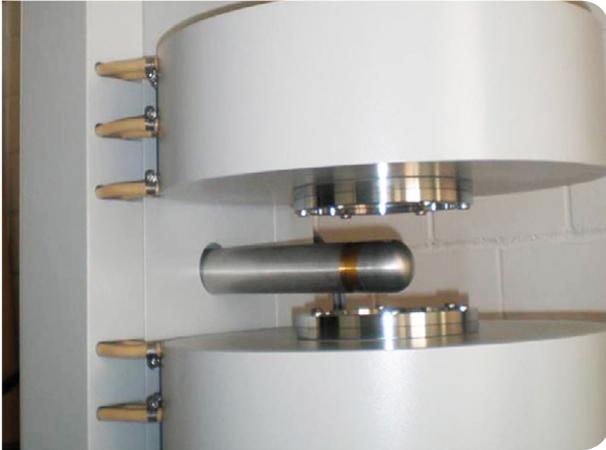

Specifications

- Two pieces available
- Maximum magnetic field: up to 1.4 T
- Temperature: 3 – 320 K
- used for small angle scattering and reflectometry

JVM1-5.0T active shielded



Specifications

- Maximum vertical magnetic field: ± 5.0 T (asymm.)
- Homogeneity of the magnetic field (\varnothing 15 mm, 30 mm in height, cylinder): 1.8%
- Beam window dimensions: horizontal and vertical angle: 30°
- Sample space: \varnothing 30 mm, height 50 mm
- Split: 30 mm
- Access angle: $\pm 5^\circ$ vertical
- Scattering angle: 330° horizontal
- VTI: 1.8 – 300 K
- Strayfield < 1 G at 1 m
- Cooling system: LHE
- Cool down: 4 hours
- Additional sample environment available: Kelvinox insert, base temperature 25 mK

Dr. Alexander Weber

Email:
Office:

al.weber@fz-juelich.de
+49.(0)89.289.10726


www.mlz-garching.de/se

HP 400 kN press



Specifications

HP 400 kN press:


- compressive force: 450 kN
- repulsive force: 220 kN
- dynamic pressure up to 10 Hz possible

Closed cycle cryostat CC-11-P adapted to be used with the HP 400 kN press.

- Temperature range: 20 K – 600 K
- Cooling power 2. stage: 1.5 W
- number of radiation shields: 1
- Heater cartridge 25 Ω / 100W
- Cool-down time RT – 20 K: 4 h
- Sample space: max \varnothing 16 mm, height 10 mm
- Pressure range: 200 MPa for \varnothing 16 mm anvil

A special challenge is to go to low temperatures at high pressures. The newly constructed high-pressure cryostat allows one to apply and vary the force in situ, even at low temperature. The sample cell is operated externally by the hydraulic HP 400 kN press. It can be mounted in a dewar and cooled down to below 10 K within four hours. With a piston-in-cylinder cell and a piston diameter of 16 mm a maximum pressure of 2.0 GPa (20 kbar) can be achieved at a sample size up to 30 mm length and 16 mm diameter. Reducing the piston diameter and the sample volume the pressure can be increased up to 100 kbar. Using anvil cells the pressure range extends to 30 GPa (300 kbar) and above. In addition, a dynamic pressure can be superimposed on the static pressure with a frequency of up to 10 Hz.

HPG-10 kbar

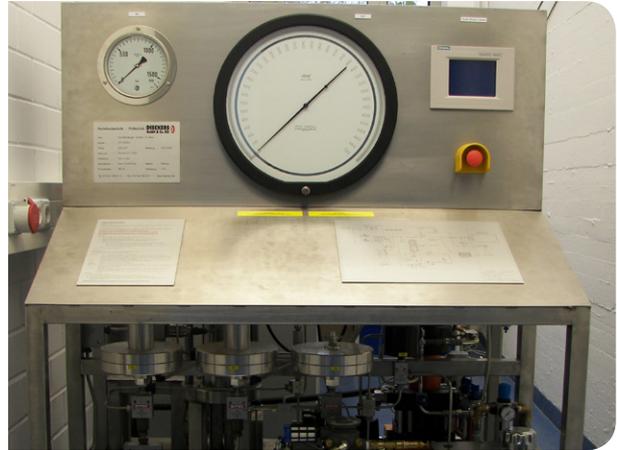

Specifications

- Pressure range $10^5 \text{ Pa} \leq p \leq 10^9 \text{ Pa}$

Pressure measurement:

- $10^5 \text{ Pa} \leq p \leq 7 \cdot 10^8 \text{ Pa}$ Heise high-precision manometer; resolution 0.1% ME
- $7 \cdot 10^8 \text{ Pa} \leq p \leq 10^9 \text{ Pa}$ transmitter; resolution $9 \cdot 10^6 \text{ Pa}$

To load gas pressure cells, the MLZ provides a gas compression unit for inert gas pressure up to 10 kbar. The unit can be operated either manually or remote controlled. The remote controlled operation allows for individually programmable pressure time profiles. The measurement system includes a high-precision manometer and a pressure transmitter that cover the range up to 10 kbar.

Dr. Jürgen Peters

Email:
Office:

juergen.peters@frm2.tum.de
+49.(0)89.289.14700



www.mlz-garching.de/se

Ultrasonic levitator (tec5 AG)



Specifications

- Frequency: 58 Hz
- Wavelength: 5.9 mm
- Optimal diameter (without drop deformation): 2.5 mm

Multiposition sample holders



Specifications

- 8 – 32 positions
- Thermalising devices:
Temperature -10°C – 120°C
(220°C with silicon oil)


Peltier furnace



Specifications

- 8 positions
- Temperature -20°C – 120°C

GISANS liquid flow cell

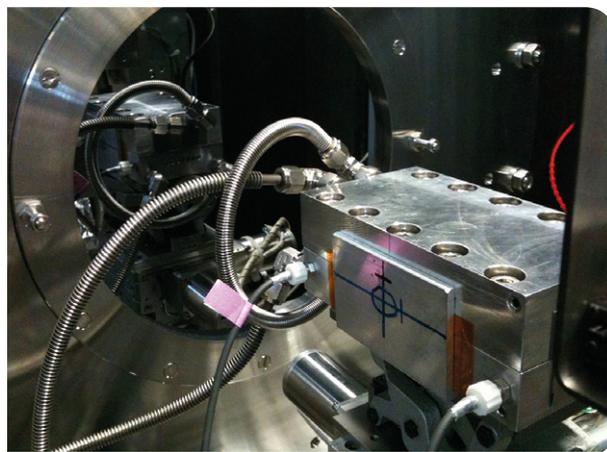


Specifications

- 15 x 5 x 2 cm³ silicon block
- Liquid volume of 7 ml
- Heatable by circulation bath

Dr. Alexander Weber

Email: al.weber@fz-juelich.de
Office: +49.(0)89.289.10726

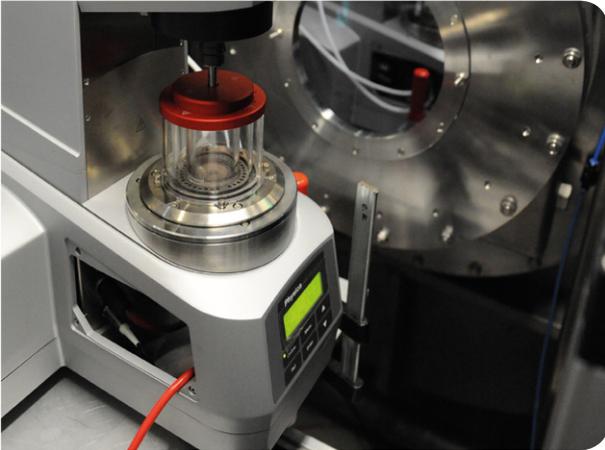


Rheometer RSA II**Specifications**

- Couette geometry for liquids
- Shear rate: 0.001 to 5000 s⁻¹
- Temperature: 20 to 200°C

Pressure cell**Specifications**

- Temperature: -40°C – 80°C
- Pressure: 5000 bar

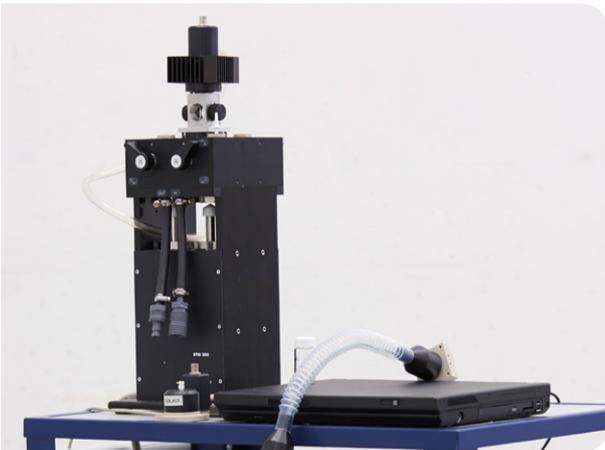Stopped-flow (Bio-Logic® SFM-300)**Specifications**

- 3 independently controlled syringes and two mixers
- Mixing ration from 1:1 to 1:100
- Exact control of flow rate
- Dead-time of mixing: < 20 msec
- Neutron-cell volume: 0.4 ml
- Recommended mixing volume: 1 ml


Dr. Alexander Weber

Email: al.weber@fz-juelich.de
Office: +49.(0)89.289.10726



www.mlz-garching.de/se


Humidity cell

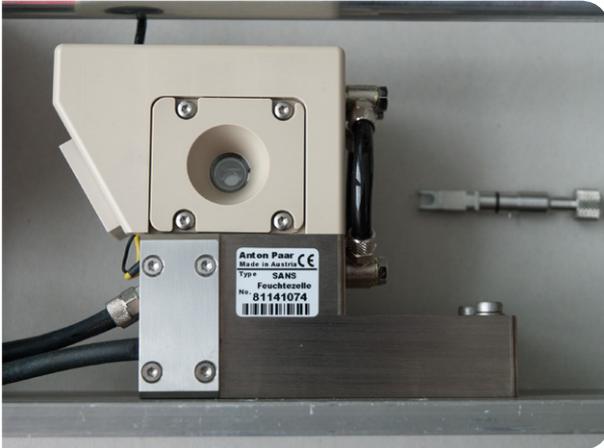

Specifications

- Temperature: 10°C – 90°C
- Humidity range: 5% – 95%
- Sample size: Ø 8 mm, thickness: 1 mm

Humidity generator

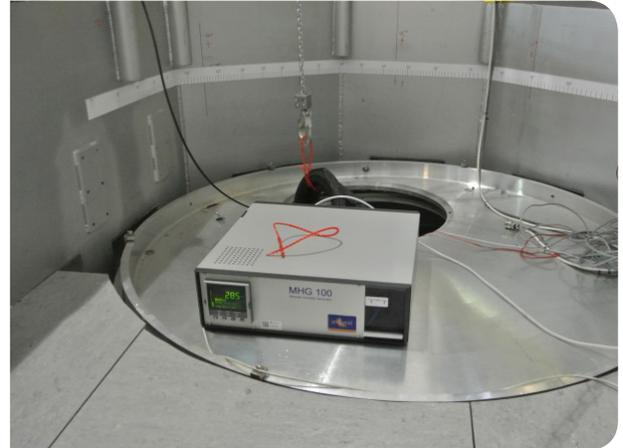

Specifications

- Massflow 0...15 l/min, saturated flow up to 10 l/min
- Rel. humidity 5% ... 95%
- Temperature +5°C ...75 (95)°C, depends on humidity
- 2 tempered pipes (1.5 m)



Measurement with customer cell at SPHERES on fuel-element parts.


Dr. Alexander Weber

Email: al.weber@fz-juelich.de
Office: +49.(0)89.289.10726



The MLZ sample environment groups supply a variety of liquid refrigerant free cryostats adapted for the different needs on the instruments. Besides standard closed cycle cryostats, either of toploding type or sample directly mounted to coldhead, adapted versions for the use inside the vertical field 7.5 T magnet system, the 40 kN press or equipped with a thermoswitch for temperatures above RT are available. For special requirements please contact the Sample Environment Group or your local contact to discuss details.

CCR - closed-cycle cryostat with sample tube

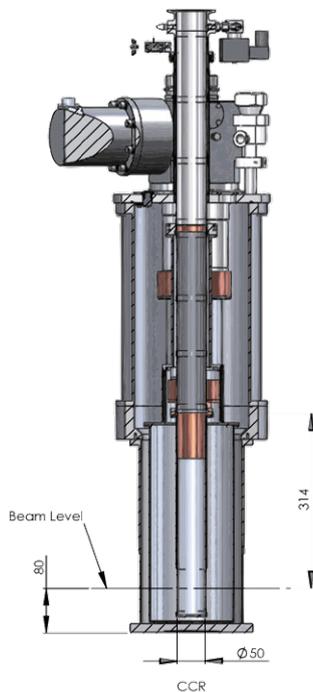

Cut through CCR cryostat.

The liquid refrigerant free closed cycle cryostats of the CCR-type are designed for a fast change of samples. They are based on pulse-tube technique refrigerators with a 6 kW water cooled compressor unit. The sample tube is connected to the pulse tube cold plate with a copper heat exchanger.


The sample tube is filled with exchange-gas. Temperature regulation is achieved by a sensor and a heater attached to the sample tube. The sample-holder is in general mounted to a sample stick (see fig. 2). Temperatures ranging from 3.5 K to RT can be regulated. To reach for temperatures above RT to 700 K the sample space has to be evacuated and a special high-temperature sample-stick has to be used. For this case temperature control is provided by a heater and sensor mounted on the sample stick.

Common features

- Based on pulse-tube refrigerator cooling power 2nd stage: 1000 mW
- 1 radiation shield connected to 1st stage
- Temperature sensor Cernox® 1.4 K – 325 K
- Heater 25 Ω / 100 W
- Temperature range: with exchange gas in sample space 2,8 K – 300 K
- Extended temperature range: evacuated sample space and HT sample stick T < 700 K
- Typical cool down times of the cryostat RT to base temperature min. 2,5 h
- Cool down time for sample change with cryostat at base temperature < 1.5 h

Dimensions

- Diameter sample space: < 50 mm and < 80 mm
- height of sample space (beam window): approx. 75 mm

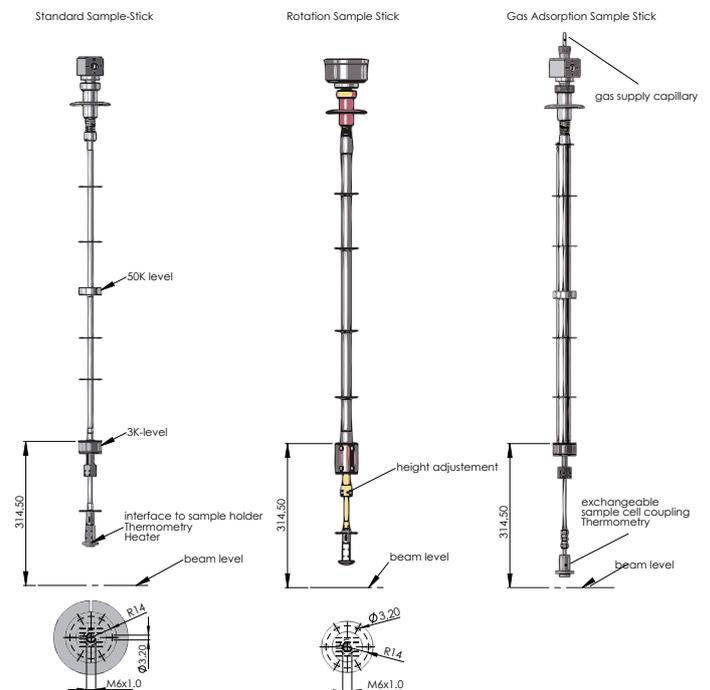

Available sample sticks provided for CCR cryostats.

Dr. Jürgen Peters

Email: juergen.peters@frm2.tum.de
Office: +49.(0)89.289.14700



www.mlz-garching.de/se

CCI - low temperature inserts

To reach temperatures below 3 K ^3He and $^3\text{He}/^4\text{He}$ insert cryostats for the CCR systems are available. Please note an additional preparative time of up to 4h is needed before first cool-down.

CCI specifications

	CCI - ^3He	CCI $^3\text{He}/^4\text{He}$
Temp. range	450 mK – 10 K	50 mK – 1 K
Cooling power	1 mW (500 mK)	15 μW (100 mK)
Heater	10 Ω	10 Ω
Cool-down time	5 h	7 h
Sample \varnothing	30 mm	30 mm
Sample space-height	145 mm	70 mm
Distance cold plate to beam level	106.5 mm	27 mm

CC- closed-cycle cryostats


The closed cycle cryostats of the CC-type are based on SHI-RDK-2025D and SHI-RDK-101D cold-heads mounted with differing isolation vacuum tails. The sample space of these liquid-cryogen free closed-cycle cryostats is inside the isolation vacuum. The thermalisation of the sample is achieved by the thermal conductivity of the sample holder and the cold-plate. The temperature sensor for temperature control is mounted on the cold-plate. To avoid temperature gradients a sample mounting with adequate thermal conductivity is needed. In case of sample with poor thermal conductivity, the usage of an exchange-gas (He) containing sample can or the use of a CCR-type cryostat should be considered.

CC-1*, CC-2-PUMA, CC-3, CC-4-PANDA

- Temperature range: 2.8 K – 300 K
- Cooling power 2. stage: 250 mW
- Radiation shields attached to 1st stage: 1
- Thermometry: Si-diode; Cernox® 1.4 K – 325 K
- Heater cartridge 25 Ω / 100 W
- Cool-down time RT – 2.8 K: 2.5 h
- Total height of sample space: 110 mm
- Diameter of sample space: 60 mm


* CC-1 prepared for the condensation of non-corrosive gases up to P = 100 bar

3He-CCI-2

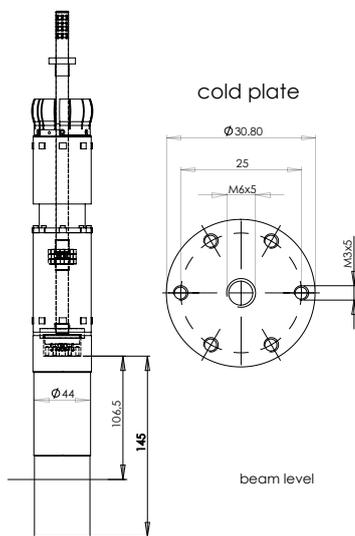


Cut through the CCI inserts.

3He/4He-CCI

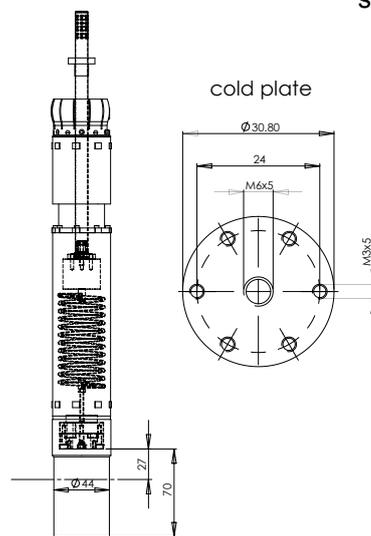

Cut through CC-3.

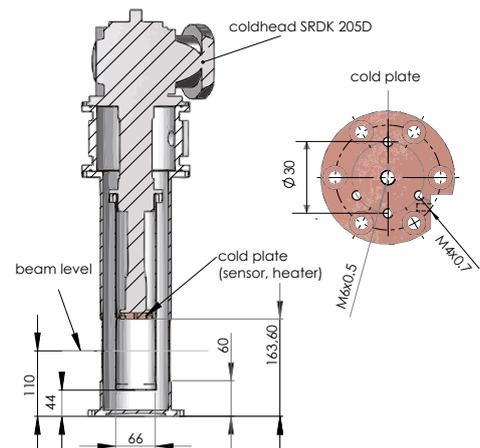


Dr. Jürgen Peters

Email: juergen.peters@frm2.tum.de
Office: +49.(0)89.289.14700



CC-5, CC-6

- Temperature range: 2.2 K – 300 K
- Cooling power 2. stage: 110 mW
- Radiation shields attached to 1st stage: 1
- Thermometry: Si-diode; Cernox® 1.4 K – 325 K
- Heater cartridge 25 Ω / 100 W
- Cool-down time RT – 2.8 K: 2 h
- Sample space: \varnothing 10 mm, height 10 mm

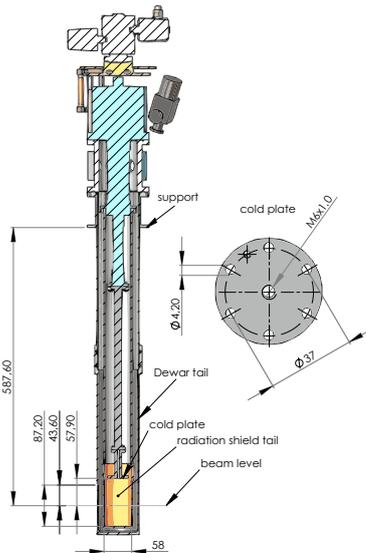


Cut through CC-5/ CC-6.

CC-7-CCM

- Temperature range: 2.2 K – 300 K
- Cooling power 2. stage: 110 mW
- Radiation shields attached to 1st stage: 1
- Thermometry: Si-diode; Cernox® 1.4 K – 325 K
- Heater cartridge 25 Ω / 100 W
- Cool-down time RT – 2.8 K: 4 h
- Sample space: \varnothing 55 mm, height 85 mm

CC-8, CC-9-RESEDA, CC-10-PUMA


- Temperature range: 2.8 K – 600 K
- Cooling power 2. stage: 110 mW
- Radiation shields attached to 1st stage: 1
- Thermometry: Pt1000 IST 10 K – 600 K
- Heater cartridge 25 Ω / 100 W
- Cool-down time RT – 2.8 K: 3 h
- Cool-down time 600 K – 2.8 K: 5 h
- Total height of sample space: < 136 mm (< 85 mm CC-9 / CC-10)
- Diameter of sample space: < 138 mm (< 55 mm CC-9 / CC-10)

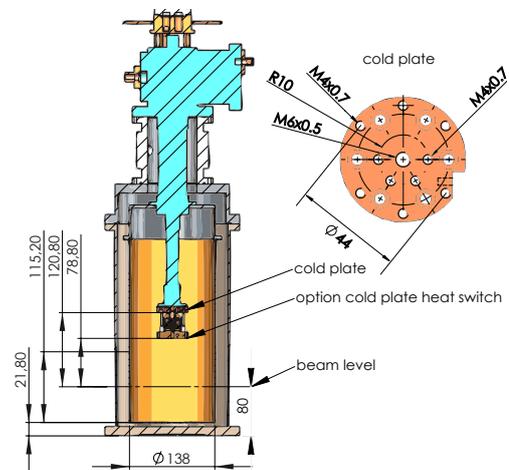

Cut through CC-8.



JTC 1-4

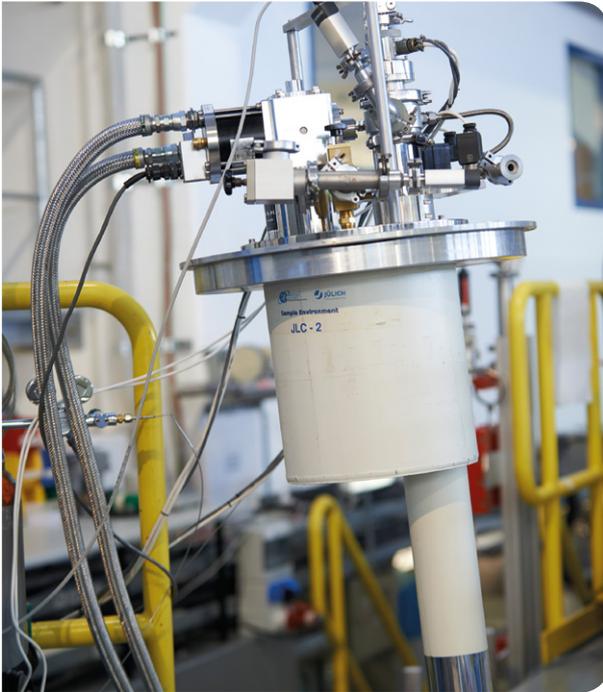

Standard cryostat.

Common features

- Heater 25 Ω / 100 W
- Temperature range:
with exchange gas in sample space
3 K – 300 K
- Extended temperature range:
evacuated sample space and HT sample stick
T < 700 K
- Typical cool down times of the cryostat
RT to base temperature
min. 2.5 h
- Cool down time for sample change with cryostat
at base temperature < 1.5 h


Dimensions

- Diameter sample space:
< 60 mm
- Height of sample space (beam window):
approx. 100 mm

Available sample sticks

- Standard stick (T 3 – 700 K)

The sample position and orientation in the beam can be adjusted by a simple height adjustment and rotation at the top part of the sample stick during the experiment.

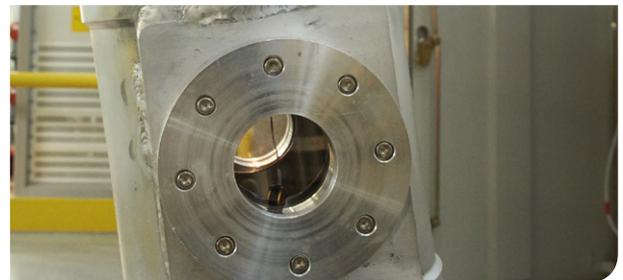

Cryostat with sapphire window.

Dr. Alexander Weber

Email: al.weber@fz-juelich.de
Office: +49.(0)89.289.10726



HTF – high temperature furnace


The heater of the HTF consist of a resistive Nb double cylinder element. Radiation shields made of Nb reduce the thermal loss. The sample is mounted on a sample rod top down, using a M8 screw. Maximum reachable temperature is 1900°C. For temperatures up to 300 °C the furnace can be filled with Ar or He exchange gas to improve regulation stability. Temperature sensors are type C thermocouples.

Further a special version of the furnace with a reduced diameter for usage inside the CCM-7.5T vertical magnet is available.

Specifications

Designation: HTF 1, 2, 3, 4 - SPODI

- Temperature range (vacuum): RT – 1900 °C
- Temperature range (exchange gas): RT – 900°C
- Thermometry: type C thermocouple
- Total height of sample space: 100 mm
- Diameter of sample space: 45 mm
- Sample rod tail: M8 (male)

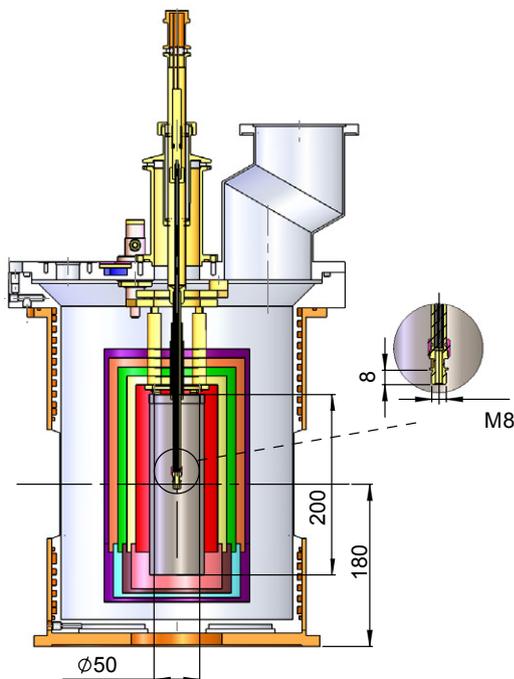

Cut through a standard HTF.

IRF – Infrared light furnace


Compact dimensions are the key feature of the light furnace. The sample is placed at the focus of four halide lamps, therefore restricting the sample volume to a few mm³. The sample can be heated in vacuum or any convenient atmosphere up to ambient pressure respectively. Of course the maximum temperature depends on pressure.

A dedicated version of this furnace for use together with a load frame is available. Twelve halogen bulb lamps heat the sample up to 900 °C with an almost homogeneous temperature distribution. The sample volume has 6-8 mm diameter and 30 mm length. By means of an additional heat shield temperatures beyond 950 °C are expected.

Specifications

Designation: IRF 1, 2

- Temperature range (vacuum): RT – 1200 °C
- Temperature range (Ar 100mbar): RT – 300 °C
- Thermometry: type K, R, S thermocouple
- Height of sample space: 15 mm
- Diameter of sample space: 10 mm
- Sample rod tail: M4 (male)

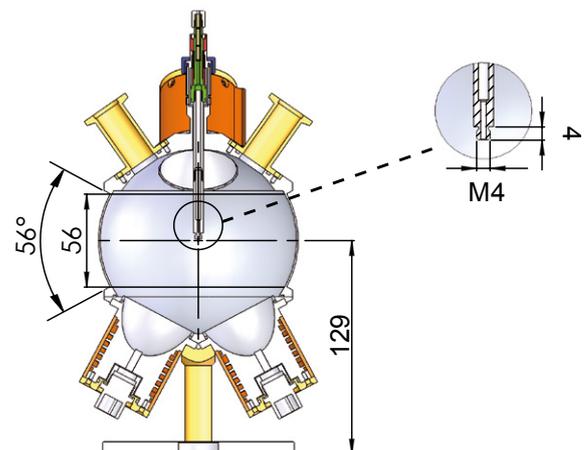

Cut through IRF light furnace.

Dr. Jürgen Peters

Email: juergen.peters@frm2.tum.de
Office: +49.(0)89.289.14700



www.mlz-garching.de/se


PF – Polarised neutron furnace

For experiments using a polarised neutron beam a furnace with special bifilar heater cartridges is available for the temperature range up to 700 °C.

Specifications

Designation: PF 1

- Temperature range (vacuum): RT – 700 °C
- Thermometry: Pt100
- Height of sample space: 80 mm
- Diameter of sample space: 45 mm
- Sample rod tail: M6 (male)

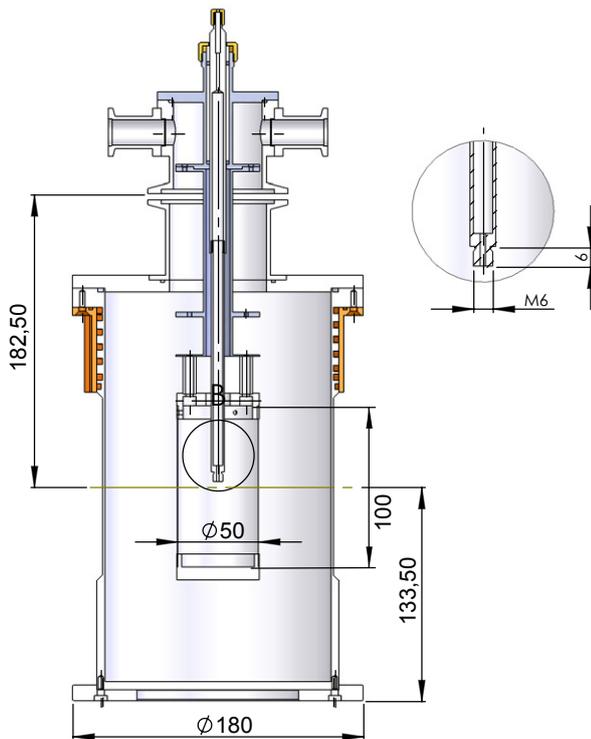

Cut through PF 1 furnace for polarised neutron experiments.


CTF1 – Circulation thermostat furnace

For the temperature range -20 °C to 200 °C a furnace using a thermalised circulating medium allows for a precise regulation of the particularly homogeneous sample temperature.

Specifications

Designation: CTF 1

- Temperature range (vacuum): -30 °C – 200 °C
- Thermometry: Pt100
- Height of sample space: 80 mm
- Diameter of sample space: 48 mm
- Sample rod tail: M8 (male)

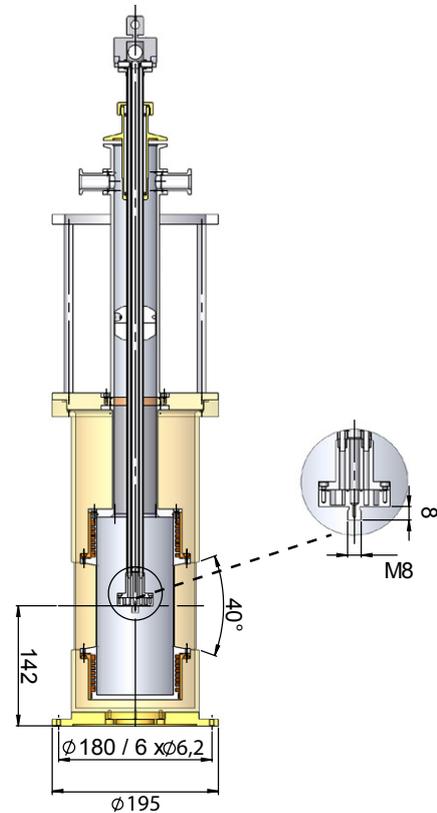

Cut through the CTF circulation thermostat furnace.

Dr. Jürgen Peters

Email: juergen.peters@frm2.tum.de
Office: +49.(0)89.289.14700



Special application furnace: battery furnace





Battery furnace set-up.

For research, e.g. on rechargeable energy storage systems such as lithium-ion batteries, dedicated heater elements are available. Depending on battery size and temperature specification the furnaces have to be designed and manufactured (for details, please get in contact with the Sample Environment Group). Heating is accomplished by an appropriate set of high performance cartridge heaters allowing for temperatures up to 700°C. A remote controlled power supply provides the heater current. Optional temperatures below room temperature (RT) are possible.

Specifications

- Temperature : RT – 700 °C (optional < RT)
- Thermometry : PT 1000
- Dimension: adapted to battery size

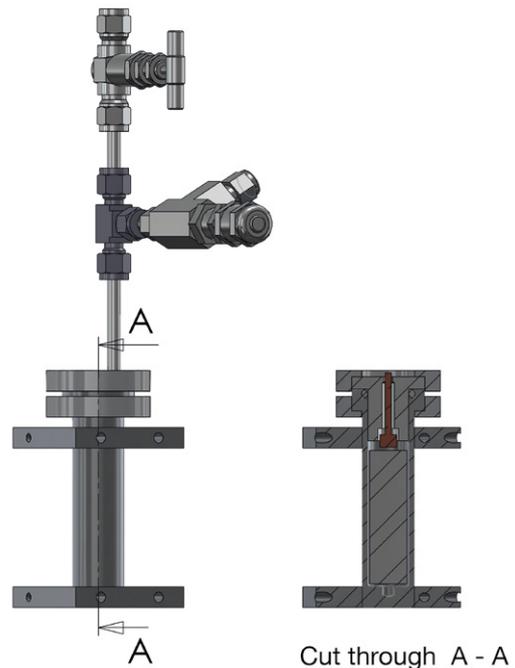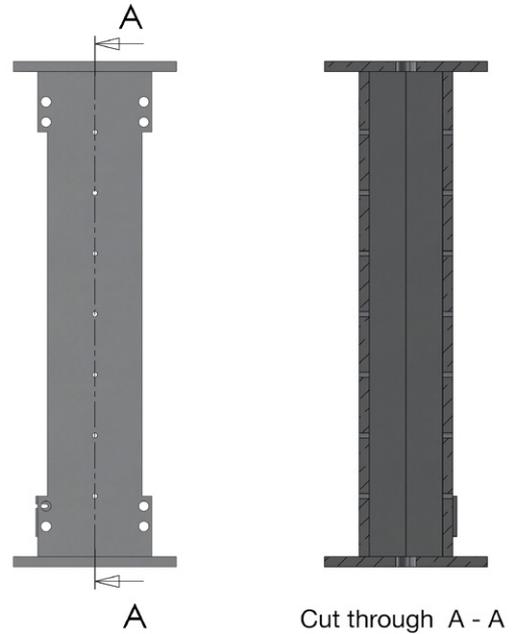

Examples for battery furnaces available at MLZ.

Dr. Jürgen Peters

Email: juergen.peters@frm2.tum.de
Office: +49.(0)89.289.14700



Biology Laboratory



The biology laboratory is a 40 m² facility which serves to prepare biological and soft matter samples for measurements at the neutron scattering instruments. It is equipped to handle samples in inert gas atmosphere within a glovebox and pre-characterise samples by UV-VIS spectroscopy and dynamic light scattering. A low-temperature freezer can be used to store samples at temperatures down to -80°C.

Equipment

- Fine balance
max. weight 405 g, resolution 1 mg
- Elga Labpure water system (18 MOhm)
- Glovebox H₂O and O₂ < 1ppm, Ar atmosphere
- pH-meter with a 3 mm diameter probe
- Freezer -80°C
- Refrigerator +4°C and freezer -18°C
- UV-VIS spectrometer NanoDrop
- Malvern Zetasizer Nano S
- Optical microscope
- Thermostated centrifuge up to 60 000 G
- Rheometer
- ÄKTA FPLC for protein purification

Chemistry Laboratory



The chemistry and sample preparation laboratory is a 45 m² facility equipped with double hood to prepare samples for subsequent measurements on the neutron scattering instruments. It offers a basic tool set to handle samples from users, prepare solutions or even perform chemical reactions resulting in specific modifications. Samples can be handled under inert gas atmosphere like N₂ or Ar.

Equipment

- Fine balance
max. weight 405 g, resolution 1 mg
- Balance
max. weight 4.2 kg, resolution 0.1 g
- Refrigerator with freeze box -18°C
- Vacuum heating cabinet
max. temperature 300°C, min. pressure 25 mbar, inertgas: N₂, Ar
- Thermostatic bath, 20 – 80°C
- Ultra-sonic bath
- Rotary evaporator
- Rheometer Rheometrics Solid Analyzer RSA II, temperature from liquid N₂ to 500°C

Dr. Tobias Schrader

Email: t.schrader@fz-juelich.de
Office: +49.(0)89.289.10743



www.mlz-garching.de/user-labs

Dr. Henrich Frielinghaus

h.frielinghaus@fz-juelich.de
+49.(0)89.289.10706



www.mlz-garching.de/user-labs

Materials Science Laboratory



In collaboration with HZG and TUM, the Materials Science Lab serves as laboratory to prepare samples (polishing, cutting, etching, and annealing) for subsequent neutron scattering measurements. It is equipped with several analytic instruments for instance X-ray Diffractometer (XRD) mainly for battery research and high-performance materials, Differential Scanning Calorimeter (DSC), Micro-Hardness and Optical Microscopes to perform complementary analytics.

Equipment

Analytic instruments:

- PANalytical Empyrean XRD (Cu, Mo, XRR, Transm. Reflection)
- Differential Scanning Calorimetry (-80°C – 600°C)
- Struers Micro Hardness DuraScan (HV 0.1 – 10)
- BioLogic – Bi-Potentiostat
- Digital Microscope & Stereomicroscope

Preparation instruments:

- GERO – Programmable High Temperature Tube furnace
- ATM BRILLANT 200 – corundum wheel saw
- ATM OPAL 410 – resin sample hot embedding
- ATM SAPHIR 520 – polishing machine
- Binder – vacuum heating cabinet (200°C)
- Well Diamond Wire Saw

Sample Preparation Laboratory



For the preparation or post-treatment of samples, the MLZ offers a sample preparation laboratory in the Neutron Guide Hall West. In addition to the equipment listed below, the sample preparation lab provides a supply of argon, helium, compressed air, and ultrapure water. Usually, small amounts of commonly used solvents and laboratory dishes are available. Additionally, you can use the official proposal form to point out your special requirements, so that your local contact and the Sample Environment Group can prepare appropriate conditions for your experiment in time.

Mobile glove boxes are installed in both, Neutron Guide and Experimental Hall, to be used close to the instruments.

Equipment

- Fume hood
- Glove box
- Refrigerator
- Freezer
- Scales
- Stirrers
- Heating oven
- Ultrasonic bath

Dr. habil. Ralph Gilles

ralph.gilles@frm2.tum.de
+49.(0)89.289.14665



www.mlz-garching.de/user-labs

Dr. Jürgen Peters

juergen.peters@frm2.tum.de
+49.(0)89.289.14700



www.mlz-garching.de/user-labs

TEM: JEOL 200 kV JEM-FS2200




Transmission electron microscopy is a complementary technique with neutron scattering in the frame of soft matter studies. With the TEM, real space investigations are performed to access information about shape, size, and size distribution of particles, self-assembly and aggregation. Users will be supported by JCNS staff (M.S. Appavou) to conduct the suitable preparation and TEM investigation.

Specifications

- JEOL 200 kV JEM-FS2200 with a field emission gun (FEG) and energy filter
- Magnification from 50x to 1 000 000x
- Image resolution of 0.2 nm in point and 0.1 nm in lattice
- Tietz CCD camera with 2048 × 2048 pixels
- Grid preparation by carbon sputtering on bare grids
- Glow discharge process to make the coated grids hydrophilic
- Plunge freezer to prepare cryo-specimen for direct imaging (Cryo-DI), for freeze fracture and direct imaging (FFDI) or freeze fracture before shadowing with Pt/C or Ta/W, so-called freeze fracture and electron microscopy (FFEM)
- Cryo-specimen investigation at low temperature around -180°C under vacuum.
- Ultramicrotome for room temperature and cryo sectioning

Thin film laboratory: MBE system

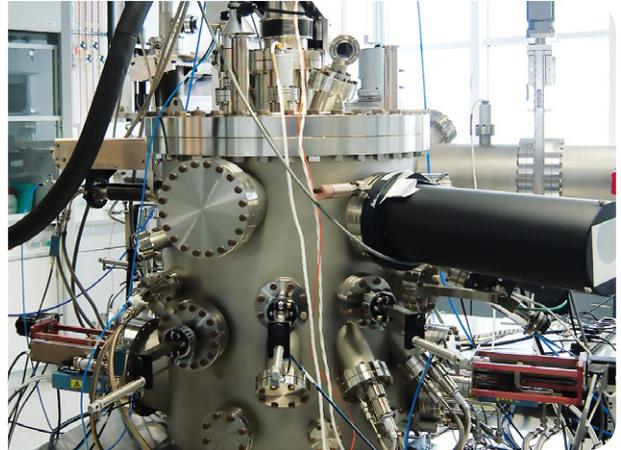

The Jülich Centre for Neutron Science operates a state-of-the-art oxide **M**olecular **B**eam **E**pitaxy (MBE) system (DCA Instruments Oy, Finland) to prepare tailored samples for the investigation with the neutron reflectometer MARIA, other neutron scattering instruments or methods.

The system is equipped with six effusion cells and two e-guns (each with four crucibles) for co-deposition. An RF plasma atom source enables atomic oxygen introduction into the chamber to grow oxides. Thin film growth analysis is performed by reflection high energy electron diffraction (RHEED). Additionally low energy electron diffraction (LEED) and Auger electron spectroscopy (AES) are offered for thin film characterisation.

Sample growth will be performed in strong collaboration with the instrument scientist.

Please ask for the feasibility of your desired sample!

Equipment and services

- Thin film growth using MBE technique
- Atomic Force Microscopy (AFM)
- Magneto-Optical Kerr Effect (MOKE) set-up, magnetic field up to 0.7 T, temperature range 4.5 – 420 K

Dr. Marie-Sousai Appavou

Email: m.s.appavou@fz-juelich.de
Office: +49.(0)89.289.10747



www.mlz-garching.de/user-labs

Dr. Sabine Pütter

s.puetter@fz-juelich.de
+49.(0)89.289.10742



www.mlz-garching.de/user-labs

X-ray Tomography: v|tome|x



In order to offer an X-ray facility complementary to the neutron tomography station ANTARES, the FRM II and the Chair of Biomedical Physics at the TUM have started to operate a high resolution computer tomography facility. The “micro CT VtomeX” is also available to users of the MLZ. Samples, that can be examined, include geo and composite materials, semiconductors and biomedical specimens.

The X-ray tomography station is featured by its high flexibility: Due to two different exchangeable X-ray tubes, both, pictures with high resolution and pictures with lower resolution but higher contrast, can be taken. The detector allows for fast and highly contrasted pictures. The reconstruction of the data is accomplished within a few minutes due to a cluster of four computers using graphic cards to calculate the images.

Possible scanning parameters

- Max. sample diameter 230 mm
- Max. sample height 420 mm
- Min. resolution < 1 micron (isotropic)
- Max. voxel size of reconstruction 2048³
- Max. X-ray energy 240 keV
- Max. sample weight 10 kg
- Typical scanning time 1–120 min

Location

Institute for Medical Engineering (IMETUM)
Boltzmannstr. 11, 85748 Garching

Dr. Klaus Achterhold

klaus.achterhold@ph.tum.de
+49.(0)89.289.12559



www.mlz-garching.de/user-labs



General information for users



Access to the MLZ

User Office

From idea to publication

General information for users


The User Office organises the scientific experiments at the MLZ, handles the proposal rounds and takes care of everything scientists visiting the neutron source in the scope of their work have to deal with.


You can call us, drop in when visiting the MLZ, or find all available information at

www.mlz-garching.de/user-office

There you learn about

- our Terms of References, for example regarding publications,
- reactor cycles and upcoming proposal deadlines,
- the online User Office System,
- workflows,
- requirements for the access to the site,
- radiation protection regulations,
- accommodations at Garching,
- the public transport system,
- possibilities regarding financial support,
- forms and templates for proposals and reports ready for download,
- and much more...

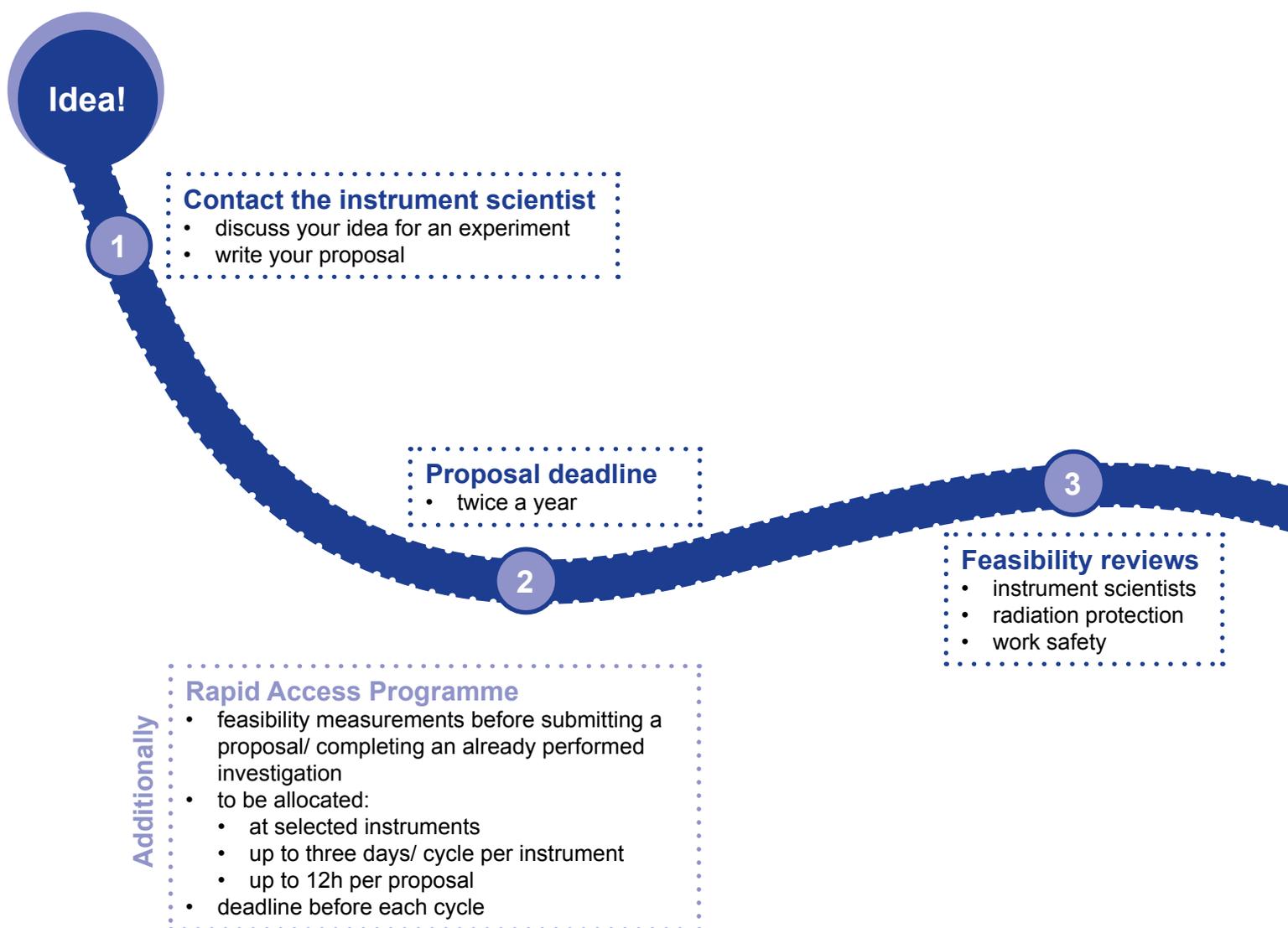

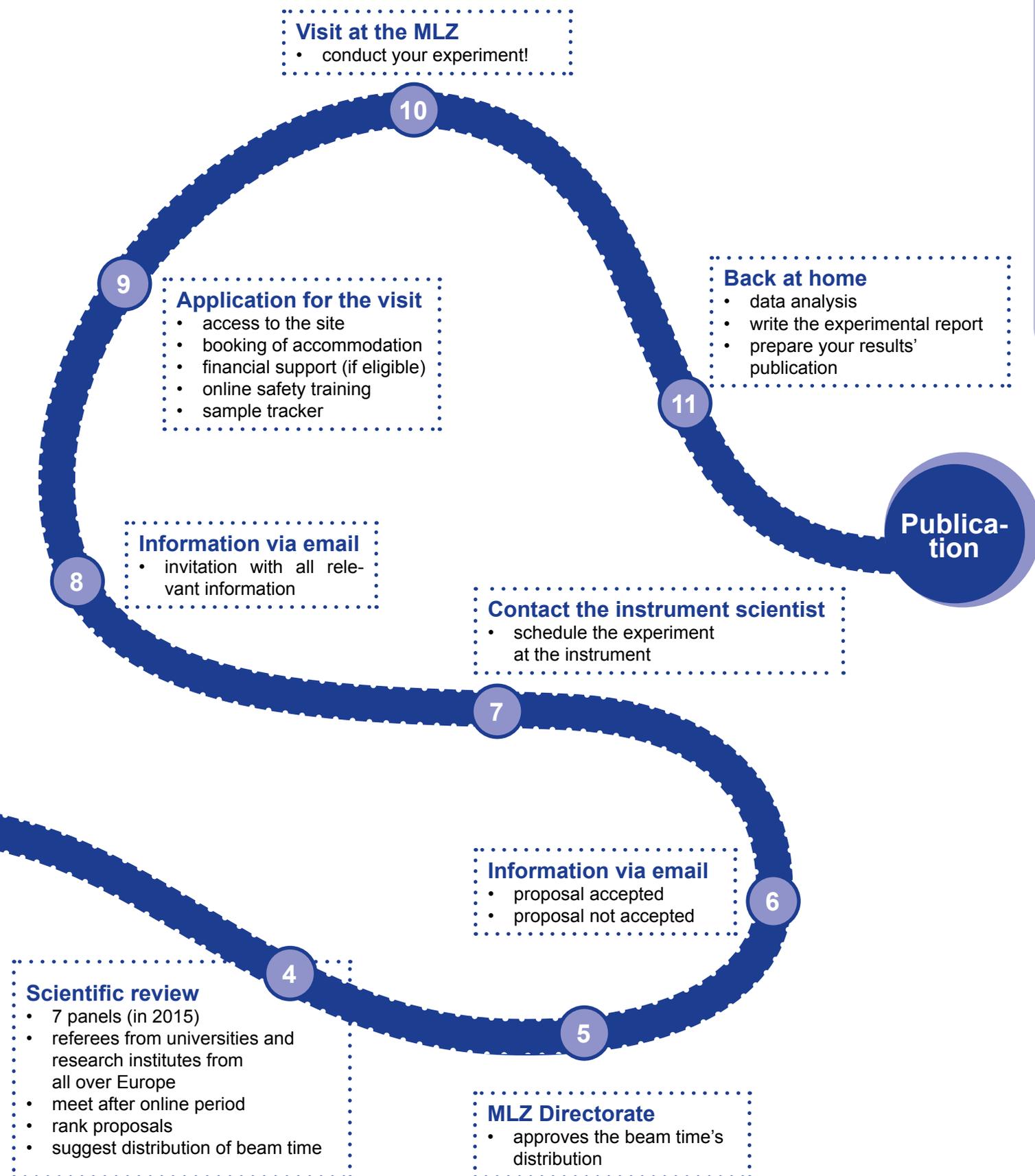

Dr. Flavio Carsughi

Dr. Ina Lommatzsch

Ramona Bucher



Email: f.carsughi@fz-juelich.de
Office: +49.(0)89.289.10703

ina.lommatzsch@frm2.tum.de
+49.(0)89.289.10794

r.bucher@fz-juelich.de
+49.(0)89.289.11751

www.mlz-garching.de/user-office

Access to the MLZ

Access to the site

- You had to apply online for the visit well in advance before your experiment starts.
- A valid passport or ID-card is required. Driving licenses or other personal documents are not sufficient!

Radiation protection requests from Users working in Germany:

- 'Strahlenpass', dose record not older that 3 months
- dosimeter (capable of detecting neutron and gamma radiation)
- the home institution needs a valid license according to §15 StrlSchV as well as an 'Abgrenzungsvertrag'

- access to the internet
- control experiments remotely
- access data


Reception

ON ARRIVAL

- getting a badge
- blue folder with process slip, experimental safety sheet, sample tracker print-out, further information
- your local contact will be informed

ON DEPARTURE

- return your badge
- return the completed process slip

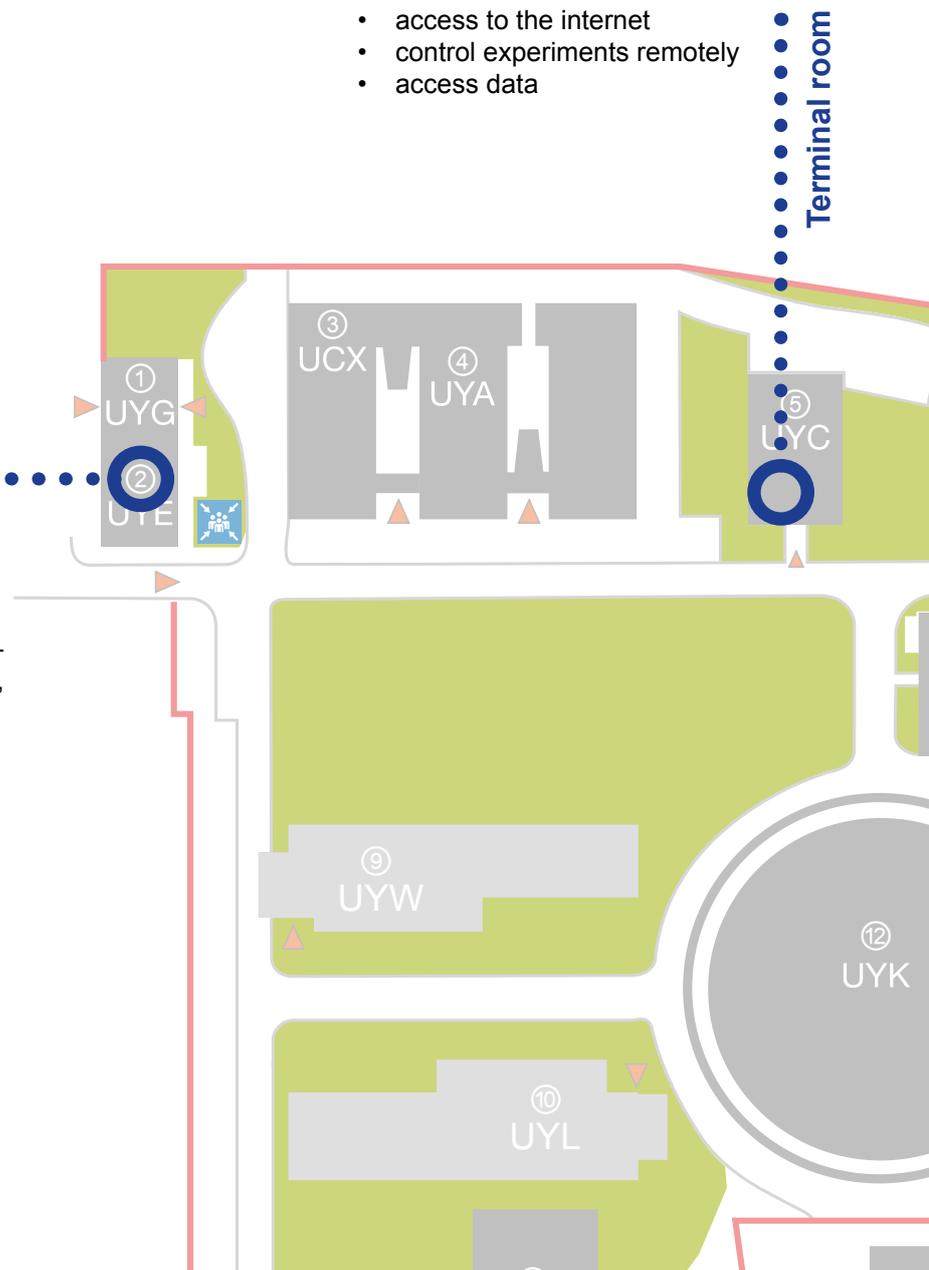

Radiation protection requests from Users working outside Germany:

- subdivided in “radiological workers” and persons, who have never worked in a controlled area
- download the appropriate form from the User Office Webpage
- complete it, sign it, and present it upon arrival


ON ARRIVAL

- register for dosimetry
- register your samples

ON DEPARTURE

- check-out of dosimetry
- discuss your sample’s clearance

Radiation protection office

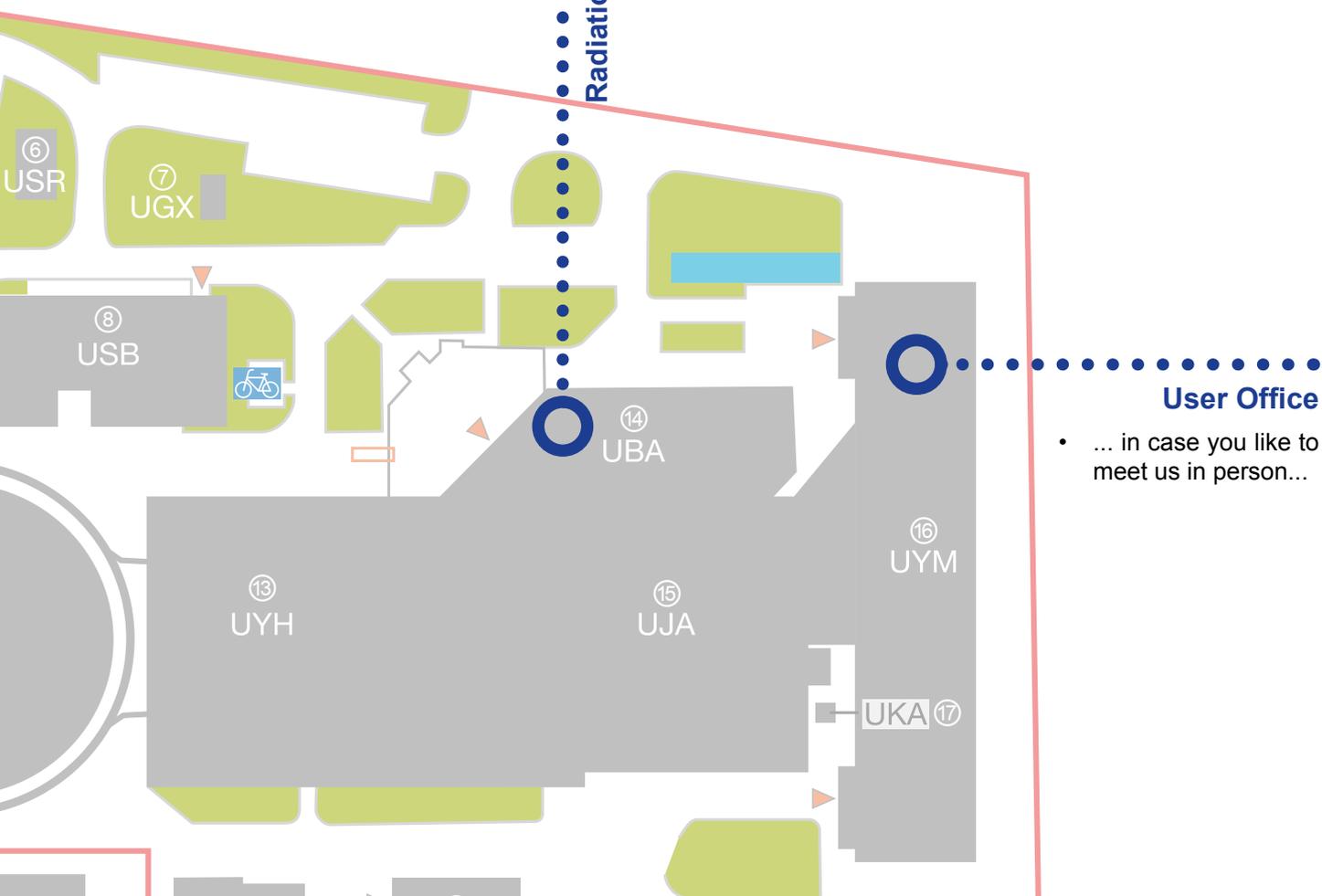

User Office

- ... in case you like to meet us in person...

Dr. Flavio Carsughi

Email: f.carsughi@fz-juelich.de
Office: +49.(0)89.289.10703

Dr. Ina Lommatzsch

ina.lommatzsch@frm2.tum.de
+49.(0)89.289.10794

Ramona Bucher

r.bucher@fz-juelich.de
+49.(0)89.289.11751



www.mlz-garching.de/user-office

Partner institutions




Bayerisches Geoinstitut
Universität Bayreuth
www.bgi.uni-bayreuth.de



Georg-August-Universität Göttingen

- Institut für Physikalische Chemie
www.uni-pc.gwdg.de/eckold
- Geowissenschaftliches Zentrum
www.uni-goettingen.de/de/125309.html

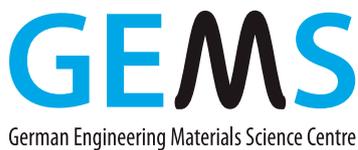

German Engineering Materials Science Centre GEMS
Helmholtz-Zentrum Geesthacht GmbH
www.hzg.de/institutes_platforms/gems/



Jülich Centre for Neutron Science JCNS
Forschungszentrum Jülich GmbH
www.jcns.info



Karlsruher Institut für Technologie

- Institut für Angewandte Materialien – Energiespeichersysteme (IAM-ESS)
www.iam.kit.edu



Ludwig-Maximilians-Universität München

- Sektion Kristallographie
www.lmu.de/kristallographie
- Sektion Physik
www.softmatter.physik.uni-muenchen.de



MAX-PLANCK-GESELLSCHAFT

Max-Planck-Institut für Festkörperforschung, Stuttgart
www.fkf.mpg.de



RWTH Aachen

- Institut für Kristallographie
www.xtal.rwth-aachen.de
- Institut für Anorganische Chemie
www.ac.rwth-aachen.de



Technische Universität Clausthal

- Institut für Werkstoffkunde und Werkstofftechnik
www.iww.tu-clausthal.de



Technische Universität Dresden

- Institut für Festkörperphysik
www.physik.tu-dresden.de/ifp



Technische Universität München

- E13 – Lehrstuhl für Funktionelle Materialien
www.e13.physik.tu-muenchen.de



Technische Universität München

- E18 – Lehrstuhl für Experimentalphysik I
www.e18.ph.tum.de



E21
Arbeitsgebiet stark
korrelierte Elektronensysteme

Technische Universität München

- E21 – Lehrstuhl für Neutronenstreuung
www.e21.ph.tum.de



Technische Universität München

- Exzellenzcluster „Origin and Structure of the Universe“
www.universe-cluster.de



Klinikum rechts der Isar

Technische Universität München

- MRI - Klinikum Rechts der Isar
www.med.tum.de



RCM
Radiochemie München

Technische Universität München

- RCM - Radiochemie München
www.rcm.tum.de



TECHNISCHE
UNIVERSITÄT
WIEN
Vienna University of Technology

Technische Universität Wien

- Neutronen- & Quantenphysik
Forschungsbereich am Atominstitut Wien
Arbeitsgruppe Abele
ati.tuwien.ac.at/research_areas/neutron_quantum_physics/

der Bundeswehr
Universität  *München*

Universität der Bundeswehr München

- Institut für Angewandte Physik und Messtechnik
www.unibw.de/lrt2

Universität zu Köln



Universität zu Köln

- Institut für Kernphysik
www.ikp.uni-koeln.de
- II. Physikalisches Institut
www.ph2.uni-koeln.de

Imprint

Publisher

Technische Universität München
Forschungs-Neutronenquelle
Heinz Maier-Leibnitz (FRM II)
Lichtenbergstraße 1
D-85747 Garching

User Office
Heinz Maier-Leibnitz Zentrum (MLZ)
++49.(0)89.289.10703 /.10794 /.11751
useroffice@mlz-garching.de
www.mlz-garching.de

Editors

Dr. Ina Lommatzsch
Dr. Connie Hesse

Design and typesetting

Ramona Bucher
Dr. Ina Lommatzsch

Third revised edition (2015)
Second revised edition (2013)
First edition (2010)

Picture credits

Bernhard Ludewig
66

Forschungszentrum Jülich:
3 (right)

Andreas Heddergott / Astrid Eckert, TUM:
8, 88 (left), 58

Wenzel Schürmann, TUM:
All instrument pictures;
Cover, Backcover, 3 (left), 88 (right), 89 (left, top),
90 (right, top), 104


Editors and authors:
other images

Schematic illustrations of the instruments:
Ramona Bucher
Julia Fridgen
Sabrina Kressierer
Dr. Peter Link
Christine Sturz



MLZ User Office
Phone: +49.(0)89.289.10703 / .10794 / .11751
Email: useroffice@mlz-garching.de

MLZ Instrument Suite

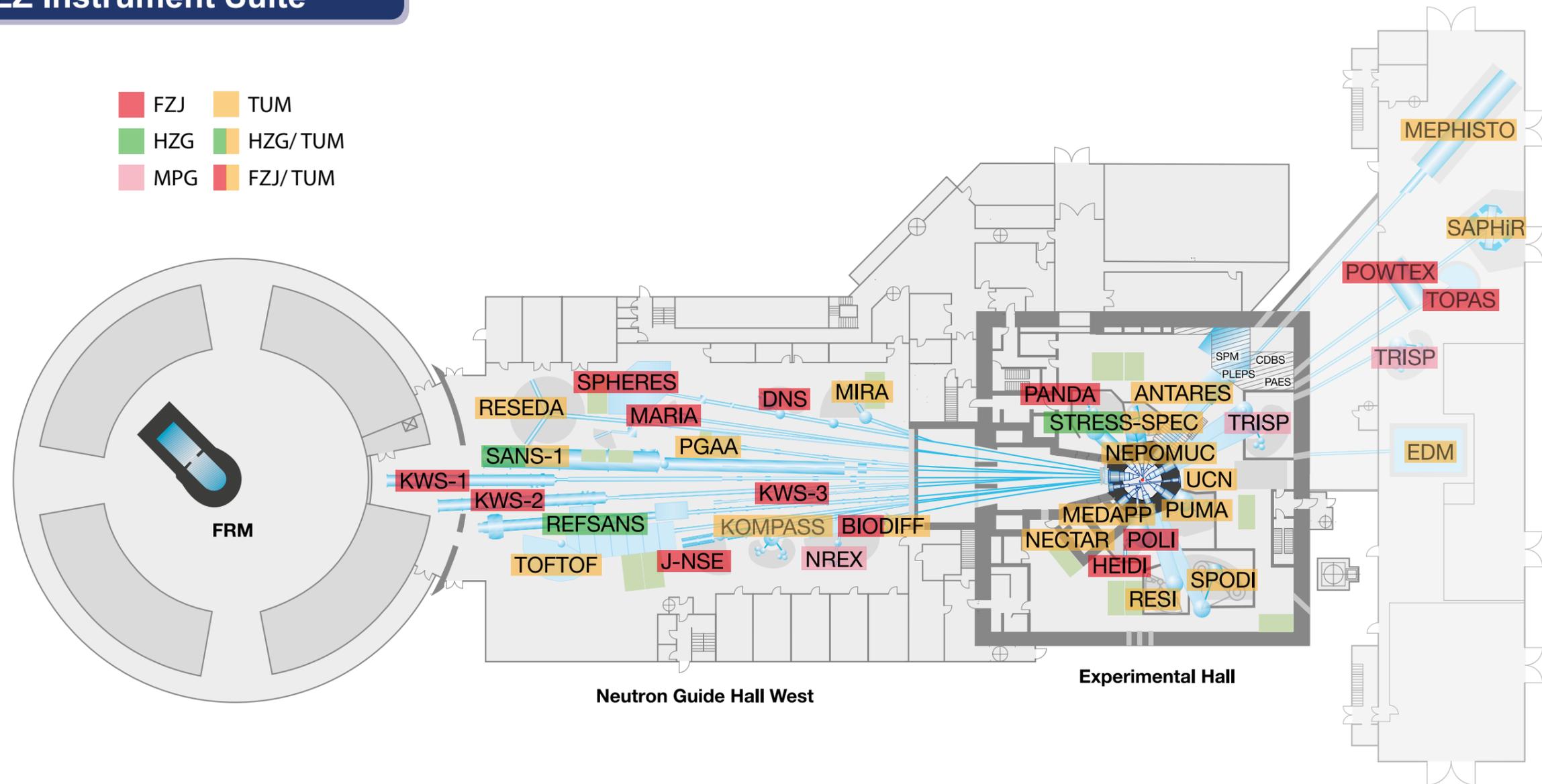

Instrument	Description	Neutrons	Status	Operated by	Funding
ANTARES	Radiography and tomography	cold	operation	TUM	TUM
BIODIFF	Diffractionmeter for large unit cells	cold	operation	TUM, JCNS	TUM, FZJ
DNS	Diffuse scattering spectrometer	cold	operation	JCNS	FZJ
HEIDI	Single crystal diffractometer	hot	operation	RWTH Aachen	FZJ
J-NSE	Spin-echo spectrometer	cold	operation	JCNS	FZJ
KOMPASS	Three axes spectrometer	cold	construction	Uni Köln, TUM	BMBF
KWS-1	Small angle scattering	cold	operation	JCNS	FZJ
KWS-2	Small angle scattering	cold	operation	JCNS	FZJ
KWS-3	Very small angle scattering	cold	operation	JCNS	FZJ
MARIA	Magnetic reflectometer	cold	operation	JCNS	FZJ
MEPHISTO	Facility for particle physics, PERC	cold	reconstruction	TUM	TUM, DFG
MIRA	Multipurpose instrument	cold	operation	TUM	TUM
MEDAPP	Medical irradiation treatment	fast	operation	TUM	TUM
NECTAR	Radiography and tomography	fast	operation	TUM	TUM
NEPOMUC	Positron source, CDBS, PAES, PLEPS, SPM	-	operation	TUM, UniBw München	TUM, BMBF
NREX	Reflectometer with X-ray option	cold	operation	MPI Stuttgart	MPG
PANDA	Three axes spectrometer	cold	operation	JCNS	FZJ

Instrument	Description	Neutrons	Status	Operated by	Funding
PGAA	Prompt gamma activation analysis	cold	operation	Uni Köln	TUM
PUMA	Three axes spectrometer	thermal	operation	Uni Göttingen, TUM	TUM
POLI	Single-crystal diffractometer polarized neutrons	hot	operation	RWTH Aachen	BMBF, FZJ
POWTEX	Time-of-flight diffractometer	thermal	construction	RWTH Aachen, Uni Göttingen, JCNS	BMBF, FZJ
REFSANS	Reflectometer	cold	operation	GEMS	HZG
RESEDA	Resonance spin-echo spectrometer	cold	operation	TUM	TUM
RESI	Single crystal diffractometer	thermal	operation	LMU	TUM
SANS-1	Small angle scattering	cold	operation	TUM, GEMS	TUM, HZG
SAPHIR	Six anvil press for radiography and diffraction	thermal	construction	BGI	BMBF
SPHERES	Backscattering spectrometer	cold	operation	JCNS	FZJ
SPODI	Powder diffractometer	thermal	operation	KIT	TUM
STRESS-SPEC	Materials science diffractometer	thermal	operation	TUM, TU Clausthal, GEMS	TUM, HZG
TOFTOF	Time-of-flight spectrometer	cold	operation	TUM	TUM
TOPAS	Time-of-flight spectrometer	thermal	construction	JCNS	FZJ
TRISP	Three axes spin-echo spectrometer	thermal	operation	MPI Stuttgart	MPG
UCN	Ultra cold neutron source, EDM	ultra-cold	construction	TUM	TUM, DFG